

# **EU PROFICIENCY TESTS EUPT-CF16, 2022**

# Pesticide Residues in Barley Kernels

**Final Report** 

Version 2

Mette Erecius Poulsen Ederina Ninga Elena Hakme

November 2022

#### Organiser:





#### **Mette Erecius Poulsen**

Director of EURL Cereals and Feedingstuff National Food Institute Research Group for Analytical Food Chemistry Technical University of Denmark Kemitorvet 1

DK-28 Lyngby Phone: +45 3588 7463 Fax: +45 3588 7448 E-Mail: eurl-cf@food.dtu.dk http://www.eurl-pesticides.eu

#### Organising Team at EURL for Pesticide Residues in Cereals and Feedingstuff:

Helen Fodnæss

Elena Hakme

Ederina Ninga

Ban M. Kadhum

Susanne Pless

Lisbet Pilhkjær

# **Quality Control Group:**

Antonio Valverde University of Almería, Spain

Paula Medina-Pastor European Food Safety Authority (EFSA), Italy

# **Advisory Group**

Amadeo R. Fernández-Alba University of Almeria, Spain

André de Kok Former WFSR - Wageningen Food Safety Research, Wageningen, The Netherlands

Bjoern Hardebush Chemisches und Veterinäruntersuchungsamt (CVUA) Freiburg, Germany

Carmen Ferrer Amate University of Almeria, Spain

Finbarr O'Regan Pesticide Control Laboratory, Co. Kildare, Ireland

Hans Mol WFSR - Wageningen Food Safety Research, Wageningen, The Netherlands

Magnus Jezussek Bavarian Authority of Health and Food Safety, Erlangen, Germany

Michelangelo Anastassiades Chemisches und Veterinäruntersuchungsamt (CVUA) Stuttgart, Germany

Patrizia Pelosi Istituto Superiore di Sanità, Roma, Italy

Ralf Lippold Chemisches und Veterinäruntersuchungsamt (CVUA) Freiburg, Germany

Hermann Unterluggauer AGES, Institute for Food Safety, Innsbruck, Austria Štěpán Radim Czech Agriculture and Food Inspection Authority

Tuija Pihlström National Food Agency, Uppsala, Sweden

# **PREFACE**

Regulation (EU) No 2017/625 [1], defines the general tasks and duties of the European Union Reference Laboratories (EURLs) for Food, Feed and Animal Health including the organisation of comparative tests. These proficiency tests (PTs) are carried out on an annual basis, and aim to improve the quality, accuracy and comparability of the analytical results generated by EU Member States within the framework of the EU multiannual co-ordinated control and national monitoring programmes. Participation in the proficiency test scheme "European Union Proficiency Tests (EUPTs) for pesticide residues" is mandatory according to Article 28 of Regulation (EC) No 396/2005 on maximum residue levels of pesticides in, or on, food and feed of plant and animal origin [2], as long as the analytical scope of the PT and the laboratory overlap.

The present EUPT was the sixteenth organized within the frame of the EURL activities with cereal or feed matrices as Test Items. The previous PTs were EUPT-C1/SRM2 on wheat, EUPT-C2 on wheat, EUPT-C3/SRM4 on hay, EUPT-C4 on rye, EUPT-C5/SRM6 on rice, EUPT-C6 on barley, EUPT-CF7 on animal feed, EUPT-CF8 on wheat, EUPT-CF9 on maize, EUPT-CF10 on rye flour, EUPT-CF11 on oat flour, EUPT-CF12 on hay flour, EUPT-CF13 on rye kernels, EUPT-CF14 on rice kernels and the EUPT-CF15 on rapeseed cake. The PTs in 2007, 2009, 2011, 2015 and 2020 were jointly organised by the EURL-CF and EURL-SRM using same cereal and focusing on both MRM and SRM pesticides. The other PTs have only focused on MRM-pesticides. The barley kernels used for EUPT-CF16 were treated both with formulations in the field and post-harvest in the laboratory.

Participation in EUPT-CF16 was compulsory for all National Reference Laboratories (NRLs) and Official Laboratories (OfLs) within the EU involved in the determination of pesticide residues in cereals for human or animal consumption using multi residue methods for their national programmes. Official laboratories from EFTA countries (Iceland, Norway and Switzerland), as well as official laboratories from EU-candidate states, were invited to take part in this EUPT. Selected laboratories from Third Countries were also allowed to take part in this exercise, but their results, together with the EU-candidate state laboratories, were not used when establishing the Assigned Values for each pesticide.

DG-SANTE will have full access to all data from EUPTs including the lab-code/lab-name key. The same will apply to all NRLs regarding data from laboratories belonging to their own country network. The results of this EUPT may be further presented to the European Commission Standing Committee for Animal Health and the Food Chain.



# **CONTENTS**

| PREFACE       |   | 7  |
|---------------|---|----|
| CONTENTS.     |   | 9  |
| 1. INTRODUC   | CTION   | 11 |
| 1.1 Analytic  | al methods  | 11 |
| 1.2 Selectio  | n of Pesticides for the Target Pesticide List                               | 11 |
| 1.3 Prepara   | tion of the Test Item   | 12 |
| 1.4 Homoge    | neity test  | 12 |
| 1.5 Stability | tests   | 13 |
| 1.6 Organisa  | ational details   | 14 |
| 2. EVALUAT    | ON OF THE RESULTS   | 15 |
| 2.1 False po  | sitives and negatives   | 15 |
| 2.2 Estimati  | on of the true concentration (x <sub>pt</sub> )                             | 15 |
| 2.3 Uncertai  | nty of the assigned value   | 15 |
| 2.4 Standar   | d deviation of the assigned value (target standard deviation)               | 15 |
| 2.5 Z scores  | 3   | 16 |
| 2.6 Categor   | y A and B classification and combined z scores (AZ²)                        | 16 |
| 3. RESULTS    |   | 17 |
| 3.1 Summai    | y of reported results   | 17 |
| 3.2 Assigne   | d values, target standard deviations and Alg A standard deviations          | 20 |
| 3.3 Assessn   | nent of laboratory performance  | 21 |
| 3.4 Trends i  | n numbers of participating laboratories and their performance               | 40 |
| 3.5 Summai    | y, conclusions and prospects for the EUPTs on pesticide residues in cereals | 41 |
| 4. ACKNOWI    | LEDGEMENTS  | 42 |
| 5. REFEREN    | CES   | 42 |
| APPENDICE     | S   | 43 |
| Appendix 1    | List of Laboratories registered to participate in the EUPT-CF14             | 43 |
| Appendix 2    | Target Pesticide List   | 47 |
| Appendix 3    | Homogeneity data  | 53 |
| Appendix 4    | Stability figures   | 55 |
| Appendix 5    | Graphical presentation of z scores  | 61 |
| Annexes       |   | 79 |
| Annex 1       | General protocol  | 79 |
| Annex 2       | Specific protocol   | 99 |

# EUROPEAN COMMISSION EURL PROFICIENCY TEST ON PESTICIDE RESIDUES IN CEREALS EUPT-CF16, 2022

#### 1. INTRODUCTION

On 1 December 2021 the announcement of the 16<sup>th</sup> European Commission's Proficiency Test on Cereals and Feed (EUPT-CF16) was published on the EURL website, together with the Calendar and the Pesticide Target List including all compounds that could potentially be present in the Test Item. The Target Pesticides List included 169 individual compulsory compounds and 53 voluntary requiring the use of multi residue methods (MRMs), along with a minimum required reporting level (MRRL) stipulated for each compound. Links to The General Protocol containing information (Annex 1) that is common to all EUPTs, the Specific protocol (Annex 2), as well as a list of labs that are obliged to take part in the EUPT-CF16, were provided via the homepage. Laboratories were able to register online from December 2021 to 28 of February 2022. In total 151 laboratories from EU and EFTA countries agreed to participate in the test as well as 9 laboratories from EU-Candidate States and Third Countries (Appendix 1).

The barley were sprayed in the field with 19 pesticides. The cultivation was performed in 2021 in Denmark by the Danish Centre for Food and Agriculture at Aarhus University. After analyses of the pesticide residues content, it was decided to additionally spike in the laboratory with 8 pesticides, which were either not included in the field treatments or where residues were too low for the evaluation.

The pesticides employed for the field treatment were selected by the EURL-CF and the EUPT quality control group. The application rates and harvest intervals chosen were based on previous experience and data from supervised residue trials. The test material was checked for homogeneity before shipping to participants. Furthermore, the stabilities of the pesticides in the Test Item were checked several times during the period of time allowed for laboratories to complete the PT exercise.

The participating laboratories were provided with 100 g portions of the rapeseed cake Test Item. The Test Items were shipped to participants on 28 March 2022 and the deadline for submission of results to the Organiser was the 2 May 2022. The deadline for submission of additional information for false negative results was the 11 May 2022. The participants were asked to analyse the Test Item and report the concentrations of any pesticide residues found that were included in the Target Pesticide List (**Appendix 2**). Submission of results was performed online via the DTU Webtool.

#### 1.1 Analytical methods

The QuEChERS method [3] was used by the organiser to test the homogeneity and stability of the Test Items. Determination was performed by GC-MS/MS and LC-MS/MS.

— QuEChERS - Citrate buffered (EN 151662): Cold water was added to a milled portion of the test item and shaken. Acetonitrile was added immediately and the tube was shaken again. A salt and buffer mixture was then added together with ceramic homogenizers and the sample was shaken vigorously for 1 min. After centrifugation, an aliquot of the supernatant was cleaned by freezing out. After additional centrifugation of the cold extract 1 ml of supernatant was filtrated and transferred in a autosampler vial for the LC/MS/MS analysis. The remaining extract supernatant was transferred to a tube containing PSA and MgSO<sub>4</sub>. After shaking and centrifugation the extract was ready for analysis by GC-MS/MS.

# 1.2 Selection of Pesticides for the Target Pesticide List

The pesticides to be included in the target pesticides list were selected by the Organiser and the Quality Control Group, taking into account the present and upcoming scope of the EU multi-annual coordinated control programme, the working document, and pesticides according to their relevance and risk-potential, as well as pesticides relevant to the specific commodity (barley kernels). The overall capacity and capability of the laboratories within the EU, as assessed from previous PTs and surveys, was also taken into account. The minimum required reporting level (MRRL) for all pesticides in the target list was in general set at 0.01 mg/kg. However, for 20 pesticides the MRRL were set at or below 0.005 mg/kg.

#### 1.3 Preparation of the Test Item

The field spraying was performed in 2021 in Denmark and organised by Danish Centre for Food and Agriculture at Aarhus University. Approximately, 40 kg of the harvested barley were used for this PT. It was decided to additionally spike in the laboratory with eight pesticides, which were either not included in the field treatments or where residues were too low for the evaluation (**Table 1**). Spiking in the laboratory was performed using formulations or pure standards. Seven portions of 1.4 kg of the non-treated barley was spiked and subsequently mixed with 20 kg of field treated barley kernels and homogenised thoroughly. One hundred gram portions of the homogenized barley kernels were then weighed out into screw-capped polyethylene plastic bottles, sealed, numbered, and stored in a freezer at about -20 °C prior to homogeneity testing and distribution to participants.

**Table 1**. Pesticides used for application in the field and/or spiked in the laboratory.

| Pesticides           | Application in field | Spike in laboratory | Formulation/standard       |
|----------------------|----------------------|---------------------|----------------------------|
| Aclonifen            |                      | X                   | Analytical Standard        |
| Azoxystrobin         | x                    |                     | Amistar                    |
| Benzovindiflupyr*    | X                    |                     | Elatus plus                |
| Boscalid             | x                    |                     | Pictor active              |
| Dieldrin             |                      | X                   | Analytical Standard        |
| Endosulfan-beta      |                      | X                   | Analytical Standard        |
| Fenpicoxamid*        | X                    |                     | Univoq/Analytical Standard |
| Fenpropathrin        |                      | X                   | Analytical Standard        |
| Fluopyram            | X                    |                     | Propulse                   |
| Fluxapyroxad         | X                    |                     | Imtrex                     |
| Lambda-cyhalothrin   | X                    |                     | Lamdex                     |
| Lindane              |                      | X                   | Analytical Standard        |
| Mefentrifluconazole* | X                    |                     | Revisol                    |
| Metconazole*         |                      | X                   | Analytical Standard        |
| Proquinazid          | X                    |                     | Talius                     |
| Prothioconazole      | X                    |                     | Proline/Propulse           |
| Pyraclostrobin       | X                    |                     | Comet Pro/Pictor active    |
| Tau-Fluvalinate      | X                    |                     | Mavrik                     |
| Teflubenzuron        |                      | X                   | Analytical standard        |

<sup>\*</sup>Voluntary pesticides.

The barley kernels were also field treated with ethephon, mepiquat, trihexapac-ethyl. These compounds were not included on the target list.

# 1.4 Homogeneity test

Ten bottles of the Test Items were randomly chosen and analyses were performed on duplicate portions taken from each bottle using the analytical methods described in section 1.1. The sequence of analyses and injections were also randomly chosen. Quantification was performed using a 5-point calibration curve constructed from matrix-matched standards.

The statistical evaluation was performed according to the International Harmonized Protocols published by IUPAC, ISO and AOAC [4]. An overview of the statistical analyses of the homogeneity test is shown in **Table 2**. The individual residues data from the homogeneity tests, as well as the results of the statistical analyses, are given in **Appendix 3**.

The homogeneity test is to show that the between-bottle variance is not greater than the within-bottle variance. The acceptance criteria to show that the Test Items were sufficiently homogeneous for the proficiency test was that:

 $S_s^2 < c$  where  $S_s$  is the between-bottle sample standard deviation and  $c = F_1 \times \sigma_{all}^2 + F_2 \times s_{an}^2$ ; F1 and F2 being constants with values of 1.83 and 0.93, respectively, from the 11 samples taken,  $\sigma_{all}^2 = 0.3 \times FFP$  RSD (25%) x the analytical sampling mean for all pesticides, and  $s_{an}$  is the estimate of the analytical standard deviation.

As all pesticides passed the homogeneity test, when the Test Item was stored at -18 °C, the Test Item was considered to be sufficiently homogeneous and suitable for the EUPT-CF16.

**Table 2**. Statistical evaluation of the homogeneity test data (n=22 analyses using a sub-sample of 5 g in each case).  $S_s$ : Between Sampling Standard Deviation.

| Pesticides              | Mean, mg/kg | Ss2     | c      | Ss2 < c |
|-------------------------|-------------|---------|--------|---------|
| Aclonifen               | 0.058       | 0.00000 | 0.0001 | Pass    |
| Azoxystrobin            | 1.24        | 0.00232 | 0.0234 | Pass    |
| Benzovindiflupyr        | 0.885       | 0.00126 | 0.0114 | Pass    |
| Boscalid                | 1.92        | 0.00566 | 0.0627 | Pass    |
| Dieldrin                | 0.025       | 0.00000 | 0.0000 | Pass    |
| Endosulfan-beta         | 0.034       | 0.00001 | 0.0000 | Pass    |
| Fenpropathrin           | 0.089       | 0.00004 | 0.0001 | Pass    |
| Fenpicoxamid            | 0.036       | 0.00000 | 0.0000 | Pass    |
| Fluopyram               | 1.41        | 0.00065 | 0.0306 | Pass    |
| Fluxapyroxad            | 0.730       | 0.00035 | 0.0089 | Pass    |
| Lambda-cyhalothrin      | 0.074       | 0.00002 | 0.0001 | Pass    |
| Lindane                 | 0.037       | 0.00000 | 0.0000 | Pass    |
| Mefentrifluconazole     | 0.750       | 0.00063 | 0.0089 | Pass    |
| Metconazole             | 0.035       | 0.00001 | 0.0000 | Pass    |
| Proquinazid             | 0.211       | 0.00012 | 0.0063 | Pass    |
| Prothioconazole-desthio | 1.22        | 0.00059 | 0.0245 | Pass    |
| Pyraclostrobin          | 4.68        | 0.02390 | 0.3342 | Pass    |
| Tau-Fluvalinate         | 0.327       | 0.00063 | 0.0016 | Pass    |
| Teflubenzuron           | 0.043       | 0.00002 | 0.0000 | Pass    |

#### 1.5 Stability tests

The analytical methods described briefly above (in section 1.1) were also used for the stability tests.

The stability test was performed according to ISO 13528, Annex B [5]. Two different storage temperatures were used; room temperature and -18 °C. Six sub-samples (analytical portions) were analysed on each test day. A pesticide is considered to be adequately stable if  $|x_1 - y_i| \le 0.3 \times \sigma$ , where  $x_1$  is the mean value of the first stability test,  $y_i$  the mean value of the last stability test and  $\sigma$  the standard deviation used for proficiency assessment (25% of the assigned value):

The dates of testing were as follows:

Day 1: 28 March 2022 Day 2: 13 April 2022 Day 3: 2 May 2022

The results of the stability test for storage temperature -18 °C are given in Table 3. All pesticides passed the test at -18 °C. At room temperature aclonifen, dieldrin, fenpicoxamid, fenpropathrin, lindane, metconazole and teflubenzuron did not pass the test. All the laboratories were instructed to store the test item at -18 degree and the stability test was consequently accepted. See the individual stability figures for all pesticides in **Appendix 4**.

Table 3. Statistical evaluation of the stability test data at -18 °C.

| Pesticides              | Mean, mg/kg | x1 - yi | 0.3×σ  | x1 - yi   ≤ 0.3×σ |
|-------------------------|-------------|---------|--------|-------------------|
| Aclonifen               | 0.060       | 0.004   | 0.004  | Pass              |
| Azoxystrobin            | 1.09        | 0.042   | 0.091  | Pass              |
| Benzovindiflupyr        | 0.937       | 0.055   | 0.063  | Pass              |
| Boscalid                | 2.63        | 0.106   | 0.164  | Pass              |
| Dieldrin                | 0.023       | 0.002   | 0.002  | Pass              |
| Endosulfan-beta         | 0.030       | 0.0007  | 0.0030 | Pass              |
| Fenpicoxamid            | 0.036       | 0.003   | 0.004  | Pass              |
| Fenpropathrin           | 0.074       | 0.007   | 0.007  | Pass              |
| Fluopyram               | 1.61        | 0.093   | 0.122  | Pass              |
| Fluxapyroxad            | 0.809       | 0.049   | 0.068  | Pass              |
| Lambda-cyhalothrin      | 0.080       | 0.003   | 0.005  | Pass              |
| Lindane                 | 0.027       | 0.002   | 0.003  | Pass              |
| Mefentrifluconazole     | 0.805       | 0.056   | 0.063  | Pass              |
| Metconazole             | 0.031       | 0.003   | 0.003  | Pass              |
| Proquinazid             | 0.223       | 0.016   | 0.020  | Pass              |
| Prothioconazole-desthio | 1.71        | 0.111   | 0.122  | Pass              |
| Pyraclostrobin          | 5.04        | 0.372   | 0.383  | Pass              |
| Tau-Fluvalinate         | 0.255       | 0.022   | 0.024  | Pass              |
| Teflubenzuron           | 0.032       | 0.005   | 0.001  | Pass              |

#### 1.6 Organisational details

#### 1.6.1 Access to documents, registration and confidentiality

In the invitation letter, all NRLs and OfLs were requested to register using the online registration link from December 2021. All documents related to this EUPT (Calendar, Target Pesticides List, Specific Protocol, General Protocol) were uploaded to the EURL website and the CIRCA platform. Laboratories that were intending not to participate were given the opportunity to explain the reasons for their non-participation. Participants from Candidate countries and third countries did also have access to another online registration link. On 14 March, the participants received a link to DTU web tool, along with login credentials and were asked to enter the web tool and to select the scope of pesticides they wanted to be evaluated on. This had to be done before the samples were shipped to the participants.

#### 1.6.2 Distribution of the Test Item

On 28 March 2022, the Test Item ( $100 \, g$ ) was shipped to all participants in insulated polystyrene boxes containing a freezer block. The laboratories were asked to check the state of the sample on receipt and to enter the web tool to report whether they accept/not accept the Test Item. No blank test material was send.

#### 1.6.3 Submission of results

The participants had to submit their results via a web tool. All participants had access to the result-submission website from a few days after shipment until the result-submission deadline (11 May 2022). Participants were asked not only to report their analytical results, but also to give information regarding accreditation, reporting limits and details regarding the methods they used to analyse the Test Item.

#### 2. EVALUATION OF THE RESULTS

The results were evaluated according to the general and specific protocols (Annex 1 and 2). However, the main points are listed below.

#### 2.1 False positives and negatives

#### 2.1.1 False positives

These are results of pesticides from the Target Pesticides List, that are reported at or above, their respective MRRLs although they were: (i) not detected by the Organiser, even after repeated analyses, and/or (ii) not detected by the overwhelming majority (e.g. > 95%) of the participating laboratories that had targeted these specific pesticides. In certain instances, case-by-case decisions by the EUPT-Panel may be necessary. Any results reported lower than the MRRL will not be considered as false positives, even though these results should not have been reported.

#### 2.1.2 False negatives

These are results for pesticides reported by the laboratories as 'analysed' but without reporting numerical values although they were: a) used by the Organiser to treat the Test Item and b) detected by the Organiser as well as the majority of the participants that had targeted these specific pesticides at, or above the respective MRRLs. Results reported as '< RL' (RL= Reporting Limit of the laboratory) will be considered as not detected and will be judged as false negatives. In certain instances, case-by-case decisions by the EUPT-Panel may be necessary. In cases of the assigned value being less than a factor of 3 times the MRRL, false negatives will typically not be assigned. The EUPT-Panel may decide to take case-by-case decisions in this respect after considering all relevant factors such as the result distribution and the reporting limits of the affected labs.

#### 2.2 Estimation of the true concentration (xpt)

In order to minimise the influence of out-lying results on the statistical evaluation, the assigned value xpt (= consensus concentration) will typically be estimated using robust estimate of the participants' mean  $(x^*)$  as described in ISO 13528:2015, taking into account the results reported by only EU and EFTA countries laboratories. In special justifiable cases, the EUPT-Panel may decide to eliminate certain results traceably associated with gross errors, or to use only the results of a subgroup consisting of laboratories that have repeatedly demonstrated good performance for the specific compound in the past.

# 2.3 Uncertainty of the assigned value

The uncertainty of the assigned values u(xpt) is calculated according to ISO 13528:2015 as:

$$u(xpt) = 1.25 \frac{s*}{\sqrt{p}}$$

where  $s^*$  is the robust standard deviation and p is the number of results.

#### 2.4 Standard deviation of the assigned value (target standard deviation)

The target standard deviation of the assigned value (FFP- $\sigma_{pt}$ ) will be calculated using a Fit-For-Purpose approach with a fixed Relative Standard Deviation (FFP-RSD) of 25% as follows:

$$FFP-\sigma_{pt} = 0.25 *x_{pt}$$

The percentage FFP-RSD is set at 25% based on experience from results of previous EUPTs. The EUPT-Panel reserves the right to also employ other approaches on a case-by-case basis considering analytical difficulties and experience gained from previous proficiency tests.

For informative purposes the robust relative standard deviation (CV\*) is calculated according to ISO 13528:2015; Chapter 7.7 (Consensus value from participant results) following Algorithm A in Annex C.

#### 2.5 Z scores

A z-score for each laboratory/pesticide combination was calculated according to the following equation:

$$z_i = \frac{\left(x_i - x_{pt}\right)}{FFP - \sigma_{pt}}$$

where  $x_i$  is the value reported by the laboratory,  $x_{pt}$  is the assigned value, and FFP- $\sigma_{pt}$  is the standard deviation using FFP approach. Z scores was rounded to one decimal place. For the calculation of combined z scores (see below) the original z scores will be used and rounded to one decimal place after calculation.

Any z scores > 5 will be typically reported as '> 5' and a value of '5' will be used to calculate combined z scores.

Z scores will be interpreted in the following way as is set in the ISO 17043:2010 [6]:

 $|z| \le 2$  Acceptable 2 < |z| < 3 Questionable  $|z| \ge 3$  Unacceptable

For results considered as false negatives, z scores will be calculated using the MRRL or RL (the laboratory's Reporting Limit) if RL < MRRL. Where, using this approach, the calculated z scores for false negatives are > -3 (still questionable), they will be fixed at -3.5 to underline that these are unacceptable results. These z-scores will typically appear in the z-score histograms and used in the calculation of combined z-scores.

#### 2.6 Category A and B classification and combined z scores (AZ2)

The EUPT-Panel will decide if and how to classify the laboratories into two categories - A or B. Currently, laboratories that are able to analyse at least 90% of the compulsory pesticides in the target pesticides list, have correctly detected and quantified a sufficiently high percentage of the pesticides present in the Test Item (at least 90%) and reported no false positives, will have demonstrated 'sufficient scope' and can therefore be classified into Category A. For the 90% criteria, the number of pesticides needed to be correctly analysed to have sufficient scope will be calculated by multiplying the number of compulsory pesticides from the Target Pesticides List by 0.9 and rounding to the nearest full number with 0.5 decimals being rounded downwards.

For evaluation of the overall performance of laboratories within Category A, the Average of the Squared z Score ( $AZ^2$ ) will be used. The  $AZ^2$  is calculated as follows:

$$AZ^{2} = \frac{\sum_{i=1}^{n} Z_{i}^{2}}{n}$$

where "n" is the number of each laboratory's z scores that were considered in this formula. For the calculation, any z-score > 5 was set at "5". Based on the AZ<sup>2</sup> achieved, the laboratories are classified as follows:

 $AZ^2 \le 2$ Good $2 < AZ^2 < 3$ Satisfactory $AZ^2 \ge 3$ Unsatisfactory

The AZ<sup>2</sup> is considered being of lesser importance than the individual z scores.

Laboratories within Category B are ranked according to the total number of pesticides that they correctly reported to be present in the Test Item. The number of acceptable z scores achieved is listed as well.

#### 3. RESULTS

#### 3.1 Summary of reported results

In total, 151 EU and EFTA laboratories, from 29 different countries (26 EU member states), agreed to participate in this proficiency test. Five EU participants did not submit results. Additionally, nine participants from non-EU Countries registered for the PT. The participating laboratories are listed in **Appendix 1**.

An overview of results submitted by laboratories from the EU and EFTA can be seen in **Table 4**. All reported analytical results for the pesticide residues are shown in **Table 10a-c** and in **Appendix 5**. However, only results submitted by laboratories from EU and EFTA countries are included in **Table 4**, **8-9** and **13** and the z scores histograms are shown in **Appendix 5**.

**Table 4.** Overview of number of results, number of not analysed (NA), number of not detected (ND = false negatives) and the percentage of laboratories that reported results for the pesticides in the Test Item. Only results submitted by laboratories from the EU and EFTA are included in this table.

| Pesticides              | No. of reported No. of NA results |     | False negatives | % of labs reporting results <sup>1</sup> |  |
|-------------------------|-----------------------------------|-----|-----------------|--|--|
| Aclonifen               | 97                                | 51  | 9               | 66                                       |  |
| Azoxystrobin            | 139                               | 9   | 0               | 94                                       |  |
| Benzovindiflupyr        | 67                                | 81  | 4               | 45                                       |  |
| Boscalid                | 135                               | 13  | 0               | 91                                       |  |
| Dieldrin                | 140                               | 8   | 12              | 95                                       |  |
| Endosulfan-beta         | 143                               | 5   | 6               | 97                                       |  |
| Fenpicoxamid            | 48                                | 100 | 7               | 32                                       |  |
| Fenpropathrin           | 127                               | 21  | 5               | 86                                       |  |
| Fluopyram               | 132                               | 16  | 0               | 89                                       |  |
| Fluxapyroxad            | 122                               | 26  | 0               | 82                                       |  |
| Lambda-cyhalothrin      | 136                               | 12  | 1               | 92                                       |  |
| Lindane                 | 139                               | 9   | 7               | 94                                       |  |
| Mefentrifluconazole     | 50                                | 98  | 7               | 34                                       |  |
| Metconazole             | 125                               | 23  | 8               | 84                                       |  |
| Proquinazid             | 112                               | 36  | 0               | 76                                       |  |
| Prothioconazole-desthio | 122                               | 26  | 4               | 82                                       |  |
| Pyraclostrobin          | 126                               | 22  | 1               | 85                                       |  |
| Tau-Fluvalinate         | 131                               | 17  | 1               | 89                                       |  |
| Teflubenzuron           | 115                               | 33  | 21              | 78                                       |  |

 $<sup>^1</sup>$  '% results' have been calculated using the number of laboratories that reported results for each particular compound and the total number of EU laboratories that submitted results (n = 148). False negatives are included in reported results.

Azoxystrobin, boscalid, dieldrin, endosulfan-beta, lambda-cyhalothrin and lindane were the most frequently analysed compounds with  $\geq 90$  % of the labs submitting results for these compounds. Fenpropathrin, fluopyram, fluxapyroxad, metconazole, proquinazid, prothioconazole-desthio, pyraclostrobin, tau-fluvalinate and teflubenzuron were analysed and reported by 78-86% of the participants. Aclonifen, benzovindiflupyr, fenpicoxamid and mefentrifluconazole were only analysed and reported by 34-66% of participants.

#### 3.1.1 False positives

Fourteen participants (all form EU and EFTA) countries reported twenty-five results for fifteen different additional pesticides above the MRRL that had not been used to treat the Test Item (**Table 5**). The pesticides were: Acrinathrin,

biphenyl, buprofezin, chlorfluazuron, cyprodinil, endosulfan alpha, endosulfan sulfate, endrin, etoxazole, fenbuconazole, HCH-alpha, HCH-beta, hexaconazole, prosulfocarb, prothiofos, thiabendazole and trifloxystrobin. In all cases the compounds were not detected either by the Organizer, or by the other participating laboratories. The reported results were therefore considered to be false positives.

Table 5. False positive results at or above 0.01 mg/kg, the concentration detected in mg/kg, the determination

| Lab code | Pesticides         | Concentration<br>mg/kg | Determination technique | RL, mg/kg | MRRL, mg/kg |
|----------|--------------------|------------------------|-------------------------|-----------|-------------|
| 35       | Acrinathrin        | 0.056                  | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 91       | Acrinathrin        | 0.029                  | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 112      | Acrinathrin        | 0.07                   | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 162      | Acrinathrin        | 0.028                  | GC-Q-Orbitrap           | 0.01      | 0.01        |
| 154      | Biphenyl           | 0.022                  |                         | 0.01      | 0.01        |
| 128      | Buprofezin         | 0.0105                 | LC-MS/MS QQQ            | 0.01      | 0.01        |
| 128      | Chlorfluazuron     | 0.04                   | LC-MS/MS QQQ            | 0.01      | 0.01        |
| 135      | Cyprodinil         | 0.048                  | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 154      | Endosulfan alpha   | 0.02                   |                         | 0.01      | 0.01        |
| 162      | Endosulfan alpha   | 0.023                  | GC-Q-Orbitrap           | 0.01      | 0.01        |
| 72       | Endosulfan sulfate | 0.03                   | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 9        | Endrin             | 0.017                  | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 153      | Endrin             | 0.0211                 | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 112      | Etoxazole          | 0.027                  | LC-MS/MS QQQ            | 0.01      | 0.01        |
| 112      | Fenbuconazole      | 0.014                  | LC-MS/MS QQQ            | 0.01      | 0.005       |
| 97       | HCH-alpha          | 0.0266                 | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 5        | HCH-beta           | 0.031                  | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 135      | HCH-beta           | 0.033                  | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 157      | HCH-beta           | 0.0379                 | GC-MS/MS (QQQ)          | 0.01      | 0.01        |
| 73       | Hexaconazole       | 0.99                   | LC-MS/MS QQQ            | 0.01      | 0.01        |
| 162      | Hexaconazole       | 1.24                   | LC-MS/MS QQQ            | 0.01      | 0.01        |
| 112      | Prosulfocarb       | 0.017                  | LC-MS/MS QQQ            | 0.01      | 0.01        |
| 154      | Prothiofos         | 0.775                  |                         | 0.01      | 0.01        |
| 154      | Thiabendazole      | 0.02                   |                         | 0.01      | 0.01        |
| 112      | Trifloxystrobin    | 0.016                  | LC-MS/MS QQQ            | 0.01      | 0.01        |

# 3.1.2 Findings of compounds below the MRRL mg/kg

Apart from the false positive results above and the results for the pesticide residues listed in Table 10, two participants reported results for three other pesticide residues, see Table 6. These results were not evaluated as false positives because the concentrations were below the MRRL at  $0.01\ mg/kg$ .

Table 6. Reported results in mg/kg at or below the MRRL

| Lab code | Pesticides   | Concentration<br>mg/kg | Determination<br>technique | RL, mg/kg | MRRL, mg/kg |
|----------|--------------|------------------------|----------------------------|-----------|-------------|
| 10       | Ethoprophos  | 0.0015                 | GC-MS/MS (QQQ)             | 0.1       | 0.01        |
| 112      | Fluopicolide | 0.007                  | LC-MS/MS QQQ               | 0.1       | 0.01        |
| 10       | Prosulfocarb | 0.0012                 | LC-MS/MS QQQ               | 0.1       | 0.01        |

# 3.1.3 False negatives

Not reported results for pesticides actually present in the Test Item were judged as false negatives. **Table 7** summarizes the number of reported false negatives for each pesticide. Forty-seven participants submitted 69 false negatives results for 13 different pesticides, which represents 3.3% of the total number of results submitted by EU and EFTA laboratories. Around 30% of the EU and EFTA participants (45 laboratories) reported false negative results...

**Table 7.** False negative results (FN)

| Table 7. False negative results (FN). |           |                  |          |                 |              |               |                        |         |                          |             |                             |                |                 |
|---------------------------------------|-----------|------------------|----------|-----------------|--------------|---------------|------------------------|---------|--------------------------|-------------|-----------------------------|----------------|-----------------|
| Lab code                              | Aclonifen | Benzovindiflupyr | Dieldrin | Endosulfan beta | Fenpicoxamid | Fenpropathrin | Lambda-<br>cyhalothrin | Lindane | Mefentriflu-<br>conazole | Metconazole | Prothioconazole-<br>desthio | Pyraclostrobin | Tau-Fluvalinate |
| 5                                     |           |                  |          |                 |              |               |                        | FN      |                          |             |                             |                |                 |
| 9                                     |           |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 23                                    |           |                  |          |                 |              | FN            |                        |         |                          |             |                             |                |                 |
| 28                                    |           |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 36                                    |           |                  |          |                 | FN           |               |                        |         |                          |             |                             |                |                 |
| 41                                    |           |                  | FN       | FN              |              |               |                        | FN      |                          |             |                             |                |                 |
| 44                                    |           |                  | FN       | FN              |              |               |                        |         |                          |             |                             |                |                 |
| 47                                    |           |                  | FN       | FN              |              |               |                        |         |                          |             |                             |                |                 |
| 49                                    |           |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 52                                    | FN        |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 53                                    |           |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 54                                    |           |                  |          |                 |              |               |                        |         |                          |             |                             | FN             |                 |
| 61                                    | FN        |                  |          |                 |              |               |                        |         |                          |             |                             |                |                 |
| 66                                    |           |                  |          |                 | FN           |               |                        |         |                          |             |                             |                |                 |
| 69                                    |           |                  | FN       | FN              |              |               |                        |         |                          |             |                             |                | FN              |
| 70                                    |           |                  |          |                 |              |               |                        | FN      |                          |             |                             |                |                 |
| 71                                    |           |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 72                                    |           |                  |          | FN              |              |               |                        |         |                          |             |                             |                |                 |
| 73                                    |           |                  |          |                 |              | FN            |                        |         |                          |             |                             |                |                 |
| 74                                    |           |                  |          |                 |              |               |                        |         |                          |             | FN                          |                |                 |
| 87                                    |           |                  |          |                 |              |               |                        |         |                          |             | FN                          |                |                 |
| 91                                    | FN        |                  |          |                 |              |               |                        |         |                          |             | FN                          |                |                 |
| 93                                    |           |                  | FN       |                 |              |               |                        |         |                          |             |                             |                |                 |
| 95                                    |           | FN               |          |                 | FN           |               |                        |         | FN                       |             |                             |                |                 |
| 97                                    |           |                  |          |                 |              |               |                        | FN      |                          |             |                             |                |                 |
| 99                                    | FN        |                  |          |                 |              |               |                        |         |                          | FN          |                             |                |                 |
| 101                                   |           |                  |          |                 |              |               |                        |         | FN                       |             |                             |                |                 |
| 103                                   |           | FN               |          |                 | FN           |               |                        |         | FN                       |             |                             |                |                 |
| 105                                   | FN        |                  | FN       | FN              |              | FN            |                        |         |                          |             |                             |                |                 |
| 106                                   |           |                  |          |                 |              |               |                        |         |                          | FN          |                             |                |                 |
| 109                                   | FN        |                  |          |                 |              |               |                        |         |                          |             |                             |                |                 |
| 112                                   |           |                  |          |                 |              |               |                        |         |                          | FN          |                             |                |                 |
| 119                                   |           |                  |          |                 |              |               |                        | FN      | FN                       | FN          |                             |                |                 |

| Lab code | Aclonifen | Benzovindiflupyr | Dieldrin | Endosulfan beta | Fenpicoxamid | Fenpropathrin | Lambda-<br>cyhalothrin | Lindane | Mefentriflu-<br>conazole | Metconazole | Prothioconazole-<br>desthio | Pyraclostrobin | Tau-Fluvalinate |
|----------|-----------|------------------|----------|-----------------|--------------|---------------|------------------------|---------|--------------------------|-------------|-----------------------------|----------------|-----------------|
| 120      |           |                  |          |                 |              |               |                        |         | FN                       |             |                             |                |                 |
| 126      |           |                  |          |                 |              |               |                        |         |                          | FN          |                             |                |                 |
| 128      |           |                  |          |                 |              |               | FN                     |         |                          |             |                             |                |                 |
| 129      |           |                  |          |                 |              |               |                        |         |                          | FN          |                             |                |                 |
| 135      |           |                  |          |                 |              |               |                        | FN      |                          | FN          |                             |                |                 |
| 138      |           |                  |          |                 |              | FN            |                        |         |                          |             |                             |                |                 |
| 141      | FN        |                  |          |                 |              |               |                        |         |                          |             |                             |                |                 |
| 142      | FN        |                  |          |                 |              |               |                        |         |                          |             |                             |                |                 |
| 145      |           |                  |          |                 |              |               |                        |         |                          | FN          |                             |                |                 |
| 153      | FN        |                  |          |                 |              |               |                        |         |                          |             |                             |                |                 |
| 154      |           |                  |          |                 |              | FN            |                        |         |                          |             |                             |                |                 |
| 157      |           |                  |          |                 |              |               |                        | FN      |                          |             |                             |                |                 |
| 158      |           |                  |          |                 | FN           |               |                        |         | FN                       |             |                             |                |                 |
| 162      |           |                  | FN       | FN              |              |               |                        |         |                          |             | FN                          |                |                 |

# 3.2 Assigned values, target standard deviations and Alg A standard deviations

# 3.2.1 Assigned values

The Assigned Values were calculated as the Algorithm A mean (Alg A mean), including the reported results submitted by laboratories from EU and EFTA countries.

All assigned values for the pesticides can be seen in **Table 8**. For the evaluated pesticides the assigned values were in the range of 0.023-4.39 mg/kg. The calculated Alg A mean for teflubenzuron was less than three times the MRRL. Consequently, the results for this compound could not be evaluated and the values for the compounds are given for informative purposes only

The uncertainty of the assigned values is calculated according to ISO 13528 [5] as:

$$\mu = 1.25 \; \frac{s *}{\sqrt{n}}$$

Where  $s^*$  is the robust standard deviation estimate and n is the number of datapoints equal to the number of results used to calculate the assigned value (number of results in **Table 9**)

# 3.2.2 Target standard deviations and Alg A standard deviations

The target standard deviation was obtained using a fixed FFP-RSD value of 25%. In parallel, the Algorithm A standard deviation (Alg A-RSD) was calculated for informative purposes only. The range of Alg A-RSD values was for the evaluated pesticide in the range of 17-27 % but on average, the Alg A-RSD was 22%, lower than 25% FFP-RSD used for the z score calculations.

**Table 8.** Assigned values and their uncertainty in mg/kg, Fit-For-Purpose Relative Standard Deviation (FFP RSD) and Robust Relative Standard Deviation (Alg A RSD) for the pesticides present in the Test Item.

| Pesticides              | MRRL,<br>mg/kg | Assigned<br>value, mg/kg | Uncertainty,<br>mg/kg | FFP RSD,<br>% | Alg A RSD,<br>% |
|-------------------------|----------------|--------------------------|-----------------------|---------------|-----------------|
| Aclonifen               | 0.01           | 0.046                    | 0.001                 | 25            | 19              |
| Azoxystrobin            | 0.01           | 1.04                     | 0.029                 | 25            | 27              |
| Benzovindiflupyr        | 0.01           | 0.718                    | 0.022                 | 25            | 20              |
| Boscalid                | 0.01           | 1.88                     | 0.049                 | 25            | 24              |
| Dieldrin                | 0.01           | 0.023                    | 0.001                 | 25            | 21              |
| Endosulfan-beta         | 0.01           | 0.034                    | 0.001                 | 25            | 20              |
| Fenpicoxamid            | 0.01           | 0.040                    | 0.001                 | 25            | 19              |
| Fenpropathrin           | 0.01           | 0.084                    | 0.002                 | 25            | 19              |
| Fluopyram               | 0.01           | 1.40                     | 0.033                 | 25            | 22              |
| Fluxapyroxad            | 0.01           | 0.774                    | 0.019                 | 25            | 22              |
| Lambda-cyhalothrin      | 0.01           | 0.061                    | 0.002                 | 25            | 24              |
| Lindane                 | 0.01           | 0.031                    | 0.001                 | 25            | 19              |
| Mefentrifluconazole     | 0.01           | 0.722                    | 0.032                 | 25            | 25              |
| Metconazole             | 0.01           | 0.035                    | 0.001                 | 25            | 26              |
| Proquinazid             | 0.01           | 0.225                    | 0.006                 | 25            | 24              |
| Prothioconazole-desthio | 0.01           | 1.36                     | 0.038                 | 25            | 24              |
| Pyraclostrobin          | 0.01           | 4.39                     | 0.131                 | 25            | 27              |
| Tau-Fluvalinate         | 0.01           | 0.269                    | 0.007                 | 25            | 24              |
| Teflubenzuron           | 0.01           | 0.016                    | 0.000                 | 25            | 22              |

 $<sup>^{1}</sup>$  The assigned values are less than 3 times the MRRL and consequently shown for informative purposes only.

#### 3.3 Assessment of laboratory performance

# 3.3.1 Z scores

Z scores have been calculated for all the quantified pesticides using the FFP RSD of 25%. **Table 9** shows an overview of the acceptable, questionable, and unacceptable z scores and **Tables 10a/c**-show the individual results and z scores for each laboratory and pesticide together with the assigned values. A graphical representation of the z scores (for EU and EFTA countries) can be seen in **Appendix 5**.

Of the reported results for the evaluated pesticides, more than 90% were acceptable azoxystrobin, benzovindiflupyr, boscalid, endosulfan-beta, fenpropathrin, fluopyram, fluxapyroxad, proquinazid, pyraclostrobin and tau-fluvalinate. For aclonifen, dieldrin, fenpicoxamid, lambda-cyhalothrin, lindane, mefentrifluconazole, metconazole, prothioconazole-desthio and teflubenzuron between 73-89% of the results were acceptable.

Table 9. Number of acceptable, questionable, unacceptable z scores, and false negatives.

| Pesticides              | No. of<br>reported<br>results | Assigned<br>values | Acceptable<br>% | Questionable<br>% | Unacceptable <sup>1</sup><br>% | False<br>negatives<br>% |
|-------------------------|-------------------------------|--------------------|-----------------|-------------------|--------------------------------|-------------------------|
| Aclonifen               | 97                            | 0.046              | 86              | 3                 | 11                             | 9                       |
| Azoxystrobin            | 139                           | 1.04               | 95              | 4                 | 1                              | 0                       |
| Benzovindiflupyr        | 67                            | 0.718              | 91              | 3                 | 6                              | 6                       |
| Boscalid                | 135                           | 1.88               | 93              | 4                 | 3                              | 0                       |
| Dieldrin                | 140                           | 0.023              | 86              | 2                 | 11                             | 9                       |
| Endosulfan-beta         | 143                           | 0.034              | 93              | 1                 | 6                              | 4                       |
| Fenpicoxamid            | 48                            | 0.040              | 73              | 4                 | 23                             | 15                      |
| Fenpropathrin           | 127                           | 0.084              | 92              | 3                 | 5                              | 4                       |
| Fluopyram               | 132                           | 1.40               | 91              | 5                 | 5                              | 0                       |
| Fluxapyroxad            | 122                           | 0.774              | 93              | 4                 | 3                              | 0                       |
| Lambda-cyhalothrin      | 136                           | 0.061              | 87              | 8                 | 5                              | 1                       |
| Lindane                 | 139                           | 0.031              | 89              | 4                 | 7                              | 5                       |
| Mefentrifluconazole     | 50                            | 0.722              | 78              | 8                 | 14                             | 14                      |
| Metconazole             | 125                           | 0.035              | 82              | 3                 | 14                             | 6                       |
| Proquinazid             | 112                           | 0.225              | 91              | 7                 | 2                              | 0                       |
| Prothioconazole-desthio | 122                           | 1.36               | 88              | 4                 | 8                              | 3                       |
| Pyraclostrobin          | 126                           | 4.39               | 93              | 6                 | 2                              | 1                       |
| Tau-Fluvalinate         | 131                           | 0.269              | 91              | 5                 | 5                              | 1                       |
| Teflubenzuron           | 115                           | 0.016              | 87              | 3                 | 25                             | 21                      |

<sup>&</sup>lt;sup>1</sup> Unacceptable z scores includes false negative results.

#### 3.3.2 Analytical methods used

More than five different analytical methods have been used by the laboratories. For the majority of the results, 73%, QuEChERS, Citrate buffered (EN 151662) was used. However, variations in the clean-up procedures were reported by the labs, e.g. some used a freezing out step (28% of the participants), centrifugation (14%), some used d-SPE with PSA/MgSO4 (34%), some used d-SPE with ODS/ MgSO4 (4.1%) and other used different combination of ODS, PSA, C18, z-sep (10%). Liquid-liquid partition was used by 6.2% of the participants. Consequently, it was not one specific method.

Other extraction methods have been used; the original QuEChERS version method (J. AOAC 86, 2003) and QuEChERS-Acetate buffered (AOAC Official method 2007.01) were used by 5.5% and 7.7% of the laboratory, respectively. The Mini-Luke method and the SweEt method were each used by 3.3% of the participants. The remaining 7.7% of the participants used other methods. More than 93% of the reported results derived from a method where water was added before extraction.

For milling, 53% of the labs used a knife mill and 27% of the labs used centrifugal mill. Moreover, 7.4% used a disk mill, 2.1% used an horizontal mill, and 0.6% used a hammer mill. Furthermore, 9.2% of the labs did not specify the type of mill used or did not mill at all.

GC instruments was used for 47% of the results, mainly GC-MS/MS (82%), but also GC-MS (7.5%) and GC-  $(\mu)$  ECD (7.5%) was used. GC-NPD was used for 1 results, GC-iontrap for 7, GC-TOF and GC-Orbitrap for 15 result. LC instruments was used for 53% of the reported results, mainly LC-MS/MS (91%) but 6.8% used high resolution instrument like LC-Orbitrap, LC-Q-Orbitrap or LC-Q-TOF. Finally, 13 results were based on LC-Iontrap. No result were analysed using specific detectors such as LC-Fluorescence, LC-UV, or LC-DAD.

 $\textbf{Table 10a.} \ \ \text{Results for the mandatory pesticides aclonifen, azoxystrobin, boscalid, dieldrin, endosulfan-beta, fenpropathrin, fluopyram and fluxapyroxad in mg/kg, the corresponding z scores, MRRLs and the assigned values.$ 

| values.         |           |                         |               |                         |              |                         |             |                         |                 |                         |               |                         |               |                         |              |                         |
|-----------------|-----------|-------------------------|---------------|-------------------------|--------------|-------------------------|-------------|-------------------------|-----------------|-------------------------|---------------|-------------------------|---------------|-------------------------|--------------|-------------------------|
| Laboratory code | Aclonifen | Z-scores (FFP RSD (25%) | Azoxystrobin  | Z-scores (FFP RSD (25%) | Boscalid     | Z-scores (FFP RSD (25%) | Dieldrin    | Z-scores (FFP RSD (25%) | Endosulfan-beta | Z-scores (FFP RSD (25%) | Fenpropathrin | Z-scores (FFP RSD (25%) | Fluopyram     | Z-scores (FFP RSD (25%) | Fluxapyroxad | Z-scores (FFP RSD (25%) |
| MRRL            | 0.01      | res                     | 0.01          | res                     | 0.01         | res                     | 0.01        | Sal                     | 0.01            | res                     | 0.01          | Tes                     | 0.01          | res                     | 0.01         | res                     |
| Assigned        |           | -800                    |               | 008-                    |              | oos-                    |             | 008-                    |                 | oos-                    |               | oos-                    |               | -SC0                    |              | 008-                    |
| value           | 0.046     | Z                       | 1.04          | Z                       | 1.88         | Z                       | 0.023       | Z                       | 0.034           | Z                       | 0.084         | Z                       | 1.40          | Z                       | 0.774        | Z                       |
| 1               | 0.027     | -1.7                    | 1.196         | 0.6                     | 2.424        | 1.2                     | 0.018       | -0.9                    | 0.025           | -1.1                    | 0.058         | -1.2                    | 1.737         | 1.0                     | 1.008        | 1.2                     |
| 2               | 0.042     | -0.4                    | 1.405         | 1.4                     | 2.015        | 0.3                     | 0.022       | -0.2                    | 0.034           | 0.0                     | 0.07          | -0.7                    | 1.465         | 0.2                     | 1.007        | 1.2                     |
| 4               | 0.0498    | 0.3                     | 0.892         | -0.6                    | 1.61         | -0.6                    | 0.0269      | 0.6                     | 0.0375          | 0.4                     | 0.0928        | 0.4                     | 1.22          | -0.5                    | 0.71         | -0.3                    |
| 5               | 0.046     | 0.0                     | 0.706         | -1.3                    | 1.48         | -0.9                    | 0.025       | 0.3                     | 0.033           | -0.1                    | 0.074         | -0.5                    | 1.11          | -0.8                    | 0.564        | -1.1                    |
| 6               |           |                         | 0.814         | -0.9                    | 1.392        | -1.0                    | 0.019       | -0.7                    | 0.03            | -0.5                    | 0.067         | -0.8                    | 1.196         | -0.6                    | 0.531        | -1.3                    |
| 8               | 0.0379    | -0.7                    | 0.91          | -0.5                    | 1.75         | -0.3                    | 0.0235      | 0.0                     | 0.035           | 0.1                     | 0.078         | -0.3                    | 1.58          | 0.5                     | 0.73         | -0.2                    |
| 9               | 0.041     | -0.5                    | 0.972         | -0.3                    | 1.98         | 0.2                     | FN          | -3.1                    | 0.028           | -0.7                    | 0.068         | -0.8                    | 1.45          | 0.1                     | 0.8          | 0.1                     |
| 10              | 0.0353    | -0.9                    | 0.54          | -1.9                    | 0.865        | -2.2                    | 0.0204      | -0.5                    | 0.0336          | 0.0                     | 0.0822        | -0.1                    | 0.589         | -2.3                    | 0.305        | -2.4                    |
| 11              |           |                         | 0.78          | -1.0                    | 1.19         | -1.5                    | 0.0171      | -1.1                    | 0.0289          | -0.6                    | 0.116         | 1.5                     | 0.405         | -2.8                    |              |                         |
| 12              |           |                         | 1.13          | 0.3                     | 1.87         | 0.0                     | 0.023       | 0.0                     | 0.033           | -0.1                    |               |                         | 1.4           | 0.0                     | 0.802        | 0.1                     |
| 13              | 0.0468    | 0.0                     | 1.07          | 0.1                     | 1.98         | 0.2                     | 0.0257      | 0.4                     | 0.0396          | 0.7                     | 0.0941        | 0.5                     | 1.41          | 0.0                     | 0.775        | 0.0                     |
| 14              | 0.064     | 1.5                     | 1.076         | 0.1                     | 1.932        | 0.1                     | 0.027       | 0.6                     | 0.033           | -0.1                    | 0.103         | 0.9                     | 1.689         | 0.8                     | 0.894        | 0.6                     |
| 15              | 0.051     | 0.4                     | 1.089         | 0.2                     | 1.838        | -0.1                    | 0.028       | 0.8                     | 0.035           | 0.1                     | 0.088         | 0.2                     | 1.376         | -0.1                    | 0.717        | -0.3                    |
| 16              | 0.051     | 0.4                     | 1.098         | 0.2                     | 1.943        | 0.1                     | 0.031       | 1.3                     | 0.039           | 0.6                     | 0.112         | 1.3                     | 1.644         | 0.7                     | 0.912        | 0.7                     |
| 17              |           |                         |               |                         |              |                         |             |                         |                 |                         |               |                         |               |                         |              |                         |
| 18              |           |                         | 1.11          | 0.3                     | 1.74         | -0.3                    | 0.0182      | -0.9                    | 0.026           | -0.9                    |               |                         | 1.01          | -1.1                    | 0.889        | 0.6                     |
| 19              |           |                         | 0.173         | -3.3                    | 0.274        | -3.4                    | 0.012       | -1.9                    | 0.018           | -1.9                    | 0.036         | -2.3                    | 0.151         | -3.6                    |              |                         |
| 20              |           |                         |               |                         |              |                         | 0.037       | 2.4                     | 0.058           | 2.8                     |               |                         |               |                         |              |                         |
| 21              | 0.055     | 0.8                     | 0.809         | -0.9                    | 2.045        | 0.3                     | 0.023       | 0.0                     | 0.037           | 0.4                     | 0.086         | 0.1                     | 1.356         | -0.1                    | 0.839        | 0.3                     |
| 22              | 0.036     | -0.9                    | 1.146         | 0.4                     | 2.086        | 0.4                     | 0.028       | 8.0                     | 0.034           | 0.0                     | 0.093         | 0.4                     | 1.556         | 0.4                     | 0.9          | 0.7                     |
| 23              |           |                         | 1.06          | 0.1                     | 1.9          | 0.0                     | 0.023       | 0.0                     | 0.049           | 1.8                     | FN            | -3.5                    | 1.43          | 0.1                     | 0.698        | -0.4                    |
| 24              | 0.053     | 0.6                     | 1.31          | 1.0                     | 2.06         | 0.4                     | 0.02        | -0.6                    | 0.033           | -0.1                    | 0.091         | 0.3                     | 1.42          | 0.1                     | 0.781        | 0.0                     |
| 25              | 0.045     | 0.4                     | 4.40          | 0.0                     | 4 50         | 0.0                     | 0.031       | 1.3                     | 0.042           | 0.9                     | 0.0046        | 0.0                     | 4.00          | 0.4                     | 0.60#        | 0.7                     |
| 26              | 0.0476    | 0.1                     | 1.12          | 0.3                     | 1.73         | -0.3                    | 0.0233      | 0.0                     | 0.0347          | 0.1                     | 0.0846        | 0.0                     | 1.38          | -0.1                    | 0.635        | -0.7                    |
| 27              | 0.04      | -0.5                    | 0.92          | -0.5                    | 1.88         | 0.0                     | 0.022       | -0.2                    | 0.034           | 0.0                     | 0.067         | -0.8                    | 1.429         | 0.1                     | 0.633        | -0.7                    |
| 28              | 0.055     | 0.8                     | 1.33<br>0.895 | -0.6                    | 2.1<br>1.796 | -0.2                    | FN<br>0.033 | -3.1<br>1.7             | 0.039           | 0.6<br>2.1              | 0.094         | 0.5                     | 1.34<br>1.139 | -0.2<br>-0.7            | 0.81         | -0.9                    |
| 30              | 0.033     | -0.3                    | 1.02          | -0.0                    | 1.73         | -0.2                    | 0.0293      | 1.0                     | 0.032           | -0.5                    | 0.0913        | 0.3                     | 1.34          | -0.7                    | 0.678        | -0.5                    |
| 31              | 0.0420    | -0.2                    | 0.893         | -0.1                    | 1.49         | -0.8                    | 0.02        | -0.6                    | 0.023           | -1.3                    | 0.06          | -1.2                    | 1.69          | 0.8                     | 0.969        | 1.0                     |
| 32              | 0.011     | 0.2                     | 0.075         | 0.0                     | 1.17         | 0.0                     | 0.0183      | -0.9                    | 0.03            | -0.5                    | 0.00          | 1.2                     | 1.07          | 0.0                     | 0.505        | 1.0                     |
| 33              |           |                         | 1.11          | 0.3                     | 1.86         | 0.0                     | 0.021       | -0.4                    | 0.026           | -0.9                    |               |                         | 1.45          | 0.1                     | 0.85         | 0.4                     |
| 34              | 0.042     | -0.4                    | 0.768         | -1.1                    | 1.283        | -1.3                    | 0.024       | 0.1                     | 0.031           | -0.4                    | 0.075         | -0.4                    | 1.166         | -0.7                    | 0.573        | -1.0                    |
| 35              | 0.035     | -1.0                    | 1.209         | 0.6                     | 1.793        | -0.2                    | 0.01        | -2.3                    | 0.031           | -0.4                    | 0.093         | 0.4                     | 1.77          | 1.1                     | 0.926        | 0.8                     |
| 36              | 0.039     | -0.6                    | 0.711         | -1.3                    | 1.354        | -1.1                    | 0.024       | 0.1                     | 0.034           | 0.0                     | 0.08          | -0.2                    | 1.157         | -0.7                    | 0.662        | -0.6                    |
| 37              | 0.036     | -0.9                    | 0.82          | -0.9                    | 1.4          | -1.0                    | 0.028       | 0.8                     | 0.034           | 0.0                     | 0.091         | 0.3                     | 0.94          | -1.3                    | 0.51         | -1.4                    |
| 38              |           |                         | 0.981         | -0.2                    | 1.8          | -0.2                    |             |                         |                 |                         | 0.072         | -0.6                    |               |                         |              |                         |
| 39              |           |                         | 1             | -0.2                    | 1.8          | -0.2                    | 0.025       | 0.3                     | 0.038           | 0.5                     | 0.103         | 0.9                     | 1.59          | 0.5                     |              |                         |
| 40              |           |                         |               |                         |              |                         | 0.0282      | 0.8                     | 0.0444          | 1.2                     |               |                         |               |                         |              |                         |
|                 |           |                         |               |                         |              |                         |             |                         |                 |                         |               |                         |               |                         |              |                         |

| Taboratory code | Aclonifen | Z-scores (FFP RSD (25%) | Azoxystrobin | Z-scores (FFP RSD (25%) | Boscalid      | Z-scores (FFP RSD (25%) | Dieldrin | Z-scores (FFP RSD (25%) | Endosulfan-beta | Z-scores (FFP RSD (25%) | Fenpropathrin | Z-scores (FFP RSD (25%) | Fluopyram     | Z-scores (FFP RSD (25%) | Fluxapyroxad   | Z-scores (FFP RSD (25%) |
|-----------------|-----------|-------------------------|--------------|-------------------------|---------------|-------------------------|----------|-------------------------|-----------------|-------------------------|---------------|-------------------------|---------------|-------------------------|----------------|-------------------------|
| Assigned        |           | -score                  |              | -score                  |               | -score                  |          | -score                  |                 | -score                  |               | -score                  |               | -score                  |                | -score                  |
| value           | 0.046     | Z                       | 1.04         | Z                       | 1.88          | Z                       | 0.023    | Z                       | 0.034           | 2                       | 0.084         | Z                       | 1.40          | Z                       | 0.774          | 2                       |
| 41              |           |                         | 1.273        | 0.9                     | 2.352         | 1.0                     | FN       | -3.1                    | FN              | -3.5                    |               |                         |               |                         |                |                         |
| 42              |           |                         | 1.39         | 1.3                     | 2.1           | 0.5                     | 0.0335   | 1.8                     | 0.0511          | 2.0                     |               |                         | 1.5           | 0.3                     | 0.893          | 0.6                     |
| 43              | 0.05      | 0.3                     | 1.46         | 1.6                     | 2.48          | 1.3                     | 0.027    | 0.6                     | 0.037           | 0.4                     | 0.084         | 0.0                     | 1.65          | 0.7                     | 0.92           | 0.8                     |
| 44<br>45        | 0.022     | -2.1<br>1.6             | 0.966        | -0.3                    | 2.11          | -0.2                    | FN 0.02  | -3.1                    | FN<br>0.05      | -3.5<br>1.9             | 0.102         | 0.8                     | 1.05          | -1.7<br>-1.0            | 0.405<br>0.845 | -1.9                    |
| 46              | 0.003     | 1.0                     | 1            | -0.2                    | 1.8           | -0.2                    | 0.0134   | -1.7                    | 0.0207          | -1.6                    | 0.0823        | -0.1                    | 1.39          | 0.0                     | 0.863          | 0.5                     |
| 47              |           |                         | 0.66         | -1.5                    | -10           | · · ·                   | FN       | -3.1                    | FN              | -3.5                    | 0.039         | -2.1                    |               |                         |                | 0.0                     |
| 48              |           |                         |              |                         |               |                         | 0.0281   | 0.8                     | 0.0437          | 1.1                     |               |                         |               |                         |                |                         |
| 49              |           |                         | 0.765        | -1.1                    | 1.714         | -0.4                    | FN       | -3.1                    | 0.032           | -0.2                    | 0.086         | 0.1                     | 1.17          | -0.7                    |                |                         |
| 50              | 0.0483    | 0.2                     | 1.04         | 0.0                     | 2.05          | 0.4                     | 0.023    | 0.0                     | 0.0358          | 0.2                     | 0.089         | 0.2                     | 1.36          | -0.1                    | 0.74           | -0.2                    |
| 51              | 0.054     | 0.7                     | 1.05         | 0.0                     | 1.94          | 0.1                     | 0.041    | 3.0                     | 0.036           | 0.2                     | 0.089         | 0.2                     | 1.4           | 0.0                     | 1.023          | 1.3                     |
| 52              | FN        | -3.1                    | 0.699        | -1.3                    | 1.353         | -1.1                    | FN       | -3.1                    | 0.034           | 0.0                     | 0.086         | 0.1                     | 0.972         | -1.2                    | 0.518          | -1.3                    |
| 53              | 0.044     | -0.2                    | 0.306        | -2.8                    | 0.616         | -2.7                    | FN       | -3.1                    | 0.059           | 2.9                     | 0.088         | 0.2                     | 0.188         | -3.5                    | 0.126          | -3.3                    |
| 54<br>55        | 0.029     | -1.5                    | 1.04         | 0.0                     | 1.94<br>2.551 | 0.1                     | 0.018    | -0.9                    | 0.028           | -0.7                    | 0.08          | -0.2                    | 1.53<br>1.866 | 1.3                     | 0.714          | -0.3                    |
| 56              | 0.029     | -0.6                    | 1.431        | 1.5                     | 2.561         | 1.4                     | 0.018    | 0.0                     | 0.028           | -0.7                    | 0.108         | 1.1                     | 1.822         | 1.2                     | 1.008          | 1.2                     |
| 57              |           | 0.0                     |              | -10                     |               |                         | 0.0206   | -0.5                    | 0.0316          | -0.3                    |               |                         |               |                         |                |                         |
| 58              | 0.036     | -0.9                    | 1.326        | 1.1                     | 2.661         | 1.7                     | 0.018    | -0.9                    | 0.029           | -0.6                    | 0.066         | -0.9                    | 1.682         | 0.8                     | 0.863          | 0.5                     |
| 59              |           |                         |              |                         |               |                         | 0.021    | -0.4                    | 0.022           | -1.4                    |               |                         |               |                         |                |                         |
| 60              |           |                         |              |                         |               |                         | 0.0232   | 0.0                     | 0.0336          | 0.0                     |               |                         |               |                         |                |                         |
| 61              | FN        | -3.1                    | 0.595        | -1.7                    | 1.52          | -0.8                    | 0.014    | -1.6                    | 0.025           | -1.1                    | 0.089         | 0.2                     | 1.21          | -0.5                    | 0.714          | -0.3                    |
| 62              | 0.049     | 0.2                     | 1.07         | 0.1                     | 1.58          | -0.6                    | 0.015    | -1.4                    | 0.046           | 1.4                     | 0.083         | -0.1                    | 1.17          | -0.7                    | 0.729          | -0.2                    |
| 63              | 0.051     | 0.4                     | 1.5          | 1.8                     | 2.54          | 1.4                     | 0.027    | 0.6                     | 0.043           | 1.1                     | 0.105         | 1.0                     | 1.83          | 1.2                     | 1.15           | 1.9                     |
| 64<br>65        | 0.046     | 0.0                     | 1.258        | 0.8                     | 2.54          | 1.4                     | 0.024    | 0.1                     | 0.033           | -0.1<br>1.1             | 0.08          | -0.2                    | 1.74          | 1.0                     | 0.956          | 0.9                     |
| 66              | 0.0385    | -0.7                    | 0.638        | -1.6                    | 1.26          | -1.3                    | 0.03     | -0.6                    | 0.043           | -0.2                    | 0.0803        | -0.2                    | 1.04          | -1.0                    | 0.61           | -0.8                    |
| 67              |           |                         | 1.36         | 1.2                     | 4.74          | >5                      | 0.0254   | 0.4                     | 0.0351          | 0.1                     | 0.0792        | -0.2                    | 3.45          | >5                      |                |                         |
| 68              | 0.042     | -0.4                    | 0.97         | -0.3                    | 3.8           | 4.1                     | 0.024    | 0.1                     | 0.036           | 0.2                     |               |                         |               |                         |                |                         |
| 69              | 0.039     | -0.6                    | 1.43         | 1.5                     | 2.3           | 0.9                     | FN       | -3.1                    | FN              | -3.5                    | 0.093         | 0.4                     | 1.47          | 0.2                     | 0.905          | 0.7                     |
| 70              | 0.0392    | -0.6                    | 1.13         | 0.3                     | 1.92          | 0.1                     | 0.0208   | -0.4                    | 0.0344          | 0.0                     | 0.0678        | -0.8                    | 0.8           | -1.7                    | 0.71           | -0.3                    |
| 71              | 0.04      | -0.5                    | 1.3          | 1.0                     | 2.1           | 0.5                     | FN       | -3.1                    | 0.034           | 0.0                     | 0.091         | 0.3                     | 1.2           | -0.6                    | 0.55           | -1.2                    |
| 72              | 0.053     | 0.6                     | 1.67         | 2.4                     | 2.75          | 1.8                     | 0.021    | -0.4                    | FN              | -3.5                    | 0.073         | -0.5                    | 1.82          | 1.2                     | 0.998          | 1.2                     |
| 73<br>74        |           |                         | 1.32         | 0.3                     | 1.91<br>1.87  | 0.1                     | 0.053    | >5<br>-0.1              | 0.033           | -0.1                    | FN            | -3.5                    | 1.28          | -0.3<br>0.1             | 0.63           | -0.7<br>-0.1            |
| 75              | 0.049     | 0.2                     | 1.12         | 1.1                     | 2.424         | 1.2                     | 0.0226   | 0.3                     | 0.035           | 0.0                     | 0.071         | -0.6                    | 1.705         | 0.1                     | 0.749          | 0.8                     |
| 76              | 0.076     | 2.6                     | 1.499        | 1.8                     | 2.422         | 1.2                     | 0.020    | 0.0                     | 0.000           | V.1                     | 0.11          | 1.2                     | 1.745         | 1.0                     | 0.879          | 0.5                     |
| 77              |           |                         |              |                         |               |                         |          |                         |                 |                         |               |                         |               |                         |                |                         |
| 78              | 0.0542    | 0.7                     | 1.27         | 0.9                     | 2.63          | 1.6                     | 0.0146   | -1.5                    | 0.0241          | -1.2                    | 0.103         | 0.9                     | 2.29          | 2.5                     | 1.39           | 3.2                     |
| 79              | 0.043     | -0.3                    | 0.862        | -0.7                    | 2.283         | 0.9                     | 0.022    | -0.2                    | 0.028           | -0.7                    | 0.076         | -0.4                    | 1.501         | 0.3                     | 0.897          | 0.6                     |
| 80              |           |                         | 1.01         | -0.1                    | 1.77          | -0.2                    | 0.0199   | -0.6                    | 0.0272          | -0.8                    |               |                         | 1.2           | -0.6                    | 0.766          | 0.0                     |

| ope             | _         | <u>.</u>                | ii           | <u>.</u>                |          | ි                       |          | 9                       | oeta            | 9                       | rin           | ල                       | e e       | ္က                      | ad           | ြွ                      |
|-----------------|-----------|-------------------------|--------------|-------------------------|----------|-------------------------|----------|-------------------------|-----------------|-------------------------|---------------|-------------------------|-----------|-------------------------|--------------|-------------------------|
| Laboratory code | Aclonifen | .25%                    | Azoxystrobin | .25%                    | Boscalid | .25%                    | Dieldrin | .25%                    | Endosulfan-beta | .25%                    | Fenpropathrin | .25%                    | Fluopyram | .25%                    | Fluxapyroxad | .25%                    |
| orato           | Aclo      | SD (                    | oxys         | SD (                    | Bosc     | SD (                    | Diel     | SD (                    | osal            | SD (                    | ıpro          | SD (                    | luop      | SD (                    | тар          | SD (                    |
| Lab             |           | FP.                     | Az           | FP F                    |          | FP F                    |          | FP F                    | End             | FP.                     | Fer           | FP F                    | E.        | FP F                    | 띮            | HP H                    |
| MRRL            | 0.01      | Z-scores (FFP RSD (25%) | 0.01         | Z-scores (FFP RSD (25%) | 0.01     | Z-scores (FFP RSD (25%) | 0.01     | Z-scores (FFP RSD (25%) | 0.01            | Z-scores (FFP RSD (25%) | 0.01          | Z-scores (FFP RSD (25%) | 0.01      | Z-scores (FFP RSD (25%) | 0.01         | Z-scores (FFP RSD (25%) |
| Assigned        | 0.01      | scor                    | 0.01         | scor                    | 0.01     | scor                    | 0.01     | scor                    | 0.01            | scor                    | 0.01          | scor                    | 0.01      | scor                    | 0.01         | scor                    |
| value           | 0.046     | Z                       | 1.04         | Ż                       | 1.88     | Ż                       | 0.023    | Ż                       | 0.034           | Z                       | 0.084         | ż                       | 1.40      | Ż                       | 0.774        | Ż                       |
| 81              | 0.053     | 0.6                     | 1.272        | 0.9                     | 2.174    | 0.6                     | 0.025    | 0.3                     | 0.039           | 0.6                     | 0.089         | 0.2                     | 2.549     | 3.3                     | 1.607        | 4.3                     |
| 82              | 0.044     | -0.2                    | 1.12         | 0.3                     | 1.5      | -0.8                    | 0.0225   | -0.1                    | 0.0261          | -0.9                    | 0.0737        | -0.5                    | 1.39      | 0.0                     | 0.811        | 0.2                     |
| 83              | 0.011     | 0.2                     | 1.184        | 0.5                     | 2.71     | 1.8                     | 0.0216   | -0.3                    | 0.0288          | -0.6                    | 0.0613        | -1.1                    | 1.786     | 1.1                     | 0.011        | 0.2                     |
| 84              | 0.047     | 0.1                     | 0.671        | -1.4                    | 1.002    | -1.9                    | 0.0210   | 0.0                     | 0.034           | 0.0                     | 0.102         | 0.8                     | 0.7       | -2.0                    | 0.57         | -1.1                    |
| 85              | 0.047     | 0.1                     | 1.22         | 0.7                     | 1.76     | -0.3                    | 0.023    | -1.3                    | 0.034           | -0.7                    | 0.102         | 0.0                     | 1.43      | 0.1                     | 0.702        | -0.4                    |
| 86              | 0.042     | -0.4                    | 0.722        | -1.2                    |          | -0.3                    | 0.010    | -1.5                    | 0.020           | -0.7                    | 0.082         | -0.1                    | 1.008     |                         |              |                         |
|                 | 0.042     | -0.4                    |              |                         | 1.729    |                         | 0.016    | 4.0                     | 0.022           | 1.4                     |               |                         |           | -1.1                    | 0.565        | -1.1                    |
| 87              | 0.0505    | 0.4                     | 0.83         | -0.8                    | 1.3      | -1.2                    | 0.016    | -1.3                    | 0.022           | -1.4                    | 0.086         | 0.1                     | 5.4       | >5                      | 2.4          | >5                      |
| 88              | 0.0507    | 0.4                     | 1.03         | 0.0                     | 1.85     | -0.1                    | 0.0199   | -0.6                    | 0.0319          | -0.2                    | 0.084         | 0.0                     | 1.26      | -0.4                    | 0.808        | 0.2                     |
| 89              | 0.062     | 1.4                     | 0.796        | -0.9                    | 1.48     | -0.9                    | 0.02     | -0.6                    | 0.033           | -0.1                    | 0.072         | -0.6                    | 1.167     | -0.7                    | 0.594        | -0.9                    |
| 90              | 0.048     | 0.2                     | 0.91         | -0.5                    | 2        | 0.3                     | 0.024    | 0.1                     | 0.03            | -0.5                    | 0.081         | -0.2                    | 1.24      | -0.5                    | 0.74         | -0.2                    |
| 91              | FN        | -3.1                    | 1.08         | 0.1                     | 2.28     | 0.8                     | 0.018    | -0.9                    | 0.022           | -1.4                    | 0.062         | -1.1                    | 1.45      | 0.1                     | 0.87         | 0.5                     |
| 92              |           |                         | 1.198        | 0.6                     | 1.82     | -0.1                    | 0.019    | -0.7                    | 0.029           | -0.6                    |               |                         | 1.426     | 0.1                     | 0.824        | 0.3                     |
| 93              | 0.04      | -0.5                    | 0.937        | -0.4                    | 1.686    | -0.4                    | FN       | -3.1                    | 0.037           | 0.4                     | 0.05          | -1.6                    | 1.399     | 0.0                     | 0.746        | -0.1                    |
| 94              | 0.0465    | 0.0                     | 0.996        | -0.2                    | 1.61     | -0.6                    | 0.0235   | 0.0                     | 0.0368          | 0.3                     | 0.111         | 1.3                     | 1.1       | -0.9                    | 0.631        | -0.7                    |
| 95              | 0.054     | 0.7                     | 1.112        | 0.3                     | 2.142    | 0.6                     | 0.02     | -0.6                    | 0.029           | -0.6                    | 0.088         | 0.2                     | 2.06      | 1.9                     | 0.991        | 1.1                     |
| 96              |           |                         | 0.95         | -0.4                    | 1.8      | -0.2                    | 0.023    | 0.0                     | 0.034           | 0.0                     | 0.092         | 0.4                     |           |                         |              |                         |
| 97              | 0.0327    | -1.2                    | 0.598        | -1.7                    | 1.05     | -1.8                    | 0.018    | -0.9                    | 0.0267          | -0.9                    | 0.0691        | -0.7                    | 0.672     | -2.1                    | 0.357        | -2.2                    |
| 98              | 0.0419    | -0.4                    | 0.547        | -1.9                    | 1.93     | 0.1                     | 0.0209   | -0.4                    | 0.0286          | -0.6                    | 0.0749        | -0.4                    | 1.46      | 0.2                     | 0.732        | -0.2                    |
| 99              | FN        | -3.1                    | 1.264        | 0.9                     | 2.186    | 0.6                     | 0.025    | 0.3                     | 0.041           | 0.8                     | 0.079         | -0.2                    | 1.272     | -0.4                    | 1.026        | 1.3                     |
| 100             | 0.044     | -0.2                    | 1.08         | 0.1                     | 1.69     | -0.4                    | 0.0218   | -0.3                    | 0.0366          | 0.3                     | 0.0855        | 0.1                     | 1.36      | -0.1                    | 0.702        | -0.4                    |
| 101             |           |                         | 1.25         | 0.8                     | 2.33     | 1.0                     | 0.0313   | 1.4                     | 0.0399          | 0.7                     | 0.0945        | 0.5                     | 1.68      | 0.8                     | 0.839        | 0.3                     |
| 102             |           |                         | 0.656        | -1.5                    | 1.152    | -1.5                    | 0.028    | 0.8                     | 0.04            | 0.7                     | 0.091         | 0.3                     | 0.869     | -1.5                    | 0.427        | -1.8                    |
| 103             | 0.021     | -2.2                    | 1.149        | 0.4                     | 1.871    | 0.0                     | 0.027    | 0.6                     | 0.039           | 0.6                     | 0.075         | -0.4                    | 1.345     | -0.2                    | 0.679        | -0.5                    |
| 105             | FN        | -3.1                    | 0.93         | -0.4                    | 1.97     | 0.2                     | FN       | -3.1                    | FN              | -3.5                    | FN            | -3.5                    | 1.49      | 0.3                     | 0.75         | -0.1                    |
| 106             | 0.047     | 0.1                     | 1.3          | 1.0                     | 1.5      | -0.8                    | 0.02     | -0.6                    | 0.019           | -1.8                    | 0.035         | -2.3                    | 1.4       | 0.0                     | 0.9          | 0.7                     |
| 107             |           |                         | 1.27         | 0.9                     | 2        | 0.3                     | 0.043    | 3.4                     | 0.068           | 4.0                     | 0.081         | -0.2                    | 1.55      | 0.4                     | 0.709        | -0.3                    |
| 108             |           |                         | 1.407        | 1.4                     | 2.604    | 1.5                     |          |                         | 0.038           | 0.5                     | 0.084         | 0.0                     | 1.864     | 1.3                     | 1.146        | 1.9                     |
| 109             | FN        | -3.1                    | 1.32         | 1.1                     | 2.11     | 0.5                     | 0.0273   | 0.7                     | 0.0355          | 0.2                     | 0.0545        | -1.4                    | 1.57      | 0.5                     | 0.859        | 0.4                     |
| 110             | 0.0539    | 0.7                     | 0.893        | -0.6                    | 1.82     | -0.1                    | 0.0193   | -0.7                    | 0.0281          | -0.7                    | 0.0658        | -0.9                    | 1.04      | -1.0                    | 0.713        | -0.3                    |
| 111             |           |                         | 1.287        | 0.9                     | 2.301    | 0.9                     | 0.021    | -0.4                    | 0.037           | 0.4                     | 0.076         | -0.4                    | 1.549     | 0.4                     | 0.893        | 0.6                     |
| 112             |           |                         | 0.626        | -1.6                    | 0.949    | -2.0                    | 0.017    | -1.1                    | 0.028           | -0.7                    | 0.075         | -0.4                    | 3.13      | 4.9                     | 0.909        | 0.7                     |
| 113             | 0.054     | 0.7                     | 0.879        | -0.6                    | 1.68     | -0.4                    | 0.024    | 0.1                     | 0.039           | 0.6                     | 0.094         | 0.5                     | 1.196     | -0.6                    | 0.644        | -0.7                    |
| 114             | 0.0436    | -0.2                    | 1.03         | 0.0                     | 2.39     | 1.1                     | 0.0293   | 1.0                     | 0.0377          | 0.4                     | 0.102         | 0.8                     | 1.63      | 0.7                     | 0.892        | 0.6                     |
| 115             | 0.048     | 0.2                     | 0.89         | -0.6                    | 1.59     | -0.6                    | 0.021    | -0.4                    | 0.034           | 0.0                     | 0.074         | -0.5                    | 1.16      | -0.7                    | 0.61         | -0.8                    |
| 116             | 0.0506    | 0.4                     | 0.722        | -1.2                    | 1.43     | -1.0                    | 0.0252   | 0.3                     | 0.0341          | 0.0                     | 0.0842        | 0.0                     | 0.811     | -1.7                    | 0.402        | -1.9                    |
| 117             | 0.0614    | 1.3                     | 1.32         | 1.1                     | 2.16     | 0.6                     | 0.0232   | -0.1                    | 0.0353          | 0.2                     | 0.107         | 1.1                     | 1.58      | 0.5                     | 0.818        | 0.2                     |
| 118             | 3.0011    | 2.0                     | 0.989        | -0.2                    | 1.852    | -0.1                    | 5.0225   | J.1                     | 3.0333          | U.L                     | 0.0837        | 0.0                     | 1.33      | -0.2                    | 0.752        | -0.1                    |
| 119             | 0.091     | 3.9                     | 1.398        | 1.4                     | 2.431    | 1.2                     | 0.028    | 0.8                     | 0.043           | 1.1                     | 0.112         | 1.3                     | 1.774     | 1.1                     | 0.752        | 0.9                     |
| 120             | 0.091     | -0.2                    | 0.91         | -0.5                    | 1.62     |                         | 0.020    | 0.0                     | 0.043           | 1.1                     | 0.112         | 1.7                     | 1.774     | -0.4                    | 0.954        | -0.6                    |
|                 |           |                         |              |                         |          | -0.6                    | 0.024    | 0.1                     |                 |                         |               |                         |           |                         |              |                         |
| 121             | 0.04      | -0.5                    | 0.73         | -1.2                    | 1.7      | -0.4                    | 0.024    | 0.1                     | 0.036           | 0.2                     | 0.073         | -0.5                    | 1.24      | -0.5                    | 0.76         | -0.1                    |

| Laboratory code | Aclonifen        | Z-scores (FFP RSD (25%) | Azoxystrobin | Z-scores (FFP RSD (25%) | Boscalid     | Z-scores (FFP RSD (25%) | Dieldrin | Z-scores (FFP RSD (25%) | Endosulfan-beta | Z-scores (FFP RSD (25%) | Fenpropathrin    | Z-scores (FFP RSD (25%) | Fluopyram | Z-scores (FFP RSD (25%) | Fluxapyroxad | Z-scores (FFP RSD (25%) |
|-----------------|------------------|-------------------------|--------------|-------------------------|--------------|-------------------------|----------|-------------------------|-----------------|-------------------------|------------------|-------------------------|-----------|-------------------------|--------------|-------------------------|
| MRRL            | 0.01             | core                    | 0.01         | core                    | 0.01         | core                    | 0.01     | core                    | 0.01            | core                    | 0.01             | core                    | 0.01      | core                    | 0.01         | core                    |
| Assigned        | 0.046            | )S-Z                    | 1.04         | Z-S(                    | 1.88         | )S-Z                    | 0.023    | Z-S(                    | 0.034           | )S-Z                    | 0.084            | S-Z                     | 1.40      | Z-s(                    | 0.774        | Z-S(                    |
| value           |                  |                         |              |                         |              |                         |          |                         |                 |                         |                  |                         |           |                         |              |                         |
| 122             | 0.0412           | -0.4                    | 0.76         | -1.1                    | 1.56         | -0.7                    | 0.0473   | 4.1                     | 0.0329          | -0.1                    | 0.096            | 0.6                     | 1.27      | -0.4                    | 0.616        | -0.8                    |
| 123             | 0.0734           | 2.3                     | 1.44         | 1.5                     | 2.03         | 0.3                     | 0.0392   | 2.7                     | 0.0466          | 1.5                     | 0.125            | 1.9                     | 1.69      | 0.8                     | 0.666        | -0.6                    |
| 124             |                  |                         | 1.2          | 0.6                     | 1.09         | -1.7                    | 0.0295   | 1.1                     | 0.0278          | -0.7                    | 0.0925           | 0.4                     | 1.32      | -0.2                    | 0.59         | -1.0                    |
| 125             |                  |                         | 0.925        | -0.4                    |              |                         | 0.021    | -0.4                    | 0.022           | -1.4                    | 0.06             | -1.2                    | 0.884     | -1.5                    |              |                         |
| 126             | 0.116            | >5                      | 0.904        | -0.5                    | 1.756        | -0.3                    | 0.032    | 1.5                     | 0.06            | 3.1                     | 0.166            | 3.9                     | 1.552     | 0.4                     | 0.911        | 0.7                     |
| 127             |                  |                         |              |                         |              |                         |          |                         | 0.044           | 1.2                     | 0.124            | 1.9                     | 1.22      | -0.5                    |              |                         |
| 128             |                  |                         | 0.655        | -1.5                    | 0.373        | -3.2                    | 0.0287   | 0.9                     | 0.0181          | -1.9                    |                  |                         |           |                         | FN           | -3.9                    |
| 129             | 0.049            | 0.2                     | 0.497        | -2.1                    | 0.835        | -2.2                    | 0.0262   | 0.5                     | 0.0358          | 0.2                     | 0.12             | 1.7                     | 0.569     | -2.4                    | 0.306        | -2.4                    |
| 130             |                  |                         | 1.13         | 0.3                     | 2.29         | 0.9                     | 0.022    | -0.2                    | 0.032           | -0.2                    | 0.086            | 0.1                     | 1.63      | 0.7                     | 0.87         | 0.5                     |
| 131             |                  |                         | 1.54         | 1.9                     | 1.77         | -0.2                    | 0.014    | -1.6                    | 0.026           | -0.9                    | 0.000            | 0.1                     | 1.65      | 0.7                     | 0.68         | -0.5                    |
| 132             | 0.049            | 0.2                     | 1.54         | -0.2                    | 1.8          | -0.2                    | 0.022    | -0.2                    | 0.032           | -0.2                    | 0.092            | 0.4                     | 1.4       | 0.0                     | 0.72         | -0.3                    |
|                 |                  |                         |              |                         |              |                         |          |                         |                 |                         |                  |                         |           |                         |              |                         |
| 133             | 0.049            | 0.2                     | 1.25         | 0.8                     | 1.96         | 0.2                     | 0.025    | 0.3                     | 0.036           | 0.2                     | 0.1              | 0.7                     | 1.44      | 0.1                     | 0.82         | 0.2                     |
| 134             | 0.054            | 0.7                     | 1.018        | -0.1                    | 1.503        | -0.8                    | 0.023    | 0.0                     | 0.034           | 0.0                     | 0.094            | 0.5                     | 1.197     | -0.6                    | 0.549        | -1.2                    |
| 135             |                  |                         | 0.924        | -0.5                    | 1.784        | -0.2                    | 0.023    | 0.0                     | 0.036           | 0.2                     | 0.074            | -0.5                    | 1.219     | -0.5                    |              |                         |
| 136             |                  |                         | 1.5          | 1.8                     |              |                         | 0.024    | 0.1                     | 0.038           | 0.5                     | 0.1              | 0.7                     |           |                         |              |                         |
| 137             |                  |                         | 1.07         | 0.1                     | 1.56         | -0.7                    | 0.0175   | -1.0                    | 0.0311          | -0.3                    | 0.0891           | 0.2                     | 1.43      | 0.1                     | 0.635        | -0.7                    |
| 138             | 0.0501           | 0.3                     | 1.01         | -0.1                    | 1.94         | 0.1                     | 0.0208   | -0.4                    | 0.0391          | 0.6                     | FN               | -3.5                    | 1.61      | 0.6                     | 0.795        | 0.1                     |
| 139             | 0.0529           | 0.6                     | 1.14         | 0.4                     | 1.96         | 0.2                     | 0.0225   | -0.1                    | 0.0321          | -0.2                    | 0.0839           | 0.0                     | 1.27      | -0.4                    | 0.645        | -0.7                    |
| 140             | 0.0555           | 0.8                     | 1.3          | 1.0                     | 2.57         | 1.5                     | 0.0265   | 0.6                     | 0.0385          | 0.5                     | 0.0915           | 0.3                     | 1.97      | 1.6                     | 0.854        | 0.4                     |
| 141             | FN               | -3.1                    | 0.83         | -0.8                    |              |                         | 0.02     | -0.6                    | 0.03            | -0.5                    | 0.067            | -0.8                    | 1.44      | 0.1                     | 0.71         | -0.3                    |
| 142             | FN               | -3.1                    | 0.789        | -1.0                    | 2.698        | 1.7                     | 0.021    | -0.4                    | 0.031           | -0.4                    | 0.072            | -0.6                    | 1.402     | 0.0                     | 0.833        | 0.3                     |
| 143             | 0.04             | -0.5                    | 0.88         | -0.6                    | 1.7          | -0.4                    | 0.028    | 0.8                     | 0.037           | 0.4                     | 0.091            | 0.3                     | 1.3       | -0.3                    | 0.72         | -0.3                    |
| 144             |                  |                         | 1.5          | 1.8                     | 2.67         | 1.7                     | 0.04     | 2.9                     | 0.051           | 2.0                     |                  |                         | 1.72      | 0.9                     | 0.85         | 0.4                     |
| 145             |                  |                         | 0.941        | -0.4                    | 1.73         | -0.3                    | 0.019    | -0.7                    | 0.017           | -2.0                    | 0.081            | -0.2                    |           |                         |              |                         |
| 146             | 0.04             | -0.5                    | 0.97         | -0.3                    | 1.8          | -0.2                    | 0.034    | 1.8                     | 0.04            | 0.7                     | 0.092            | 0.4                     | 1.22      | -0.5                    | 0.67         | -0.5                    |
| 147             | 0.045            | -0.1                    | 1.05         | 0.0                     | 2.05         | 0.4                     | 0.024    | 0.1                     | 0.034           | 0.0                     | 0.095            | 0.5                     | 1.82      | 1.2                     | 0.993        | 1.1                     |
| 148             | 0.0526           | 0.6                     | 0.78         | -1.0                    | 1.63         | -0.5                    | 0.0256   | 0.4                     | 0.04            | 0.7                     | 0.0816           | -0.1                    | 1.37      | -0.1                    | 0.778        | 0.0                     |
| 149             |                  |                         | 1.18         | 0.5                     | 2.09         | 0.4                     | 0.019    | -0.7                    | 0.029           | -0.6                    | 0.088            | 0.2                     | 1.44      | 0.1                     | 1.07         | 1.5                     |
| 150             |                  |                         | 1.283        | 0.9                     | 2.359        | 1.0                     |          |                         | 0.038           | 0.5                     | 0.091            | 0.3                     | 1.699     | 0.8                     | 0.912        | 0.7                     |
| 151             | 0.055            | 0.8                     | 0.96         | -0.3                    | 1.989        | 0.2                     | 0.022    | -0.2                    | 0.041           | 0.8                     | 0.065            | -0.9                    | 1.486     | 0.2                     | 0.741        | -0.2                    |
| 152             | m.               | 0.4                     | 4.00         | 0.4                     | 0.6          |                         | 0.0040   |                         | 0.000           | 0.6                     | 0.0550           |                         | 4.05      | 4.0                     | 0.064        | 4.0                     |
| 153             | FN               | -3.1                    | 1.02         | -0.1                    | 2.6          | 1.5                     | 0.0213   | -0.3                    | 0.0293          | -0.6                    | 0.0553           | -1.4                    | 1.85      | 1.3                     | 0.961        | 1.0                     |
| 154             | 0.023            | -2.0                    | 0.383        | -2.5                    | 0.802        | -2.3                    | 0.02     | -0.6                    | 0.021           | -1.5                    | FN               | -3.5                    | 0.58      | -2.3                    | 0.367        | -2.1                    |
| 155             | 0.04             | -0.5                    | 1.79         | 2.9                     | 1.96         | 0.2                     | 0.03     | 1.2                     | 0.04            | 0.7                     | 0.08             | -0.2                    | 1.31      | -0.3                    | 1.23         | 2.4                     |
| 156             | 0.0506           | 0.4                     | 0.96         | -0.3<br>0.7             | 1.76         | -0.3                    | 0.0216   | 0.2                     | 0.0222          | -0.1                    | 0.0004           | 0.7                     | 1.6       | 0.6                     | 0.69         | -0.4                    |
| 157             | 0.0506<br>0.0404 | -0.5                    | 1.23         | 0.7                     | 2.28<br>1.78 | -0.2                    | 0.0216   | -0.3<br>0.2             | 0.0332          | 0.2                     | 0.0984<br>0.0812 | -0.1                    | 1.6       | 0.6                     | 0.878        | -0.3                    |
| 158<br>160      | 0.0404           | -0.5                    | 1.07         | 0.1                     | 1./8         | -0.2                    | 0.0247   | 0.2                     | 0.036           | 0.2                     | 0.0012           | -0.1                    | 1.43      | 0.1                     | 0.72         | -0.3                    |
| 161             | 0.0306           | -1.4                    | 1.14         | 0.4                     | 0.747        | -2.4                    | 0.024    | 0.1                     | 0.0348          | 0.1                     | 0.0846           | 0.0                     | 1.404     | 0.0                     | 0.81         | 0.2                     |
| 162             | 0.0300           | -1.4                    | 1.14         | 3.6                     | 2.55         | 1.4                     | FN       | -3.1                    | FN              | -3.5                    | 0.03             | -2.6                    | 1.67      | 0.8                     | 0.98         | 1.1                     |
| 163             | 0.028            | 0.1                     | 1.31         | 1.0                     | 2.13         | 0.5                     | 0.029    | 1.0                     |                 | 0.9                     | 0.03             | 0.6                     | 1.62      |                         | 0.98         | 0.0                     |
| 103             | 0.047            | 0.1                     | 1.51         | 1.0                     | 2.13         | 0.5                     | 0.029    | 1.0                     | 0.042           | 0.9                     | 0.096            | 0.6                     | 1.62      | 0.6                     | 0.//         | 0.0                     |

**Table 10b**. Results for the mandatory pesticides lambda-cyhalothrin, lindane, metconazole, proquinazid, prothioconazole-desthio, pyraclostrobin, tau-fluvalinate and teflubenzuron. The data for teflubenzuron is only for informative purposes because the level was below 3\*MRRL and consequently too low to evaluate.

| informat        |                    | rposes                  | s becau | se the                  | level v     | vas be                  | low 3*      | MKKL                    | and c            | onseq             | uently         | too lo                  | w to ev         | aluate                  | е.            |                         |
|-----------------|--------------------|-------------------------|---------|-------------------------|-------------|-------------------------|-------------|-------------------------|------------------|-------------------|----------------|-------------------------|-----------------|-------------------------|---------------|-------------------------|
| Laboratory code | Lambda-cyhalothrin | Z-scores (FFP RSD (25%) | Lindane | Z-scores (FFP RSD (25%) | Metconazole | Z-scores (FFP RSD (25%) | Proquinazid | Z-scores (FFP RSD (25%) | Prothioconazole- | es (FFP RSD (25%) | Pyraclostrobin | Z-scores (FFP RSD (25%) | Tau-Fluvalinate | Z-scores (FFP RSD (25%) | Teffubenzuron | Z-scores (FFP RSD (25%) |
| MRRL            | 0.01               | COL                     | 0.01    | COL                     | 0.01        | COL                     | 0.01        | COL                     | 0.01             | Z-scores (        | 0.01           | COL                     | 0.01            | COL                     | 0.01          | COL                     |
| Assigned value  | 0.061              | Z-s                     | 0.031   | S-Z                     | 0.035       | s-Z                     | 0.225       | s-Z                     | 1.37             | S-Z               | 4.39           | s-Z                     | 0.269           | s-Z                     | 0.016         | S-Z                     |
| 1               | 0.047              | -0.9                    | 0.027   | -0.6                    | 0.022       | -1.5                    | 0.235       | 0.2                     | 1.811            | 1.3               | 5.7            | 1.2                     | 0.329           | 0.9                     | 0.011         | -1.2                    |
| 2               | 0.067              | 0.4                     | 0.035   | 0.5                     | 0.037       | 0.2                     | 0.222       | -0.1                    | 1.492            | 0.4               | 4.229          | -0.2                    | 0.296           | 0.4                     | 0.018         | 0.5                     |
| 4               | 0.0704             | 0.6                     | 0.0349  | 0.4                     | 0.0356      | 0.1                     | 0.231       | 0.1                     | 1.43             | 0.2               | 3.25           | -1.0                    | 0.277           | 0.1                     | 0.0193        | 0.8                     |
| 5               | 0.059              | -0.1                    | FN      | -3.5                    | 0.036       | 0.1                     | 0.63        | >5                      | 1.09             | -0.8              | 3.85           | -0.5                    | 0.279           | 0.1                     | 0.014         | -0.5                    |
| 6               | 0.058              | -0.2                    | 0.03    | -0.2                    |             |                         | 0.178       | -0.8                    |                  |                   | 3.567          | -0.8                    | 0.256           | -0.2                    | 0.012         | -1.0                    |
| 8               | 0.068              | 0.5                     | 0.028   | -0.4                    | 0.031       | -0.5                    | 0.245       | 0.4                     | 1.45             | 0.2               | 4.25           | -0.1                    | 0.268           | 0.0                     | 0.0187        | 0.7                     |
| 9               | 0.0839             | 1.5                     | 0.022   | -1.2                    | 0.033       | -0.2                    | 0.21        | -0.3                    | 1.34             | -0.1              | 4.78           | 0.4                     | 0.286           | 0.3                     | 0.015         | -0.2                    |
| 10              | 0.0676             | 0.4                     | 0.0303  | -0.1                    | 0.0248      | -1.2                    | 0.13        | -1.7                    | 0.683            | -2.0              | 2.43           | -1.8                    | 0.279           | 0.1                     | 0.0138        | -0.5                    |
| 11              | 0.0458             | -1.0                    | 0.0168  | -1.9                    |             |                         |             |                         |                  |                   | 2.43           | -1.8                    | 0.235           | -0.5                    |               |                         |
| 12              | 0.072              | 0.7                     | 0.03    | -0.2                    | 0.027       | -0.9                    |             |                         | 1.45             | 0.2               |                |                         | 0.285           | 0.2                     |               |                         |
| 13              | 0.0688             | 0.5                     | 0.0384  | 0.9                     | 0.0382      | 0.4                     | 0.251       | 0.5                     | 1.49             | 0.4               | 5.01           | 0.6                     | 0.313           | 0.7                     | 0.0151        | -0.2                    |
| 14              | 0.075              | 0.9                     | 0.033   | 0.2                     | 0.028       | -0.8                    | 0.223       | 0.0                     | 1.556            | 0.6               | 4.379          | 0.0                     | 0.299           | 0.4                     | 0.016         | 0.0                     |
| 15              | 0.06               | -0.1                    | 0.034   | 0.3                     | 0.037       | 0.2                     | 0.212       | -0.2                    | 1.319            | -0.1              | 4.31           | -0.1                    | 0.299           | 0.4                     | 0.013         | -0.7                    |
| 16              | 0.085              | 1.6                     | 0.037   | 0.7                     | 0.028       | -0.8                    | 0.25        | 0.4                     | 1.511            | 0.4               | 4.563          | 0.2                     | 0.307           | 0.6                     | 0.017         | 0.3                     |
| 17              |                    |                         |         |                         |             |                         |             |                         |                  |                   |                |                         |                 |                         |               |                         |
| 18              | 0.0691             | 0.5                     | 0.0299  | -0.2                    | 0.0278      | -0.8                    |             |                         | 1.38             | 0.0               |                |                         | 0.327           | 0.9                     |               |                         |
| 19              | 0.021              | -2.6                    | 0.015   | -2.1                    |             |                         |             |                         | 0.217            | -3.4              |                |                         | 0.075           | -2.9                    |               |                         |
| 20              |                    |                         | 0.047   | 2.0                     |             |                         |             |                         |                  |                   |                |                         |                 |                         |               |                         |
| 21              | 0.064              | 0.2                     | 0.034   | 0.3                     | 0.033       | -0.2                    | 0.347       | 2.2                     | 1.269            | -0.3              | 5.503          | 1.0                     | 0.306           | 0.5                     | 0.026         | 2.5                     |
| 22              | 0.061              | 0.0                     | 0.035   | 0.5                     | 0.032       | -0.3                    | 0.236       | 0.2                     | 1.65             | 0.8               | 5.92           | 1.4                     | 0.313           | 0.7                     | 0.012         | -1.0                    |
| 23              | 0.07               | 0.6                     | 0.038   | 0.8                     | 0.048       | 1.5                     | 0.254       | 0.5                     | 1.45             | 0.2               | 4.68           | 0.3                     | 0.362           | 1.4                     | 0.019         | 0.8                     |
| 24              | 0.046              | -1.0                    | 0.036   | 0.6                     | 0.039       | 0.4                     | 0.194       | -0.6                    | 1.12             | -0.7              | 4.9            | 0.5                     | 0.221           | -0.7                    | 0.015         | -0.2                    |
| 25              |                    |                         | 0.036   | 0.6                     |             |                         |             |                         |                  |                   |                |                         |                 |                         |               |                         |
| 26              | 0.0716             | 0.7                     | 0.0358  | 0.6                     | 0.0363      | 0.1                     | 0.167       | -1.0                    | 1.08             | -0.8              | 5.71           | 1.2                     | 0.277           | 0.1                     | 0.0159        | 0.0                     |
| 27              | 0.054              | -0.4                    | 0.03    | -0.2                    | 0.03        | -0.6                    | 0.223       | 0.0                     | 1.414            | 0.1               | 4.489          | 0.1                     | 0.289           | 0.3                     | 0.012         | -1.0                    |
| 28              | 0.073              | 0.8                     | 0.038   | 0.8                     | 0.023       | -1.4                    | 0.244       | 0.3                     | 1.47             | 0.3               | 5.16           | 0.7                     | 0.286           | 0.3                     | 0.024         | 2.0                     |
| 29              | 0.07               | 0.6                     | 0.047   | 2.0                     | 0.087       | >5                      | 0.222       | -0.1                    | 1.2              | -0.5              | 3.939          | -0.4                    | 0.32            | 8.0                     | 0.02          | 1.0                     |
| 30              | 0.0643             | 0.2                     | 0.035   | 0.45                    | 0.036       | 0.1                     | 0.195       | -0.5                    | 1.45             | 0.2               | 4.55           | 0.1                     | 0.264           | -0.1                    | 0.0103        | -1.4                    |
| 31              | 0.042              | -1.2                    | 0.021   | -1.3                    | 0.076       | 4.7                     | 0.176       | -0.9                    | 1.05             | -0.9              | 4.2            | -0.2                    | 0.292           | 0.3                     | 0.011         | -1.2                    |
| 32              |                    |                         | 0.0206  | -1.4                    |             |                         |             |                         |                  |                   |                |                         |                 |                         |               |                         |
| 33              | 0.071              | 0.7                     | 0.027   | -0.6                    | 0.053       | 2.0                     |             |                         | 1.59             | 0.7               |                |                         | 0.29            | 0.3                     |               |                         |
| 34              | 0.065              | 0.3                     | 0.033   | 0.2                     | 0.03        | -0.6                    | 0.189       | -0.6                    | 0.852            | -1.5              | 3.026          | -1.2                    | 0.16            | -1.6                    | 0.014         | -0.5                    |
| 35              | 0.062              | 0.1                     | 0.048   | 2.1                     | 0.046       | 1.2                     | 0.222       | -0.1                    | 1.769            | 1.2               | 4.474          | 0.1                     | 0.3             | 0.5                     | 0.018         | 0.5                     |
| 36              | 0.06               | -0.1                    | 0.027   | -0.6                    | 0.028       | -0.8                    | 0.188       | -0.7                    | 1.048            | -0.9              | 2.953          | -1.3                    | 0.229           | -0.6                    | 0.013         | -0.7                    |
| 37              | 0.035              | -1.7                    | 0.024   | -0.9                    | 0.033       | -0.2                    | 0.16        | -1.2                    | 1.1              | -0.8              | 2.9            | -1.4                    | 0.22            | -0.7                    | 0.017         | 0.3                     |
| 38              | 0.079              | 1.2                     |         |                         |             |                         |             |                         |                  |                   | 4.14           | -0.2                    |                 |                         | 0.016         | 0.0                     |
| 39              | 0.057              | -0.3                    | 0.034   | 0.3                     |             |                         |             |                         |                  |                   | 6.08           | 1.5                     | 0.275           | 0.1                     | 0.013         | -0.7                    |
| 40              | 0.0684             | 0.5                     | 0.0345  | 0.4                     |             |                         |             |                         |                  |                   |                |                         |                 |                         |               |                         |

| Laboratory code | 10.0 Lambda-cyhalothrin | Z-scores (FFP RSD (25%) | Tindane | Z-scores (FFP RSD (25%) | Metconazole | Z-scores (FFP RSD (25%) | Proquinazid | Z-scores (FFP RSD (25%) | Prothioconazole- | Z-scores (FFP RSD (25%) | Pyraclostrobin | Z-scores (FFP RSD (25%) | 10.0 Tau-Fluvalinate | Z-scores (FFP RSD (25%) | 0.01   | Z-scores (FFP RSD (25%) |
|-----------------|-------------------------|-------------------------|---------|-------------------------|-------------|-------------------------|-------------|-------------------------|------------------|-------------------------|----------------|-------------------------|----------------------|-------------------------|--------|-------------------------|
| Assigned value  | 0.061                   | ż                       | 0.031   | ·Z                      | 0.035       | ·Z                      | 0.225       | ·Z                      | 1.37             | ·z                      | 4.39           | ż                       | 0.269                | ·Z                      | 0.016  | ż                       |
| 41              | 0.064                   | 0.2                     | FN      | -3.5                    | 0.059       | 2.7                     | 0.292       | 1.2                     |                  |                         | 4.389          | 0.0                     | 0.391                | 1.8                     | 0.006  | -2.5                    |
| 42              | 0.0926                  | 2.1                     | 0.0458  | 1.8                     | 0.0362      | 0.1                     | 0.054       |                         | 1.6              | 0.7                     | 5.26           | 0.8                     | 2.000                |                         | 0.014  | -0.5                    |
| 43              | 0.071                   | -3.0                    | 0.032   | 0.1                     | 0.034       | -0.1                    | 0.274       | 0.9                     | 1.75<br>0.968    | 1.1                     | 6.83           | -2.2                    | 0.299                | 0.4                     | 0.016  | 0.0                     |
| 45              | 0.013                   | -0.6                    | 0.01    | 0.2                     | 0.018       | -0.9                    | 0.280       | -2.8                    | 1.115            | -0.7                    | 2.4            | -1.8                    | 0.22                 | -0.7                    | 0.102  | >5                      |
| 46              | 0.0219                  | -2.6                    | 0.0265  | -0.6                    | 0.0174      | -2.0                    | 0.338       | 2.0                     | 1.28             | -0.3                    | 4.24           | -0.1                    | 0.0988               | -2.5                    | 0.021  | 1.3                     |
| 47              | 0.032                   | -1.9                    | 0.03    | -0.2                    |             |                         |             |                         |                  |                         |                |                         | 0.138                | -1.9                    |        |                         |
| 48              |                         |                         | 0.0323  | 0.1                     |             |                         |             |                         |                  |                         |                |                         |                      |                         |        |                         |
| 49              | 0.041                   | -1.3                    | 0.028   | -0.4                    | 0.025       | -1.1                    |             |                         |                  |                         | 1.986          | -2.2                    | 0.21                 | -0.9                    |        |                         |
| 50              | 0.045                   | -1.0                    | 0.0312  | 0.0                     | 0.0315      | -0.4                    | 0.225       | 0.0                     | 1.48             | 0.3                     | 4.67           | 0.3                     | 0.316                | 0.7                     | 0.0162 | 0.1                     |
| 51              | 0.09                    | 1.9                     | 0.038   | 0.8                     | 0.0294      | -0.6                    | 0.267       | 0.7                     | 1.16             | -0.6                    | 4.43           | 0.0                     | 0.508                | 3.6                     | 0.0119 | -1.0                    |
| 52              | 0.052                   | -0.6                    | 0.036   | 0.6                     | 0.043       | 0.9                     | 0.198       | -0.5                    | 0.978            | -1.1                    | 2.868          | -1.4                    | 0.206                | -0.9                    |        |                         |
| 53              | 0.036                   | -1.6                    | 0.03    | -0.2                    | 0.000       | 0.0                     | 0.064       | -2.9                    | 0.46             | -2.7                    | 2.36           | -1.9                    | 0.068                | -3.0                    |        |                         |
| 54              | 0.045                   | -1.0                    | 0.024   | -0.9                    | 0.032       | -0.3                    | 0.298       | 1.3                     | 1.11             | -0.8                    | FN<br>5.409    | -4.0<br>0.9             | 0.246                | -0.3                    | 0.022  | 1.5                     |
| 55<br>56        | 0.045                   | -0.1                    | 0.024   | -0.9                    | 0.041       | 0.7                     | 0.298       | 1.7                     | 1.61<br>1.812    | 1.3                     | 6.888          | 2.3                     | 0.194                | -0.3                    | 0.022  | 0.3                     |
| 57              | 0.039                   | -0.1                    | 0.028   | -0.4                    | 0.037       | 0.2                     | 0.322       | 1./                     | 1.012            | 1.5                     | 0.000          | 2.3                     | 0.194                | -1.1                    | 0.017  | 0.5                     |
| 58              | 0.037                   | -1.6                    | 0.03    | -0.2                    | 0.041       | 0.7                     | 0.319       | 1.7                     | 1.039            | -1.0                    | 5.159          | 0.7                     | 0.253                | -0.2                    | 0.017  | 0.3                     |
| 59              |                         |                         | 0.024   | -0.9                    |             |                         |             |                         |                  |                         |                |                         |                      |                         |        |                         |
| 60              |                         |                         | 0.0244  | -0.9                    |             |                         |             |                         |                  |                         |                |                         |                      |                         |        |                         |
| 61              | 0.065                   | 0.3                     | 0.024   | -0.9                    | 0.028       | -0.8                    | 0.201       | -0.4                    | 0.97             | -1.2                    | 4.69           | 0.3                     | 0.223                | -0.7                    | 0.019  | 0.8                     |
| 62              | 0.067                   | 0.4                     | 0.038   | 0.8                     | 0.034       | -0.1                    | 0.18        | -0.8                    | 1.34             | -0.1                    | 4.81           | 0.4                     | 0.243                | -0.4                    | 0.013  | -0.7                    |
| 63              | 0.064                   | 0.2                     | 0.035   | 0.5                     | 0.038       | 0.3                     | 0.26        | 0.6                     | 1.86             | 1.4                     | 6.25           | 1.7                     | 0.312                | 0.6                     |        |                         |
| 64              | 0.072                   | 0.7                     | 0.029   | -0.3                    | 0.028       | -0.8                    | 0.204       | -0.4                    | 1.732            | 1.1                     | 3.864          | -0.5                    | 0.278                | 0.1                     | 0.009  | -1.7                    |
| 65              | 0.0555                  | 0.0                     | 0.038   | 0.8                     | 0.0454      | 4.0                     | 0.46        | 4.0                     | 0.04             | 4.0                     | 2.0            | 4.5                     | 0.224                | 0.6                     | 0.015  | 0.2                     |
| 66              | 0.0555                  | -0.3                    | 0.0293  | -0.3                    | 0.0454      | 1.2                     | 0.16        | -1.2                    | 0.91<br>1.81     | -1.3<br>1.3             | 2.8<br>4.7     | -1.5                    | 0.231                | -0.6                    | 0.015  | -0.2                    |
| 68              | 0.0743                  | 0.7                     | 0.035   | 0.5                     |             |                         |             |                         | 1.01             | 1.5                     | 4.9            | 0.5                     | 0.39                 | 1.8                     |        |                         |
| 69              | 0.082                   | 1.4                     | 0.033   | 0.2                     | 0.05        | 1.7                     | 0.345       | 2.1                     | 45.18            | >5                      | 6.35           | 1.8                     | FN                   | -3.9                    | 0.024  | 2.0                     |
| 70              | 0.0596                  | -0.1                    | FN      | -3.5                    | 0.0138      | -2.4                    | 0.194       | -0.6                    | 1.9              | 1.6                     | 3.84           | -0.5                    | 0.28                 | 0.2                     | 0.018  | 0.5                     |
| 71              | 0.068                   | 0.5                     | 0.035   | 0.5                     | 0.011       | -2.7                    | 0.17        | -1.0                    | 0.9              | -1.4                    | 4              | -0.4                    | 0.33                 | 0.9                     | 0.016  | 0.0                     |
| 72              | 0.063                   | 0.1                     | 0.03    | -0.2                    | 0.026       | -1.0                    | 0.229       | 0.1                     | 1.71             | 1.0                     | 5.6            | 1.1                     | 0.364                | 1.4                     |        |                         |
| 73              | 0.025                   | -2.4                    |         |                         | 0.02        | -1.7                    | 0.16        | -1.2                    |                  |                         | 3.46           | -0.9                    |                      |                         |        |                         |
| 74              | 0.0687                  | 0.5                     | 0.0334  | 0.2                     | 0.0335      | -0.2                    |             |                         | FN               | -4.0                    |                |                         | 0.306                | 0.5                     | 0.0144 | -0.4                    |
| 75              | 0.047                   | -0.9                    | 0.034   | 0.3                     | 0.043       | 0.9                     | 0.248       | 0.4                     | 1.705            | 1.0                     | 6.123          | 1.6                     | 0.273                | 0.1                     | 0.017  | 0.3                     |
| 76              | 0.097                   | 2.4                     |         |                         | 0.049       | 1.6                     | 0.268       | 0.8                     | 0.502            | -2.5                    | 5.322          | 0.8                     | 0.376                | 1.6                     | 0.018  | 0.5                     |
| 77              | 0.0526                  | -0.5                    | 0.0327  | 0.2                     | 0.0684      | 3.8                     | 0.28        | 1.0                     | 2.45             | 3.2                     | 4.98           | 0.5                     | 0.204                | -1.0                    | 0.0492 | >5                      |
| 79              | 0.0326                  | 0.5                     | 0.0327  | -0.9                    | 0.038       | 0.3                     | 0.214       | -0.2                    | 1.641            | 0.8                     | 4.341          | 0.0                     | 0.233                | -0.5                    | 0.0492 | 0.0                     |
| 80              | 0.0753                  | 1.0                     | 0.034   | 0.3                     | 0.0466      | 1.3                     |             |                         | 1.25             | -0.3                    |                |                         | 0.253                | -0.2                    |        |                         |
|                 |                         |                         |         |                         |             |                         |             |                         |                  |                         |                |                         |                      |                         |        |                         |

|                 | <u> </u>           |                         |         |                         |             |                         |             |                         |                            |                         |                |                         |                 |                         |               |                         |
|-----------------|--------------------|-------------------------|---------|-------------------------|-------------|-------------------------|-------------|-------------------------|----------------------------|-------------------------|----------------|-------------------------|-----------------|-------------------------|---------------|-------------------------|
| Laboratory code | Lambda-cyhalothrin | (25%)                   | Lindane | (25%)                   | Metconazole | (25%)                   | Proquinazid | (25%)                   | Prothioconazole-<br>destio | (25%)                   | Pyraclostrobin | (25%)                   | Tau-Fluvalinate | (25%)                   | Teflubenzuron | (25%)                   |
| Labora          | Lambda-            | Z-scores (FFP RSD (25%) | Lir     | Z-scores (FFP RSD (25%) | Metc        | Z-scores (FFP RSD (25%) | Proq        | Z-scores (FFP RSD (25%) | Prothio<br>de              | Z-scores (FFP RSD (25%) | Pyrac          | Z-scores (FFP RSD (25%) | Tau-Fl          | Z-scores (FFP RSD (25%) | Teflub        | Z-scores (FFP RSD (25%) |
| MRRL            | 0.01               | cores                   | 0.01    | cores                   | 0.01        | sores                   | 0.01        | cores                   | 0.01                       | cores                   | 0.01           | ores                    | 0.01            | cores                   | 0.01          | cores                   |
| Assigned value  | 0.061              | Z-sc                    | 0.031   | Z-sc                    | 0.035       | Z-sc                    | 0.225       | Z-sc                    | 1.37                       | Z-sc                    | 4.39           | Z-sc                    | 0.269           | Z-sc                    | 0.016         | Z-sc                    |
| 81              | 0.074              | 0.9                     | 0.048   | 2.1                     | 0.039       | 0.4                     | 0.392       | 3.0                     | 3.283                      | >5                      | 6.083          | 1.5                     | 0.335           | 1.0                     |               |                         |
| 82              | 0.0643             | 0.2                     | 0.0349  | 0.4                     | 0.0291      | -0.7                    | 0.187       | -0.7                    | 1.51                       | 0.4                     | 3.46           | -0.9                    | 0.279           | 0.1                     | 0.0122        | -0.9                    |
| 83              | 0.0267             | -2.2                    | 0.0232  | -1.0                    | 0.0264      | -1.0                    |             |                         |                            |                         | 1.329          | -2.8                    | 0.0786          | -2.8                    |               |                         |
| 84              | 0.071              | 0.7                     | 0.025   | -0.8                    | 0.029       | -0.7                    | 0.207       | -0.3                    | 0.999                      | -1.1                    | 5.04           | 0.6                     | 0.181           | -1.3                    | 0.017         | 0.3                     |
| 85              | 0.065              | 0.3                     | 0.025   | -0.8                    | 0.069       | 3.9                     | 0.00=       |                         | 1.36                       | 0.0                     | 0.404          |                         | 0.241           | -0.4                    | 0.00          | 4.0                     |
| 86              | 0.067              | 0.4                     | 0.010   | 1.7                     | 0.033       | -0.2                    | 0.205       | -0.4                    | 1.111                      | -0.8                    | 3.484          | -0.8                    | 0.26            | -0.1                    | 0.02          | 1.0                     |
| 87<br>88        | 0.067              | -0.2                    | 0.018   | -1.7                    | 0.022       | -1.5                    | 0.215       | -0.2                    | FN<br>1.28                 | -4.0                    | 2.7<br>4.62    | -1.5<br>0.2             | 0.75            | >5                      | 0.0145        | -0.4                    |
| 89              | 0.0383             | -0.2                    | 0.0310  | -0.1                    | 0.0414      | 0.7                     | 0.213       | 0.4                     | 1.23                       | -0.3                    | 3.56           | -0.8                    | 0.272           | -0.9                    | 0.0143        | 0.3                     |
| 90              | 0.048              | -0.8                    | 0.031   | -0.1                    | 0.036       | -0.5                    | 0.247       | -0.6                    | 1.23                       | 0.4                     | 4.2            | -0.8                    | 0.207           | 0.8                     | 0.017         | 0.3                     |
| 91              | 0.028              | -2.2                    | 0.031   | -0.1                    | 0.039       | 0.4                     | 0.26        | 0.6                     | FN                         | -4.0                    | 5.44           | 1.0                     | 0.173           | -1.4                    | 0.017         | 0.0                     |
| 92              | 0.064              | 0.2                     | 0.026   | -0.7                    | 0.068       | 3.8                     |             |                         | 1.286                      | -0.2                    |                |                         | 0.265           | -0.1                    |               |                         |
| 93              | 0.04               | -1.4                    | 0.038   | 0.8                     | 0.03        | -0.6                    | 0.198       | -0.5                    | 1.254                      | -0.3                    | 3.9            | -0.5                    | 0.23            | -0.6                    |               |                         |
| 94              | 0.0713             | 0.7                     | 0.0338  | 0.3                     | 0.0361      | 0.1                     | 0.218       | -0.1                    | 1.24                       | -0.4                    | 3.77           | -0.6                    | 0.341           | 1.1                     | 0.0162        | 0.1                     |
| 95              | 0.059              | -0.1                    | 0.03    | -0.2                    | 0.033       | -0.2                    | 0.184       | -0.7                    | 1.664                      | 0.9                     | 4.662          | 0.2                     | 0.276           | 0.1                     | 0.021         | 1.3                     |
| 96              | 0.072              | 0.7                     | 0.033   | 0.2                     |             |                         |             |                         |                            |                         | 3.5            | -0.8                    |                 |                         |               |                         |
| 97              | 0.0436             | -1.1                    | FN      | -2.7                    | 0.0288      | -0.7                    | 0.111       | -2.0                    | 0.844                      | -1.5                    | 2.42           | -1.8                    | 0.235           | -0.5                    | 0.0145        | -0.4                    |
| 98              | 0.0568             | -0.3                    | 0.0347  | 0.4                     | 0.0262      | -1.0                    | 0.182       | -0.8                    | 1.204                      | -0.5                    | 4.55           | 0.1                     | 0.215           | -0.8                    | 0.0131        | -0.7                    |
| 99              | 0.059              | -0.1                    | 0.039   | 1.0                     | FN          | -3.5                    | 0.277       | 0.9                     | 1.643                      | 0.8                     | 6.156          | 1.6                     | 0.313           | 0.7                     | 0.02          | 1.0                     |
| 100             | 0.055              | -0.4                    | 0.0308  | -0.1                    | 0.0214      | -1.6                    | 0.196       | -0.5                    | 1.17                       | -0.6                    | 4.29           | -0.1                    | 0.236           | -0.5                    | 0.0127        | -0.8                    |
| 101             | 0.0649             | 0.3                     | 0.0361  | 0.6                     | 0.0383      | 0.4                     | 0.365       | 2.5                     | 1.7                        | 1.0                     | 6.67           | 2.1                     | 0.29            | 0.3                     | 0.0149        | -0.3                    |
| 102             | 0.051              | -0.6                    | 0.032   | 0.07                    | 0.024       | -1.3                    | 0.162       | -1.1                    | 1.003                      | -1.1                    | 2.851          | -1.4                    | 0.378           | 1.6                     | 0.016         | 0.0                     |
| 103             | 0.052              | -0.6                    | 0.019   | -1.6                    | 0.038       | 0.3                     | 0.377       | 2.7                     | 1.581                      | 0.6                     | 3.868          | -0.5                    | 0.259           | -0.1                    | 0.015         | -0.2                    |
| 105             | 0.062              | 0.1                     | 0.035   | 0.5                     | DN          | 2.5                     | 0.25        | 0.4                     | 4.2                        | 0.2                     | 4.89           | 0.5                     | 0.27            | 0.0                     | 0.014         | -0.5                    |
| 106             | 0.03               | -2.0                    | 0.02    | -1.5                    | FN          | -3.5                    | 0.21        | -0.3                    | 1.3                        | -0.2                    | 4.5            | 0.1                     | 0.15            | -1.8                    | 0.016         | 0.0                     |
| 107<br>108      | 0.064              | 0.2<br>4.6              | 0.047   | 0.6                     | 0.04        | 0.6<br>>5               | 0.252       | 0.5<br>1.5              | 1.65<br>1.795              | 0.8                     | 6.13<br>5.752  | 1.6                     | 0.311           | 0.6<br>1.7              | 0.016         | 0.0                     |
| 109             | 0.0627             | 0.1                     | 0.0716  | >5                      | 0.0507      | 1.8                     | 0.244       | 0.3                     | 1.7 75                     | 1.2                     | 3.45           | -0.9                    | 0.337           | 1.0                     | 0.0184        | 0.6                     |
| 110             | 0.0561             | -0.3                    | 0.0291  | -0.3                    | 0.0308      | -0.5                    | 0.217       | -0.1                    | 1.54                       | 0.5                     | 4.58           | 0.2                     | 0.258           | -0.2                    | 0.0129        | -0.8                    |
| 111             | 0.054              | -0.4                    |         |                         | 0.033       | -0.2                    | 0.208       | -0.3                    | 1.451                      | 0.2                     | 5.26           | 0.8                     | 0.328           | 0.9                     | 0.021         | 1.3                     |
| 112             | 0.097              | 2.4                     | 0.021   | -1.3                    | FN          | -3.5                    |             |                         | 2.09                       | 2.1                     | 4.36           | 0.0                     | 0.478           | 3.1                     |               |                         |
| 113             | 0.091              | 2.0                     | 0.084   | >5                      | 0.031       | -0.5                    | 0.183       | -0.8                    | 1.078                      | -0.8                    | 4.135          | -0.2                    | 0.295           | 0.4                     |               |                         |
| 114             | 0.0663             | 0.4                     | 0.0324  | 0.1                     | 0.0352      | 0.0                     | 0.263       | 0.7                     | 1.52                       | 0.4                     | 4.67           | 0.3                     | 0.366           | 1.4                     | 0.0142        | -0.4                    |
| 115             | 0.057              | -0.3                    | 0.032   | 0.1                     | 0.033       | -0.2                    | 0.18        | -0.8                    | 1.21                       | -0.5                    | 3.82           | -0.5                    | 0.322           | 0.8                     | 0.012         | -1.0                    |
| 116             | 0.0692             | 0.6                     | 0.0308  | -0.1                    | 0.0329      | -0.2                    | 0.142       | -1.5                    | 1.31                       | -0.2                    | 2.83           | -1.4                    | 0.302           | 0.5                     | 0.0248        | 2.2                     |
| 117             | 0.0699             | 0.6                     | 0.0319  | 0.1                     | 0.037       | 0.2                     | 0.255       | 0.5                     | 1.43                       | 0.2                     | 5.09           | 0.6                     | 0.324           | 0.8                     | 0.0189        | 0.7                     |
| 118             |                    |                         |         |                         | 0.0335      | -0.2                    | 0.206       | -0.3                    | 1.71                       | 1.0                     | 4.202          | -0.2                    | 0.289           | 0.3                     | 0.013         | -0.7                    |
| 119             | 0.062              | 0.1                     | FN      | -3.5                    | FN          | -3.5                    | 0.278       | 0.9                     | 1.717                      | 1.0                     | 6.184          | 1.6                     | 0.367           | 1.5                     | 0.031         | 3.8                     |
| 120             | 0.067              | 0.4                     |         |                         | 0.033       | -0.2                    | 0.2         | -0.4                    | 1.33                       | -0.1                    | 4.14           | -0.2                    | 0.26            | -0.1                    | 0.012         | -1.0                    |
| 121             | 0.073              | 0.8                     | 0.035   | 0.5                     | 0.029       | -0.7                    | 0.31        | 1.5                     | 1.13                       | -0.7                    | 3.46           | -0.9                    | 0.202           | -1.0                    |               |                         |

| Laboratory code | Lambda-cyhalothrin | Z-scores (FFP RSD (25%) | Lindane | Z-scores (FFP RSD (25%) | Metconazole | Z-scores (FFP RSD (25%) | Proquinazid    | Z-scores (FFP RSD (25%) | Prothioconazole-<br>destio | Z-scores (FFP RSD (25%) | Pyraclostrobin | Z-scores (FFP RSD (25%) | Tau-Fluvalinate | Z-scores (FFP RSD (25%) | Teflubenzuron | Z-scores (FFP RSD (25%) |
|-----------------|--------------------|-------------------------|---------|-------------------------|-------------|-------------------------|----------------|-------------------------|----------------------------|-------------------------|----------------|-------------------------|-----------------|-------------------------|---------------|-------------------------|
| MRRL            | 0.01               | core                    | 0.01    | COTE                    | 0.01        | core                    | 0.01           | core                    | 0.01                       | core                    | 0.01           | core                    | 0.01            | core                    | 0.01          | core                    |
| Assigned value  | 0.061              | s-Z                     | 0.031   | s-z                     | 0.035       | s-Z                     | 0.225          | z-s                     | 1.37                       | s-z                     | 4.39           | s-Z                     | 0.269           | s-Z                     | 0.016         | s-Z                     |
| 122             | 0.137              | 5.0                     | 0.0112  | -2.6                    | 0.034       | -0.1                    | 0.221          | -0.1                    | 1.2                        | -0.5                    | 3.49           | -0.8                    | 0.0292          | -3.6                    | 0.149         | >5                      |
| 123             | 0.118              | 3.8                     | 0.0315  | 0.0                     | 0.0311      | -0.5                    | 0.151          | -1.3                    | 1.71                       | 1.0                     | 5.67           | 1.2                     | 0.281           | 0.2                     | 0.0131        | -0.7                    |
| 124             | 0.0547             | -0.4                    | 0.0269  | -0.6                    | 0.0429      | 0.9                     | 0.238          | 0.2                     | 1.16                       | -0.6                    | 4.33           | -0.1                    | 0.302           | 0.5                     | 0.0147        | -0.3                    |
| 125             |                    |                         | 0.033   | 0.2                     |             |                         |                |                         |                            |                         |                |                         |                 |                         |               |                         |
| 126             | 0.152              | >5                      | 0.058   | 3.4                     | FN          | -3.5                    | 0.258          | 0.6                     | 1.282                      | -0.3                    | 3.801          | -0.5                    | 0.287           | 0.3                     |               |                         |
| 127             | 0.051              | -0.6                    |         |                         |             |                         | 0.212          | -0.2                    | 1.08                       | -0.8                    | 3.92           | -0.4                    | 0.285           | 0.2                     |               |                         |
| 128             |                    |                         | 0.0274  | -0.5                    | 0.07        | 4.0                     |                |                         |                            |                         | 0.856          | -3.2                    |                 |                         |               |                         |
| 129             | 0.108              | 3.1                     | 0.0289  | -0.3                    | FN          | -3.5                    | 0.212          | -0.2                    | 0.665                      | -2.1                    | 1.69           | -2.5                    | 0.241           | -0.4                    |               |                         |
| 130             | 0.066              | 0.3                     | 0.034   | 0.3                     | 0.032       | -0.3                    | 0.271          | 0.8                     | 1.42                       | 0.2                     | 5.57           | 1.1                     |                 |                         |               |                         |
| 131             | 0.066              | 0.3                     | 0.023   | -1.1                    | 0.072       | 4.2                     |                |                         | 1.54                       | 0.5                     |                |                         | 0.22            | -0.7                    |               |                         |
| 132             | 0.048              | -0.8                    | 0.032   | 0.1                     | 0.068       | 3.8                     | 0.19           | -0.6                    | 1.02                       | -1.0                    | 3.2            | -1.1                    | 0.32            | 0.8                     | 0.011         | -1.2                    |
| 133             | 0.059              | -0.1                    | 0.036   | 0.6                     | 0.042       | 8.0                     | 0.27           | 0.8                     | 1.29                       | -0.2                    | 4.81           | 0.4                     | 0.28            | 0.2                     | 0.013         | -0.7                    |
| 134             | 0.061              | 0.0                     | 0.033   | 0.2                     | 0.034       | -0.1                    | 0.115          | -2.0                    | 0.86                       | -1.5                    | 3.917          | -0.4                    | 0.192           | -1.1                    | 0.015         | -0.2                    |
| 135             | 0.045              | -1.0                    | FN      | -3.5                    | FN          | -3.5                    |                |                         | 3.87                       | >5                      | 4.58           | 0.2                     | 0.255           | -0.2                    |               |                         |
| 136             | 0.061              | 0.0                     | 0.038   | 0.8                     |             |                         |                |                         |                            |                         |                |                         | 0.167           | -1.5                    |               |                         |
| 137             | 0.0616             | 0.1                     | 0.0332  | 0.2                     | 0.0357      | 0.1                     | 0.179          | -0.8                    | 1.05                       | -0.9                    | 2.92           | -1.3                    | 0.301           | 0.5                     | 0.0164        | 0.1                     |
| 138             | 0.0556             | -0.3                    | 0.0265  | -0.6                    | 0.0414      | 0.7                     | 0.235          | 0.2                     | 1.22                       | -0.4                    | 4.36           | 0.0                     | 0.225           | -0.7                    | 0.0161        | 0.0                     |
| 139             | 0.0654             | 0.3                     | 0.0341  | 0.3                     | 0.0318      | -0.4                    | 0.249          | 0.4                     | 1.34                       | -0.1                    | 5.24           | 0.8                     | 0.309           | 0.6                     | 0.0183        | 0.6                     |
| 140             | 0.0628             | 0.1                     | 0.0365  | 0.6                     | 0.0289      | -0.7                    | 0.359          | 2.4                     | 1.65                       | 0.8                     | 5.18           | 0.7                     | 0.297           | 0.4                     | 0.0141        | -0.5                    |
| 141             | 0.05               | -0.7                    | 0.029   | -0.3                    | 0.031       | -0.5                    | 0.2            | -0.4                    | 1.27                       | -0.3                    | 3.46           | -0.9                    | 0.17            | -1.5                    |               |                         |
| 142             | 0.04               | -1.4                    | 0.028   | -0.4                    | 0.042       | 0.8                     | 0.22           | -0.1                    | 1.886                      | 1.5                     | 7.531          | 2.9                     | 0.218           | -0.8                    | 0.012         | -1.0                    |
| 143             | 0.056              | -0.3                    | 0.025   | -0.8                    | 0.044       | 1.0                     | 0.2            | -0.4                    | 1.2                        | -0.5                    | 4.4            | 0.0                     | 0.22            | -0.7                    | 0.018         | 0.5                     |
| 144             | 0.105              | 2.9                     | 0.042   | 1.3                     |             |                         | 0.294          | 1.2                     |                            |                         |                |                         | 0.39            | 1.8                     |               |                         |
| 145             |                    |                         | 0.013   | -2.3                    | FN          | -3.5                    |                |                         |                            |                         |                |                         |                 |                         |               |                         |
| 146             | 0.063              | 0.1                     | 0.032   | 0.1                     | 0.043       | 0.9                     | 0.21           | -0.3                    | 1.2                        | -0.5                    | 4.5            | 0.1                     | 0.24            | -0.4                    | 0.017         | 0.3                     |
| 147             | 0.071              | 0.7                     | 0.033   | 0.2                     | 0.03        | -0.6                    | 0.348          | 2.2                     | 1.9                        | 1.6                     | 7.25           | 2.6                     | 0.328           | 0.9                     | 0.01          | -1.5                    |
| 148             | 0.0602             | 0.0                     | 0.0285  | -0.4                    | 0.0335      | -0.2                    | 0.25           | 0.4                     | 1.52                       | 0.4                     | 3.51           | -0.8                    | 0.239           | -0.4                    | 0.019         | 0.8                     |
| 149<br>150      | 0.051              | -0.6<br>0.5             | 0.028   | -0.4<br>1.3             | 0.074       | 0.6                     | 0.225<br>0.237 | 0.0                     | 1.3                        | -0.2<br>0.5             | 4.39<br>5.547  | 1.0                     | 0.281           | 0.2                     | 0.015         | -0.2<br>0.3             |
| 151             | 0.053              | -0.5                    | 0.042   | -0.2                    | 0.035       | 0.0                     | 0.237          | -0.4                    | 1.39                       | 0.3                     | 4.9            | 0.5                     | 0.287           | 0.3                     | 0.017         | -0.7                    |
| 152             | 0.033              | 0.5                     | 0.03    | 0.2                     | 0.033       | 0.0                     | 0.2            | 0.1                     | 1.57                       | 0.1                     | 1.7            | 0.5                     | 0.207           | 0.5                     | 0.015         | 0.7                     |
| 153             | 0.0297             | -2.0                    | 0.0363  | 0.6                     | 0.0347      | 0.0                     | 0.21           | -0.3                    | 1.6                        | 0.7                     | 4.78           | 0.4                     | 0.102           | -2.5                    |               |                         |
| 154             | 0.015              | -3.0                    | 0.021   | -1.3                    | 0.015       | -2.3                    |                |                         | 0.556                      | -2.4                    | 1.839          | -2.3                    | 0.114           | -2.3                    |               |                         |
| 155             | 0.02               | -2.7                    | 0.03    | -0.2                    | 0.03        | -0.6                    | 0.03           | -3.5                    | 2.49                       | 3.3                     | 5.64           | 1.1                     | 0.03            | -3.6                    | 0.05          | >5                      |
| 156             |                    |                         |         |                         | 0.04        | 0.6                     |                |                         |                            |                         | 3.97           | -0.4                    |                 |                         |               |                         |
| 157             | 0.0516             | -0.6                    | FN      | -3.5                    | 0.0335      | -0.2                    | 0.243          | 0.3                     | 1.62                       | 0.7                     | 5              | 0.6                     | 0.284           | 0.2                     |               |                         |
| 158             | 0.0611             | 0.0                     | 0.034   | 0.3                     | 0.0707      | 4.1                     | 0.204          | -0.4                    | 1.4                        | 0.1                     | 4.43           | 0.0                     | 0.296           | 0.4                     | 0.017         | 0.3                     |
| 160             |                    |                         |         |                         |             |                         |                |                         |                            |                         |                |                         |                 |                         |               |                         |
| 161             | 0.0746             | 0.9                     | 0.0361  | 0.6                     | 0.0254      | -1.1                    | 0.314          | 1.6                     | 1.41                       | 0.1                     | 4.49           | 0.1                     | 0.235           | -0.5                    |               |                         |
| 162             | 0.02               | -2.7                    | 0.027   | -0.6                    | 0.027       | -0.9                    | 0.157          | -1.2                    | FN                         | -4.0                    | 6.15           | 1.6                     | 0.207           | -0.9                    |               |                         |
| 163             | 0.062              | 0.1                     | 0.046   | 1.9                     | 0.04        | 0.6                     | 0.3            | 1.3                     | 1.76                       | 1.1                     | 6.28           | 1.7                     | 0.2             | -1.0                    | 0.016         | 0.0                     |

 $\textbf{Table 10c}. \ \ \textbf{Results for voluntary compound benzovindiflupyr, fenpicoxamid and mefentrifluconazole, in mg/kg, and the corresponding z scores, MRRLs and the assigned values.}$ 

| code           | lupyr | 9,                     | mid          |                        | nazole | 9,                     |
|----------------|-------|------------------------|--------------|------------------------|--------|------------------------|
| ory            | ndif  | (25                    | oxa          | (25                    | luco   | (25)                   |
| oorat          | zovi  | RSD                    | Fenpicoxamid | RSD                    | ntrif  | RSD                    |
| Lat            | Ben   | FFP                    | l 🖁          | FFP                    | Mefe   | FFP                    |
| MRRL           | 0.01  | Z-scores (FFP RSD (25% | 0.01         | Z-scores (FFP RSD (25% | 0.01   | Z-scores (FFP RSD (25% |
| Assigned value | 0.718 | -Z                     | 0.040        | ž                      | 0.722  | -Z                     |
| 1              | 0.928 | 1.2                    | 0.035        | -0.5                   | 0.883  | 0.9                    |
| 2              | 0.693 | -0.1                   |              |                        | 0.913  | 1.1                    |
| 4              | 0.557 | -0.9                   |              |                        |        |                        |
| 5              |       |                        |              |                        |        |                        |
| 6              |       |                        |              |                        |        |                        |
| 8              |       |                        |              |                        |        |                        |
| 9              |       |                        |              |                        |        |                        |
| 10             | 0.368 | -1.9                   |              |                        | 0.345  | -2.1                   |
| 11             |       |                        |              |                        |        |                        |
| 12             |       |                        |              |                        |        |                        |
| 13             | 0.804 | 0.5                    | 0.0418       | 0.2                    | 0.693  | -0.2                   |
| 14             | 0.852 | 0.7                    | 0.037        | -0.3                   | 0.719  | 0.0                    |
| 15             | 0.705 | -0.1                   | 0.052        | 1.2                    | 0.722  | 0.0                    |
| 16             | 0.876 | 0.9                    | 0.039        | -0.1                   | 0.686  | -0.2                   |
| 17             |       |                        |              |                        |        |                        |
| 18             |       |                        |              |                        |        |                        |
| 19             |       |                        |              |                        |        |                        |
| 20             |       |                        |              |                        |        |                        |
| 21             | 0.886 | 0.9                    | 0.04         | 0.0                    | 0.852  | 0.7                    |
| 22             | 0.804 | 0.5                    | 0.039        | -0.1                   | 0.86   | 0.8                    |
| 23             |       |                        |              |                        |        |                        |
| 24             | 0.703 | -0.1                   |              |                        | 0.687  | -0.2                   |
| 25             |       |                        |              |                        |        |                        |
| 26             | 0.636 | -0.5                   |              |                        | 0.675  | -0.3                   |
| 27             | 0.731 | 0.1                    | 0.031        | -0.9                   | 0.792  | 0.4                    |
| 28             | 0.827 | 0.6                    |              |                        |        |                        |
| 29             |       |                        |              |                        |        |                        |
| 30             | 0.657 | -0.3                   | 0.0357       | -0.4                   | 0.74   | 0.1                    |
| 31             |       |                        |              |                        |        |                        |
| 32             |       |                        |              |                        |        |                        |
| 33             |       |                        |              |                        |        |                        |
| 34             | 0.528 | -1.1                   | 0.033        | -0.7                   | 0.504  | -1.2                   |
| 35             | 0.746 | 0.2                    | 0.068        | 2.8                    | 1.109  | 2.1                    |
| 36             | 0.582 | -0.8                   | FN           | -3.0                   |        |                        |
| 37             |       |                        |              |                        |        |                        |
| 38             |       |                        |              |                        |        |                        |
| 39             |       |                        |              |                        |        |                        |
| 40             |       |                        |              |                        |        |                        |

| TANA Laboratory code | 0.01  | Z-scores (FFP RSD (25%) | Fenpicoxamid | Z-scores (FFP RSD (25%) | 0.01  | Z-scores (FFP RSD (25%) |
|----------------------|-------|-------------------------|--------------|-------------------------|-------|-------------------------|
| Assigned value       | 0.718 | S-Z                     | 0.040        | s-Z                     | 0.722 | s-z                     |
| 41                   |       |                         |              |                         |       |                         |
| 42                   |       |                         |              |                         |       |                         |
| 43                   |       |                         |              |                         |       |                         |
| 44                   |       |                         |              |                         | 0.45  | -1.5                    |
| 45<br>46             |       |                         |              |                         |       |                         |
| 47                   |       |                         |              |                         |       |                         |
| 48                   |       |                         |              |                         |       |                         |
| 49                   |       |                         |              |                         |       |                         |
| 50                   |       |                         |              |                         |       |                         |
| 51                   |       |                         |              |                         |       |                         |
| 52                   |       |                         |              |                         |       |                         |
| 53                   | 0.23  | -2.7                    |              |                         | 0.24  | -2.7                    |
| 54<br>55             |       |                         |              |                         |       |                         |
| 56                   | 0.943 | 1.3                     | 0.034        | -0.6                    | 0.758 | 0.2                     |
| 57                   | 0.743 | 1.5                     | 0.034        | -0.0                    | 0.730 | 0.2                     |
| 58                   | 0.715 | 0.0                     | 0.09         | 5.0                     | 0.561 | -0.9                    |
| 59                   |       |                         |              |                         |       |                         |
| 60                   |       |                         |              |                         |       |                         |
| 61                   |       |                         |              |                         |       |                         |
| 62                   | 0.631 | -0.5                    | 0.039        | -0.1                    | 0.641 | -0.4                    |
| 63                   | 0.500 | 0.5                     | 0.004        | 0.6                     | 0.660 | 0.0                     |
| 64<br>65             | 0.599 | -0.7                    | 0.034        | -0.6                    | 0.669 | -0.3                    |
| 66                   |       |                         | FN           | -3.0                    |       |                         |
| 67                   | 0.812 | 0.5                     |              |                         |       |                         |
| 68                   |       |                         |              |                         |       |                         |
| 69                   |       |                         |              |                         |       |                         |
| 70                   |       |                         |              |                         |       |                         |
| 71                   |       |                         | 0.035        | -0.5                    | 0.48  | -1.3                    |
| 72                   | 0.48  | -1.3                    | 0.051        | 1.1                     | 1.01  | 1.6                     |
| 73<br>74             |       |                         |              |                         |       |                         |
| 75                   | 0.978 | 1.4                     | 0.049        | 0.9                     | 0.936 | 1.2                     |
| 76                   | 2.7.0 |                         | 2.017        | 3.7                     | 2.730 |                         |
| 77                   |       |                         |              |                         |       |                         |
| 78                   | 0.909 | 1.1                     |              |                         |       |                         |
| 79                   | 0.53  | -1.0                    | 0.035        | -0.5                    | 0.951 | 1.3                     |
| 80                   |       |                         |              |                         |       |                         |

|                 |                  |                         |              |                         |                     | _                       |
|-----------------|------------------|-------------------------|--------------|-------------------------|---------------------|-------------------------|
| Laboratory code | Benzovindiflupyr | Z-scores (FFP RSD (25%) | Fenpicoxamid | Z-scores (FFP RSD (25%) | Mefentrifluconazole | Z-scores (FFP RSD (25%) |
| MRRL            | 0.01             | ores                    | 0.01         | ores                    | 0.01                | ores                    |
| Assigned value  | 0.718            | Z-sc                    | 0.040        | Z-sc                    | 0.722               | Z-sc                    |
| 81              | 1.141            | 2.4                     |              |                         | 1.217               | 2.7                     |
| 82              |                  |                         |              |                         |                     |                         |
| 83              |                  |                         |              |                         |                     |                         |
| 84              | 0.617            | -0.6                    | 0.016        | -2.4                    | 0.615               | -0.6                    |
| 85              |                  |                         |              |                         |                     |                         |
| 86              |                  |                         |              |                         |                     |                         |
| 87              |                  |                         | 0.147        | >5                      |                     |                         |
| 88              |                  |                         |              |                         |                     |                         |
| 89              |                  |                         |              |                         |                     |                         |
| 90              | 0.62             | -0.5                    |              |                         |                     |                         |
| 91              |                  |                         |              |                         |                     |                         |
| 92              |                  |                         |              |                         |                     |                         |
| 93              | 0.732            | 0.1                     | 0.029        | -1.1                    | 0.637               | -0.5                    |
| 94              | 0.709            | 0.0                     | 0.0442       | 0.4                     | 0.522               | -1.1                    |
| 95              |                  |                         | FN           | -3.0                    | FN                  | -3.9                    |
| 96              |                  |                         |              |                         |                     |                         |
| 97              | 0.386            | -1.8                    |              |                         |                     |                         |
| 98              | 0.75             | 0.2                     | 0.0323       | -0.8                    | 0.707               | -0.1                    |
| 99              |                  |                         |              |                         |                     |                         |
| 100             | 0.716            | 0.0                     | 0.0331       | -0.7                    | 0.725               | 0.0                     |
| 101             | 0.729            | 0.1                     |              |                         | FN                  | -3.9                    |
| 102             | 0.513            | -1.1                    | 0.042        | 0.2                     | 0.483               | -1.3                    |
| 103             |                  |                         | FN           | -3.0                    | FN                  | -3.9                    |
| 105             |                  |                         |              |                         |                     |                         |
| 106             | 0.74             | 0.1                     |              |                         |                     |                         |
| 107             | 0.936            | 1.2                     |              |                         |                     |                         |
| 108             |                  |                         |              |                         |                     |                         |
| 109             | 0.814            | 0.5                     |              |                         |                     |                         |
| 110             | 0.669            | -0.3                    |              |                         |                     |                         |
| 111             |                  |                         |              |                         |                     |                         |
| 112             |                  |                         |              |                         |                     |                         |
| 113             |                  |                         | 0.042        | 0.2                     | 0.723               | 0.0                     |
| 114             | 0.843            | 0.7                     | 0.101        | >5                      | 0.957               | 1.3                     |
| 115             | 0.45             |                         | 0.0:-:       |                         | 0.4                 |                         |
| 116             | 0.454            | -1.5                    | 0.0409       | 0.1                     | 0.466               | -1.4                    |
| 117             | 0.5              |                         | 0.0407       | 0.1                     | 0.758               | 0.2                     |
| 118             | 0.727            | 0.1                     | 0.054        | 2.6                     | FILE                | 0.0                     |
| 119             | 0.999            | 1.6                     | 0.076        | 3.6                     | FN                  | -3.9                    |
| 120             | 0.61             | -0.6                    | 0.041        | 0.1                     | FN                  | -3.9                    |
| 121             |                  |                         |              |                         |                     |                         |

| ge             | pyr    |                        | 멸            |                         | zole  |                        |
|----------------|--------|------------------------|--------------|-------------------------|-------|------------------------|
| 00 A           | liflu  | 22%                    | ami          | 22%                     | cons  | 22%                    |
| ato            | vinc   |                        | Fenpicoxamid | 9                       | riffu | 9                      |
| ароі           | enzo   | P. R.                  | Fenj         | 7<br>23                 | fent  | F R                    |
|                | m<br>M | s (FF                  |              | s (FF                   | Me    | s (FF                  |
| MRRL           | 0.01   | Z-scores (FFP RSD (25% | 0.01         | Z-scores (FFP RSD (25%) | 0.01  | Z-scores (FFP RSD (259 |
| Assigned value | 0.718  | 3-Z                    | 0.040        | Z-:                     | 0.722 | -Z                     |
| 122            | 0.601  | -0.7                   |              |                         |       |                        |
| 123            | 0.757  | 0.2                    |              |                         |       |                        |
| 124            | 0.549  | -0.9                   | 0.0357       | -0.4                    |       |                        |
| 125            |        |                        |              |                         |       |                        |
| 126            | 0.802  | 0.5                    | 0.038        | -0.2                    | 0.775 | 0.3                    |
| 127            |        |                        |              |                         |       |                        |
| 128            |        |                        |              |                         |       |                        |
| 129            |        |                        |              |                         |       |                        |
| 130            | 0.773  | 0.3                    |              |                         |       |                        |
| 131            |        |                        |              |                         |       |                        |
| 132            | 0.75   | 0.2                    |              |                         |       |                        |
| 133            | 0.71   | 0.0                    |              |                         | 0.76  | 0.2                    |
| 134            | 0.617  | -0.6                   | 0.036        | -0.4                    | 0.591 | -0.7                   |
| 135            |        |                        |              |                         |       |                        |
| 136            |        |                        |              |                         |       |                        |
| 137            |        |                        |              |                         |       |                        |
| 138            |        |                        |              |                         |       |                        |
| 139            | 0.717  | 0.0                    | 0.0415       | 0.1                     |       |                        |
| 140            |        |                        |              |                         |       |                        |
| 141            | 0.82   | 0.6                    |              |                         |       |                        |
| 142            | 0.704  | -0.1                   | 0.033        | -0.7                    | 0.594 | -0.7                   |
| 143            |        |                        |              |                         |       |                        |
| 144            |        |                        |              |                         |       |                        |
| 145            |        |                        |              |                         |       |                        |
| 146            | 0.641  | -0.4                   | 0.036        | -0.4                    | 0.615 | -0.6                   |
| 147            |        |                        | 0.046        | 0.6                     |       |                        |
| 148<br>149     |        |                        | 0.046        | 0.6                     |       |                        |
| 150            |        |                        |              |                         |       |                        |
| 151            | 0.818  | 0.6                    | 0.042        | 0.2                     | 0.764 | 0.2                    |
| 152            | 0.010  | 0.0                    | 0.012        | Ų. <u>L</u>             | o r   | Ų.Z                    |
| 153            | 0.731  | 0.1                    |              |                         |       |                        |
| 154            |        |                        |              |                         |       |                        |
| 155            |        |                        |              |                         |       |                        |
| 156            |        |                        |              |                         |       |                        |
| 157            | 0.767  | 0.3                    | 0.0407       | 0.1                     | 0.87  | 0.8                    |
| 158            | 0.723  | 0.0                    | FN           | -3.0                    | FN    | -3.9                   |
| 160            |        |                        |              |                         |       |                        |
| 161            |        |                        |              |                         |       |                        |
| 162            |        |                        |              |                         |       |                        |
| 163            | 0.79   | 0.4                    |              |                         | 0.87  | 0.8                    |

# 3.3.3 Sum of Weighted Z scores (AZ2) – Category A

To be classified into Category A, the laboratories had to submit quantitative results for at least 90% of the compulsory pesticides present in the Test Item ( $\geq$ 13 pesticide residues, exclusive of any false negatives results), analyse for more than 90% of the compulsory pesticides on the target list and also report no false positive results. For the 87 EU and EFTA laboratories in Category A (59%), the results were additionally evaluated by calculating the Average of the Squared -Score (AZ²). Of the 87 participants 72 participants (83%) obtained AZ2 score at or below 2 (good), 6 participants (7%) obtained AZ2 values between 2-3 (satisfactory) and 9 participants (10%) obtained AZ2 values  $\geq$ 3 (unsatisfactory). An additional five laboratories from Third Countries were evaluated and classified into Category A. The AZ2 scores achieved by the labs can be seen in **Table 11**.

**Table 11.** Sum of Weighted z scores (AZ<sup>2</sup>) for laboratories in Category A, the number of pesticides detected and quantified by the laboratories, the number of false negatives reported and the classification as good, satisfactory and unsatisfactory. The table includes data for both EU and non-EU participants.

| and unsatisfactory. The table includes data for both EU and non-EU participants. |   |                                       |  |     |                |                |        |  |  |  |
|--|---|---------------------------------------|--|-----|----------------|----------------|--------|--|--|--|
| Lab code   | No. of<br>detected<br>mandatory<br>pesticides | Mandatory<br>pesticides<br>detected % | No. of<br>detected<br>volutary<br>pesticides | AZ2 | False negative | Classification | NRL    |  |  |  |
| 1  | 15  | 100                                   | 3  | 1.2 | 0              | Good           |        |  |  |  |
| 2  | 15  | 100                                   | 2  | 0.3 | 0              | Good           |        |  |  |  |
| 4  | 15  | 100                                   | 1  | 0.2 | 0              | Good           | NRL-CE |  |  |  |
| 8  | 15  | 100                                   | 0  | 0.1 | 0              | Good           |        |  |  |  |
| 13   | 15  | 100                                   | 3  | 0.2 | 0              | Good           |        |  |  |  |
| 14   | 15  | 100                                   | 3  | 0.4 | 0              | Good           |        |  |  |  |
| 15   | 15  | 100                                   | 3  | 0.1 | 0              | Good           | NRL-CF |  |  |  |
| 16   | 15  | 100                                   | 3  | 0.6 | 0              | Good           | NRL-CE |  |  |  |
| 21   | 15  | 100                                   | 2  | 0.5 | 0              | Good           |        |  |  |  |
| 22   | 15  | 100                                   | 3  | 0.4 | 0              | Good           |        |  |  |  |
| 23   | 13  | 87                                    | 0  | 1.5 | 1              | Good           |        |  |  |  |
| 24   | 15  | 100                                   | 2  | 0.3 | 0              | Good           |        |  |  |  |
| 26   | 15  | 100                                   | 2  | 0.3 | 0              | Good           |        |  |  |  |
| 27   | 15  | 100                                   | 2  | 0.2 | 0              | Good           |        |  |  |  |
| 28   | 14  | 93                                    | 1  | 1.1 | 1              | Good           | NRL-CF |  |  |  |
| 29   | 15  | 100                                   | 0  | 2.7 | 0              | Satisfactory   |        |  |  |  |
| 30   | 15  | 100                                   | 2  | 0.2 | 0              | Good           |        |  |  |  |
| 31   | 15  | 100                                   | 0  | 2.2 | 0              | Satisfactory   |        |  |  |  |
| 34   | 15  | 100                                   | 3  | 0.8 | 0              | Good           |        |  |  |  |
| 36   | 15  | 100                                   | 1  | 0.6 | 1              | Good           |        |  |  |  |
| 37   | 15  | 100                                   | 0  | 1.0 | 0              | Good           |        |  |  |  |
| 43   | 15  | 100                                   | 0  | 0.9 | 0              | Good           |        |  |  |  |
| 44   | 13  | 87                                    | 1  | 3.8 | 2              | Unsatisfactor  |        |  |  |  |
| 45   | 15  | 100                                   | 0  | 1.5 | 0              | Good           | NRL-CF |  |  |  |
| 46   | 14  | 93                                    | 0  | 1.9 | 0              | Good           |        |  |  |  |
| 50   | 15  | 100                                   | 0  | 0.1 | 0              | Good           | NRL-CE |  |  |  |
| 51   | 15  | 100                                   | 0  | 2.0 | 0              | Good           | NRL-CF |  |  |  |

| Lab code | No. of<br>detected<br>mandatory<br>pesticides | Mandatory<br>pesticides<br>detected % | No. of<br>detected<br>volutary<br>pesticides | AZ <sup>2</sup> | False negative | Classification | NRL    |
|----------|---|---------------------------------------|--|-----------------|----------------|----------------|--------|
| 52       | 13  | 87                                    | 1  | 2.1             | 2              | Satisfactory   |        |
| 55       | 15  | 100                                   | 0  | 1.0             | 0              | Good           | NRL-CF |
| 56       | 15  | 100                                   | 2  | 1.3             | 0              | Good           |        |
| 58       | 15  | 100                                   | 3  | 1.0             | 0              | Good           |        |
| 61       | 14  | 93                                    | 0  | 1.4             | 1              | Good           | NRL-FE |
| 62       | 15  | 100                                   | 2  | 0.5             | 0              | Good           |        |
| 63       | 15  | 100                                   | 1  | 1.3             | 0              | Good           |        |
| 64       | 15  | 100                                   | 2  | 0.5             | 0              | Good           |        |
| 66       | 15  | 100                                   | 1  | 0.9             | 1              | Good           |        |
| 70       | 14  | 93                                    | 0  | 1.4             | 1              | Good           |        |
| 71       | 14  | 93                                    | 3  | 1.7             | 1              | Good           |        |
| 75       | 15  | 100                                   | 3  | 0.7             | 0              | Good           |        |
| 78       | 15  | 100                                   | 0  | 3.4             | 0              | Unsatisfactor  |        |
| 79       | 15  | 100                                   | 2  | 0.3             | 0              | Good           |        |
| 81       | 15  | 100                                   | 1  | 4.9             | 0              | Unsatisfactor  | NRL-CF |
| 82       | 15  | 100                                   | 0  | 0.3             | 0              | Good           |        |
| 84       | 15  | 100                                   | 2  | 1.1             | 0              | Good           |        |
| 88       | 15  | 100                                   | 0  | 0.1             | 0              | Good           |        |
| 89       | 15  | 100                                   | 1  | 0.5             | 0              | Good           | NRL-CE |
| 90       | 15  | 100                                   | 1  | 0.2             | 0              | Good           |        |
| 93       | 14  | 93                                    | 3  | 1.1             | 1              | Good           | NRL-CF |
| 94       | 15  | 100                                   | 3  | 0.4             | 0              | Good           |        |
| 95       | 15  | 100                                   | 3  | 0.5             | 3              | Good           |        |
| 98       | 15  | 100                                   | 3  | 0.5             | 0              | Good           | NRL-CE |
| 99       | 13  | 87                                    | 0  | 1.8             | 2              | Good           |        |
| 100      | 15  | 100                                   | 2  | 0.3             | 0              | Good           | NRL-CF |
| 101      | 14  | 93                                    | 2  | 1.2             | 1              | Good           | NRL-CF |
| 102      | 14  | 93                                    | 3  | 1.5             | 0              | Good           | NRL-CF |
| 103      | 15  | 100                                   | 3  | 1.1             | 3              | Good           |        |
| 106      | 14  | 93                                    | 1  | 1.9             | 1              | Good           | NRL-CF |
| 107      | 14  | 93                                    | 0  | 2.6             | 0              | Satisfactory   |        |
| 108      | 13  | 87                                    | 0  | 5.0             | 0              | Unsatisfactor  |        |
| 110      | 15  | 100                                   | 1  | 0.3             | 0              | Good           | NRL-CF |
| 111      | 13  | 87                                    | 0  | 0.3             | 0              | Good           |        |
| 113      | 15  | 100                                   | 2  | 2.2             | 0              | Satisfactory   |        |
| 114      | 15  | 100                                   | 3  | 0.5             | 0              | Good           |        |
| 115      | 15  | 100                                   | 1  | 0.3             | 0              | Good           |        |
| 116      | 15  | 100                                   | 3  | 0.9             | 0              | Good           |        |

| Lab code | No. of<br>detected<br>mandatory<br>pesticides | Mandatory<br>pesticides<br>detected % | No. of<br>detected<br>volutary<br>pesticides | AZ <sup>2</sup> | False negative | Classification | NRL    |
|----------|---|---------------------------------------|--|-----------------|----------------|----------------|--------|
| 117      | 15  | 100                                   | 2  | 0.4             | 0              | Good           | NRL-FE |
| 119      | 13  | 87                                    | 3  | 3.1             | 3              | Unsatisfactor  |        |
| 121      | 15  | 100                                   | 0  | 0.5             | 0              | Good           |        |
| 122      | 15  | 100                                   | 1  | 4.4             | 0              | Unsatisfactor  |        |
| 123      | 15  | 100                                   | 0  | 2.7             | 0              | Satisfactory   |        |
| 124      | 14  | 93                                    | 1  | 0.6             | 0              | Good           |        |
| 126      | 14  | 93                                    | 3  | >5              | 1              | Unsatisfactor  |        |
| 129      | 14  | 93                                    | 1  | 3.5             | 1              | Unsatisfactor  |        |
| 130      | 13  | 87                                    | 1  | 0.3             | 0              | Good           |        |
| 132      | 15  | 100                                   | 0  | 1.2             | 0              | Good           |        |
| 133      | 15  | 100                                   | 1  | 0.2             | 0              | Good           | NRL-FE |
| 134      | 15  | 100                                   | 2  | 0.7             | 0              | Good           |        |
| 137      | 14  | 93                                    | 1  | 0.4             | 0              | Good           |        |
| 138      | 14  | 93                                    | 1  | 1.0             | 1              | Good           | NRL-CE |
| 139      | 15  | 100                                   | 1  | 0.2             | 0              | Good           | NRL-CF |
| 140      | 15  | 100                                   | 0  | 1.0             | 0              | Good           |        |
| 141      | 13  | 87                                    | 0  | 1.1             | 1              | Good           |        |
| 142      | 14  | 93                                    | 3  | 1.9             | 1              | Good           |        |
| 143      | 15  | 100                                   | 0  | 0.3             | 0              | Good           |        |
| 146      | 15  | 100                                   | 2  | 0.4             | 0              | Good           |        |
| 147      | 15  | 100                                   | 1  | 1.2             | 0              | Good           |        |
| 148      | 15  | 100                                   | 1  | 0.2             | 0              | Good           | NRL-CF |
| 149      | 14  | 93                                    | 0  | 1.7             | 0              | Good           |        |
| 150      | 13  | 87                                    | 1  | 0.6             | 0              | Good           |        |
| 151      | 15  | 100                                   | 2  | 0.2             | 0              | Good           |        |
| 155      | 15  | 100                                   | 0  | 4.0             | 0              | Unsatisfactor  |        |
| 158      | 15  | 100                                   | 2  | 1.2             | 2              | Good           |        |

The 67 laboratories that did not fulfil the requirements described above, were classified in Category B. The number of reported quantitative results, analysed compounds from the Target List and acceptable z scores as well as information on false negative and positive results are shown in **Table 12**. Five laboratories was moved from Category A to B due to false positive results. Four participants fulfilled the criteria of detecting 90% of the compulsory pesticides in the Test Item but did not fulfil the criteria of analysing for 90% of the compulsory pesticides on the Target List. Four participants analysed more than 90% of the pesticides on the Target List but reported <13 pesticides in the Test Item.

**Table 12.** Number and percentage of compulsory pesticides detected and quantified, number of compulsory compounds analysed from the Target List, number of voluntary pesticides detected and quantified, number of acceptable z scores, false negative and positive results, and NRL status for the laboratories in Category B.

| acceptable | z scores, false                                | e negative and  |   |   | status for t                    | ne laborator             | ies in catego            | гу в.  |
|------------|--|---|---|---|---------------------------------|--------------------------|--------------------------|--------|
| Lab code   | No. of<br>compulsory<br>pesticides<br>detected | Compulsory<br>pesticides<br>detected in<br>test item, % | Analysed of<br>compulsory<br>pesticides<br>on Target<br>List, % | No. Of<br>voluntary<br>pesticides<br>detected | No. of<br>acceptable<br>z score | No. of false<br>negative | No. of false<br>positive | NRL    |
| 51)        | 14   | 93  | 92  | 0   | 13                              | 1                        | 1                        | NRL-CF |
| 6          | 12   | 80  | 80  | 0   | 12                              | 0                        | 0                        |        |
| 91)        | 14   | 93  | 95  | 0   | 14                              | 1                        | 1                        | NRL-FE |
| 10         | 15   | 100   | 100   | 2   | 11                              | 0                        | 2                        |        |
| 11         | 10   | 67  | 61  | 0   | 9                               | 0                        | 0                        |        |
| 12         | 11   | 73  | 47  | 0   | 11                              | 0                        | 0                        |        |
| 17         | 0  | 0   | 96  | 0   | 0                               | 0                        | 0                        |        |
| 18         | 11   | 73  | 46  | 0   | 11                              | 0                        | 0                        |        |
| 19         | 10   | 67  | 58  | 1   | 2                               | 0                        | 0                        |        |
| 20         | 3  | 20  | 4   | 1   | 1                               | 0                        | 0                        |        |
| 25         | 3  | 20  | 5   | 1   | 3                               | 0                        | 0                        |        |
| 32         | 3  | 20  | 4   | 1   | 3                               | 0                        | 0                        |        |
| 33         | 11   | 73  | 46  | 1   | 10                              | 0                        | 0                        |        |
| 351)       | 15   | 100   | 100   | 2   | 13                              | 0                        | 1                        | NRL-CF |
| 38         | 5  | 33  | 59  | 0   | 5                               | 0                        | 0                        |        |
| 39         | 10   | 67  | 75  | 0   | 10                              | 0                        | 0                        |        |
| 40         | 4  | 27  | 27  | 0   | 4                               | 0                        | 0                        | NRL-CF |
| 41         | 7  | 47  | 59  | 0   | 6                               | 3                        | 0                        |        |
| 42         | 11   | 73  | 71  | 0   | 9                               | 0                        | 0                        | NRL-CF |
| 47         | 5  | 33  | 34  | 0   | 4                               | 2                        | 0                        |        |
| 48         | 3  | 20  | 4   | 0   | 3                               | 0                        | 0                        |        |
| 49         | 10   | 67  | 68  | 1   | 9                               | 1                        | 0                        |        |
| 53         | 13   | 87  | 78  | 1   | 5                               | 1                        | 0                        |        |
| 54         | 6  | 40  | 60  | 1   | 6                               | 1                        | 0                        | NRL-CE |
| 57         | 3  | 20  | 4   | 0   | 3                               | 0                        | 0                        | NRL-FE |
| 59         | 3  | 20  | 4   | 0   | 3                               | 0                        | 0                        |        |
| 60         | 3  | 20  | 4   | 1   | 3                               | 0                        | 0                        |        |
| 65         | 3  | 20  | 5   | 0   | 3                               | 0                        | 0                        |        |
| 67         | 11   | 73  | 56  | 0   | 9                               | 0                        | 0                        | NRL-FE |
| 68         | 9  | 60  | 59  | 1   | 8                               | 0                        | 0                        |        |
| 69         | 12   | 80  | 95  | 0   | 10                              | 3                        | 0                        |        |
| 721)       | 14   | 93  | 100   | 2   | 13                              | 1                        | 1                        |        |

| Lab code          | No. of<br>compulsory<br>pesticides<br>detected | Compulsory<br>pesticides<br>detected in<br>test item, % | Analysed of<br>compulsory<br>pesticides<br>on Target<br>List, % | No. Of<br>voluntary<br>pesticides<br>detected | No. of<br>acceptable<br>z score | No. of false<br>negative | No. of false<br>positive | NRL    |
|-------------------|--|---|---|---|---------------------------------|--------------------------|--------------------------|--------|
| 73                | 10   | 67  | 81  | 0   | 8                               | 1                        | 1                        |        |
| 74                | 10   | 67  | 76  | 1   | 10                              | 1                        | 0                        |        |
| 76                | 12   | 80  | 91  | 0   | 9                               | 0                        | 0                        |        |
| 77                | 0  | 0   | 8   | 1   | 0                               | 0                        | 0                        |        |
| 80                | 11   | 73  | 44  | 1   | 11                              | 0                        | 0                        |        |
| 83                | 11   | 73  | 76  | 0   | 8                               | 0                        | 0                        |        |
| 85                | 11   | 73  | 45  | 0   | 10                              | 0                        | 0                        |        |
| 86                | 11   | 73  | 83  | 1   | 11                              | 0                        | 0                        |        |
| 87                | 12   | 80  | 74  | 1   | 9                               | 1                        | 0                        |        |
| 91                | 13   | 87  | 100   | 0   | 12                              | 2                        | 1                        |        |
| 92                | 11   | 73  | 45  | 0   | 10                              | 0                        | 0                        |        |
| 96                | 8  | 53  | 50  | 1   | 8                               | 0                        | 0                        |        |
| 97                | 14   | 93  | 99  | 1   | 11                              | 1                        | 1                        | NRL-CF |
| 105               | 9  | 60  | 79  | 1   | 9                               | 4                        | 0                        |        |
| 109               | 13   | 87  | 88  | 0   | 12                              | 1                        | 0                        | NRL-FE |
| 112               | 12   | 80  | 86  | 1   | 8                               | 1                        | 6                        | NRL-CF |
| 118               | 10   | 67  | 73  | 1   | 10                              | 0                        | 0                        |        |
| 120               | 13   | 87  | 88  | 3   | 13                              | 1                        | 0                        |        |
| 125               | 6  | 40  | 24  | 0   | 6                               | 0                        | 0                        |        |
| 127               | 8  | 53  | 50  | 0   | 8                               | 0                        | 0                        |        |
| 128               | 7  | 47  | 76  | 1   | 4                               | 1                        | 2                        |        |
| 131               | 11   | 73  | 43  | 0   | 10                              | 0                        | 0                        |        |
| 135               | 10   | 67  | 79  | 1   | 9                               | 2                        | 2                        |        |
| 136               | 7  | 47  | 24  | 0   | 7                               | 0                        | 0                        |        |
| 144               | 10   | 67  | 53  | 0   | 8                               | 0                        | 0                        |        |
| 145               | 6  | 40  | 36  | 1   | 4                               | 1                        | 0                        |        |
| 152               | 0  | 0   | 44  | 0   | 0                               | 0                        | 0                        |        |
| 153               | 14   | 93  | 98  | 0   | 12                              | 1                        | 1                        |        |
| 154               | 13   | 87  | 81  | 1   | 3                               | 1                        | 4                        | NRL-CF |
| 156               | 5  | 33  | 48  | 0   | 5                               | 0                        | 0                        |        |
| 157 <sup>1)</sup> | 14   | 93  | 100   | 2   | 14                              | 1                        | 1                        |        |
| 160               | 0  | 0   | 51  | 0   | 0                               | 0                        | 0                        |        |
| 161               | 15   | 100   | 100   | 0   | 14                              | 4                        | 0                        | NRL-CE |
| 162               | 12   | 80  | 100   | 0   | 9                               | 7                        | 3                        |        |
| 163               | 15   | 100   | 100   | 2   | 15                              | 1                        | 0                        |        |

<sup>1)</sup> Laboratories that reported false positive results and consequently were moved from Category A to Category B

#### 3.4 Trends in numbers of participating laboratories and their performance

The number of EU and EFTA laboratories participating in the EUPTs on cereals has increased steadily until EUPT-CF10 where the highest number of laboratories participated. After this, the number has settled at around 150, unless the Test Item is a feed. Then the number of participants drops. The numbers from EUPT-CF10 and forward can be seen in **Table 13**.

**Table 13.** Overall trends in participation of laboratories, pesticides in the target list and test item, and performance of laboratories in the 7 latest EUPTs cereals (excluding EUPT-CF7 on feed and EUPT-CF12 on hay).

| PT and types of test item                   | EUPT-<br>CF10<br>Rye<br>flour | EUPT-<br>CF11<br>Oat<br>flour | EUPT-<br>CF12<br>Hay<br>flour | EUPT-<br>CF13<br>Oat<br>kernels | EUPT-<br>CF14<br>Rice<br>kernels | EUPT-<br>CF15<br>Rapeseed<br>cake | EUPT-<br>CF16<br>Barley<br>Kernels |
|---|-------------------------------|-------------------------------|-------------------------------|---------------------------------|----------------------------------|-----------------------------------|------------------------------------|
| Participants submitting results (EU+EFTA)   | 160                           | 149                           | 111                           | 149                             | 156                              | 129                               | 151                                |
| MRM pesticides in the Target Pesticide List | 134/7                         | 153/9                         | 155/23                        | 160/32                          | 164/38                           | 172/41                            | 169/53                             |
| MRM pesticides in the test material         | 16                            | 18                            | 8                             | 18                              | 19                               | 22                                | 19                                 |
| No. of results for MRM pesticides           | 2012                          | 2172                          | 808                           | 2007                            | 2298                             | 1315                              | 2206                               |
| Average of 'reported results', %            | 79                            | 83                            | 74                            | 75                              | 80                               | 83                                | 78                                 |
| Range of 'reported results', %              | 58-90                         | 65-93                         | 40-91                         | 44-94                           | 26-93                            | 57-93                             | 32-97                              |
| Acceptable z scores, %                      | 95                            | 89                            | 93                            | 93                              | 91                               | 87                                | 89                                 |
| Questionable z scores, %                    | 2                             | 3                             | 3                             | 3.1                             | 3                                | 7                                 | 4                                  |
| Unacceptable z scores, %                    | 2                             | 8                             | 3                             | 3.4                             | 6                                | 6                                 | 6                                  |
| False negatives, %                          | 2                             | 4                             | 1                             | 2.3                             | 3.4                              | 2.0                               | 3.3                                |
| Number of false positives                   | 0                             | 19                            | 7                             | 3                               | 14                               | 9                                 | 25                                 |
| Category A, % of participating laboratories | 53                            | 45                            | 51                            | 57                              | 57                               | 57                                | 59                                 |
| Good AZ <sup>2</sup> , %                    | 93                            | 92                            | 92                            | 91                              | 91                               | 67                                | 83                                 |
| Satisfactory AZ <sup>2</sup> , %            | 5                             | 1.5                           | 3.4                           | 5.7                             | 6.7                              | 12                                | 7                                  |
| Unsatisfactory AZ², %                       | 2                             | 6.2                           | 5.1                           | 3.4                             | 2.2                              | 22                                | 10                                 |
| Alg A RSD%                                  | 17                            | 17                            | 20                            | 18                              | 19                               | 30                                | 22                                 |

The number of pesticides included in the Target Pesticide List has also increased during this 16-year period, from 43 to 169 compulsory compounds and 53 voluntary compounds. Thus, the demands put on the participating laboratories has increased every year. Many laboratories have a limited scope and are therefore not able to cover all pesticides in the PT. In this EUPT, 23% of the laboratories were not able to analyse and detect more than 70% of pesticides present in the Test Item. Last EUPT the number was 18% and the year before it was also 25%. So no improvement was seen on this issue. The analytical scope was in average 82%.

The overall analytical performance (accuracy of measurement) if looking at the percentage of acceptable, questionable, unacceptable z scores has decreased during the last 2 EUPTs, and in EUPT-CF 16 only 89% of the results were acceptable and more questionable/unacceptable z scores were seen. The average percent of reported results in the last seven EUPT-CF has been between 74-83%. The false negative results have fluctuated between 1-4%. No decrease in false positive results has been seen and the highest number was seen in EUPT-CF16.

The percentage of Category A laboratories has increased slightly over the years and 59% of the participants were evaluated as Category A in this EUPT. For Category A the percentage of participant with AZ $^2$  was <2 (good) has been >90% for many year. However, for the rapeseed cake EUPT-CF15 this dropped significantly to 67% and in this EUPT it was 83%. This might be due to the relatively high content of some of the pesticides, giving challenges for the quantitation if the result are outside the calibration curve and the extract must be diluted.

#### 3.5 Summary, conclusions and prospects for the EUPTs on pesticide residues in cereals

The EUPT-CF16 Test Item was barley kernels containing incurred and spiked pesticides. The barley kernels have been sprayed in the field with commercially available pesticide formulations and additionally spiked post-harvest in the laboratory. The final Test Item contained the following pesticides: azoxystrobin, benzovindiflupyr, boscalid, dieldrin, endosulfan-beta, fenpicoxamid, fenpropathrin, fluopyram, fluxapyroxad, lambda-cyhalothrin, lindane, mefentrifluconazole, metconazole, proquinazid, prothioconazole-desthio, pyraclostrobin, tau-fluvalinate and teflubenzuron.

One hundred fifty-one EU and EFTA laboratories, from 29 different countries agreed to participate in this proficiency test. Three of them did not report any results due to different reasons. An additional 9 laboratories from EU candidate states and Third Countries registered for the PT and all submitted results. The Target Pesticide List distributed to the laboratories prior to the test contained 169 individual compulsory and 53 voluntary compounds.

The number of false positives and false negatives has varied between the EUPTs. Twenty-fiv false positive results were reported and the number of false negatives represented 3.3% of the total number of results. This is at levels typically seen in EUPT-CFs. The average Alg A-RSD was at 22%, lower than the FFP-RSD of 25%.

Of the reported results for the evaluated pesticides, more than 90% acceptable azoxystrobin, benzovindiflupyr, boscalid, endosulfan-beta, fenpropathrin, fluopyram, fluxapyroxad, proquinazid, pyraclostrobin and taufluvalinate. For aclonifen, dieldrin, fenpicoxamid, lambda-cyhalothrin, lindane, mefentrifluconazole, metconazole, prothioconazole-desthio and teflubenzuron between 73-89% of the results were acceptable.

The Test Item for EUPT-CF17 will be wheat kernel, and will planned to be shipped to the laboratories in February 2023. The selection of pesticides will continue to be focused on pesticides included in the scope of the EU multi-annual coordinated control programme, the working document as well as additional pesticides of relevance to feed and/or cereal production in Europe and in other parts of the world from where significant quantities of feed and cereals are imported.

#### 4. ACKNOWLEDGEMENTS

The Organisers wish to thank the members of the EURL Quality Control and Advisory Groups for their valuable advice.

#### 5. REFERENCES

- [1] Regulation (EU) No 2017/625 of the European Parliament and of the Council on official controls and other official activities performed to ensure the application of food and feed law, rules on animal health and welfare, plant health and plant protection products. Published at 0J of the EU L 95/1 of 07.04.2017
- [2] Regulation (EC) No 396/2005, published at OJ of the EU L70 of 16.03.2005, as last amended by Regulation 839/2008 published at OJ of the EU L234 of 30.08.2008.
- [3] CEN EN 15662 Foods of plant origin Determination of pesticide residues using GC-MS and/or LC-MS/MS following acetonitrile extraction/partitioning and clean-up by dispersive SPE QuEChERS-method
- [4] Thompson M., Ellison S. L. R. and Wood R., The International Harmonized Protocol for the Proficiency Testing of Analytical Chemistry Laboratories. Pure & Appl. Chem., Vol.78, No. 1, pp. 145-196, 2006.
- [5] ISO 13528:2015 Statistical methods for use in proficiency testing by interlaboratory comparison
- [6] ISO 17043:2010 Conformity assessment -- General requirements for proficiency testing

### **APPENDICES**

**Appendix 1** List of laboratories registered to participate in the EUPT-CF15

### Participating labs from EU and EFTA member states

| Country        | Institution  | City               | NRL    | Results reported |
|----------------|--|--------------------|--------|------------------|
| Austria        | AGES Innsbruck - Institute for Food Safety (PLMA)                                  | Innsbruck          | NRL-CF |                  |
| Belgium        | Sciensano - Pesticide Lab  | Brussels           | NRL-CF |                  |
| Belgium        | Primoris Belgium   | Gent - Zwijnaarde  |        |                  |
| Belgium        | FLVV - Tervuren  | Tervuren           |        |                  |
| Bulgaria       | Primoris - Bulgaria, Plovdiv   | Plovdiv            |        |                  |
| Bulgaria       | CLCTC - Sofia   Pesticide Lab  | Sofia              |        |                  |
| Croatia        | Bioinstitut d.o.o., Cakovec  | Cakovec            |        |                  |
| Croatia        | Teaching Institute of Public Health of Primorsko-goranska<br>County                | Rijeka             |        |                  |
| Croatia        | Sample Control d.o.o.  | Lucko              |        |                  |
| Croatia        | INSPECTO d.o.o. Laboratorij (Osijek)   | Osijek             |        |                  |
| Croatia        | CROATIAN VETERINARY INSTITUTE  | Zagreb             | NRL-FE |                  |
| Croatia        | Dr. Andrija Štampar - Pesticide Lab  | Zagreb             | NRL-CE |                  |
| Croatia        | Eurofins Croatiakontrola   | Zagreb             |        |                  |
| Cyprus         | Pesticide Lab - Department of Agriculture  | Nicosia            | NRL-FE |                  |
| Cyprus         | SGL - Pesticide Lab (Nicosia)  | Nicosia            | NRL-CE |                  |
| Czech Republic | UKZUZ - Czech Republic, Brno   | Brno               | NRL-FE |                  |
| Czech Republic | Czech Agriculture and Food Inspection Authority (CAFIA)                            | Praha              | NRL-CE |                  |
| Czech Republic | Metrological and Testing Laboratory, University of Chemistry and Technology Prague | Praha              |        |                  |
| Denmark        | Laboratoriet Ringsted - Pesticide Lab  | Ringsted           | NRL-FE |                  |
| Estonia        | Agricultural Research Center - Estonia, Saku                                       | Saku               | NRL-CF |                  |
| Finland        | Finnish Customs Laboratory   | Espoo              | NRL-CE |                  |
| Finland        | Chemistry Unit / Finnish Food Authority  | Helsinki           | NRL-FE |                  |
| France         | GIRPA  | Beaucouzé          |        |                  |
| France         | CERECO (GARONS)  | GARONS             |        |                  |
| France         | CAPINOV  | Landerneau         |        |                  |
| France         | INOVALYS   | Le Mans            |        |                  |
| France         | Laboratoire SCL de PARIS   | Massy Cedex        | NRL-CF |                  |
| France         | SCL (Montpellier)  | Montpellier        |        |                  |
| France         | Phytocontrol (Nimes) - Pesticide Lab   | Nimes              |        |                  |
| France         | CAMP Méditerrannée (Perpignan)   | PERPIGNAN          |        |                  |
| Germany        | BVL Unit 504 NRL for Pesticide Residues  | Berlin             | NRL-CF |                  |
| Germany        | LUA Sachsen - Pesticide Lab, Dresden   | Dresden            |        |                  |
| Germany        | Bayerisches Landesamt für Gesundheit und Lebensmittelsi-<br>cherheit               | Erlangen           |        |                  |
| Germany        | Central Institute of the Bundeswehr Medical Service Munich                         | Garching-Hochbrück |        |                  |
| Germany        | Landesamt für Verbraucherschutz St Fachbereich 3                                   | Halle/Saale        |        |                  |
| Germany        | LLG  | Halle/Saale        |        |                  |
| Germany        | Eurofins Dr. Specht Laboratorien GmbH  | Hamburg            |        |                  |
| Germany        | TLLLR  | Jena               |        |                  |
| Germany        | LTZ Augustenberg   | Karlsruhe          |        |                  |
| Germany        | LHL - Pesticide Lab (Kassel)   | Kassel             |        |                  |

| Country | Institution  | City             | NRL    | Results reported |
|---------|--|------------------|--------|------------------|
| Germany | AGROLAB LUFA GmbH  | Kiel             |        |                  |
| Germany | Chemisches und Veterinäruntersuchungsamt Rhein-Ruhr-<br>Wupper                       | Krefeld          |        |                  |
| Germany | Landeslabor Schleswig-Holstein   | Neumünster       |        |                  |
| Germany | Staatliche Betriebsgesellschaft für Umwelt und Landwirtschaft - FB42                 | Nossen           |        |                  |
| Germany | Niedersächsisches Landesamt für Verbraucherschutz und Le-<br>bensmittelsicherheit    | Oldenburg        |        |                  |
| Germany | Landeslabor Berlin-Brandenburg, Potsdam  | Potsdam          |        |                  |
| Germany | Landesamt für Landwirtschaft, Lebensmittelsicherheit und Fischerei                   | Rostock          |        |                  |
| Germany | LUA Saarland - Pesticide Lab   | Saarbrücken      |        |                  |
| Germany | LUA Rheinland-Pfalz, Institut für LM-Chemie Speyer                                   | Speyer           |        |                  |
| Germany | LUFA Speyer  | Speyer           |        |                  |
| Germany | LAVES - Pesticide Lab (Stade)  | Stade            |        |                  |
| Germany | Labor Friedle - Germany, Tegernheim  | Tegernheim       |        |                  |
| Germany | Landesbetrieb Hessisches Landeslabor, Standort Wiesbaden                             | Wiesbaden        |        |                  |
| Greece  | General Chemical State Laboratory  | Athens           | NRL-CF |                  |
| Greece  | Laboratory of pesticide residues   | Kifissia         | NRL-CF |                  |
| Greece  | Pesticide Residue Laboratory of RCPPQPC of Thessaloniki                              | Thessaloniki     |        |                  |
| Hungary | Food Chain Safety Centre Non-profit Ltd. Pesticide Residue<br>Analytical Laboratory  | Hódmezovásárhely |        |                  |
| Hungary | FCSCN Ltd., Pesticide Residue Analytical Laboratory, Miskolc                         | Miskolc          |        |                  |
| Hungary | Food Chain Safety Centre Non-profit Ltd., Pesticide Residue<br>Analytical Laboratory | Szolnok          |        |                  |
| Hungary | NFCSO Pesticide NRL, Velence   | Velence          | NRL-CF |                  |
| Iceland | Matis ohf  | Reykjavík        | NRL-CF |                  |
| Ireland | The Food Chemistry Laboratories - DAFM   | Co. Kildare      | NRL-CF |                  |
| Italy   | Polo di specializzazione alimenti-DAP Bari-ARPA Puglia                               | Bari             |        |                  |
| Italy   | laboratorio di Prevenzione ATS Bergamo   | Bergamo          |        |                  |
| Italy   | APPA Bolzano   | Bolzano          |        |                  |
| Italy   | Istituto Zooprofilattico Sperimentale Lombardia Emilia Ro-<br>magna                  | Brescia          |        |                  |
| Italy   | ARPA-ER - Pesticide Lab  | Ferrara          |        |                  |
| Italy   | Azienda Sanitaria Locale di Firenze  | Firenze          |        |                  |
| Italy   | IZS PB - Foggia  | Foggia           |        |                  |
| Italy   | IZS PLVA - Lab. Chimico Liguri   | Genova           | NRL-CF |                  |
| Italy   | ARPAL - Dipartimento Laboratorio - UO Analisi Chimiche e Fisiche                     | La Spezia        |        |                  |
| Italy   | ARPA Lazio (sez. Latina) - Pesticide Lab   | Latina           |        |                  |
| Italy   | IZSVe - Pesticide Lab  | Legnaro (Padova) |        |                  |
| Italy   | ARPAM - Pesticide Lab  | Macerata         |        | No               |
| Italy   | ATS Milano - Laboratorio di Prevenzione  | Milano           |        |                  |
| Italy   | IZS Sicilia - Pesticide Lab  | Palermo          |        |                  |
| Italy   | Istituto Zooprofilattico Sperimentale Umbria Marche "Togo<br>Rosati"                 | Perugia          |        |                  |
| Italy   | IZS Mezzogiorno  | Portici (NA)     |        | No               |
| Italy   | Istituto Superiore di Sanità Dip. DAMSA Rep. ECASS/AN                                | Roma             |        |                  |
| Italy   | IZSLT - Chemistry Department   | Roma             |        |                  |
| Italy   | ARPA VDA - Pesticide Lab   | Saint Christophe |        |                  |
| Italy   | IZS Sardegna - Pesticide Lab   | Sassari          |        |                  |
| Italy   | IZSAM - Pesticide Lab  | Teramo           |        |                  |

| Country     | Institution   | City                     | NRL    | Results reported |
|-------------|---|--------------------------|--------|------------------|
| Italy       | APPA Trento   | Trento                   |        |                  |
| Italy       | ARPA FVG - Pesticide Lab (Udine)  | Udine                    |        |                  |
| Latvia      | BIOR (Riga) - Pesticide Lab   | Riga                     | NRL-CF |                  |
| Lithuania   | National Food and Veterinary Risk Assesssment Institute   | Vilnius                  | NRL-CF |                  |
| Luxembourg  | LNS Food lab  | Dudelange                | NRL-CE |                  |
| Netherlands | Groen Agro Control  | Delfgauw                 |        |                  |
| Netherlands | Eurofins Lav Zeeuw-Vlaanderen B.v.  | Graauw                   |        |                  |
| Netherlands | Dr. A. Verwey B.V.  | Rotterdam                |        |                  |
| Netherlands | NofaLab B.V.  | Schiedam                 |        |                  |
| Netherlands | Wageningen Food Safety Research   | Wageningen               | NRL-CF |                  |
| Norway      | NIBIO, Pesticides and Natural Products Chemistry  | ÅS                       | NRL-CF |                  |
| Poland      | Laboratory of Food & Feed Safety in Bialystok   | Bialystok                |        |                  |
| Poland      | WIW ZHW (Bialystok) - Pesticide Lab   | Bialystok                |        |                  |
| Poland      | Hamilton UO-Technologia   | Grójec                   |        |                  |
| Poland      | WIW ZHW (Katowice) - Pesticide Lab  | Katowice                 |        |                  |
| Poland      | VSES Opole - Pesticide Lab  | Opole                    |        |                  |
| Poland      | Zaklad Higieny Weterynaryjnej w Opolu   | Opole                    |        |                  |
| Poland      | IPP-NRI - Pesticide Lab (Poznan)  | Poznan                   | NRL-FE |                  |
| Poland      | WIW ZHW Poznan Pesticide Lab  | Poznan                   |        |                  |
| Poland      | VSES Rzeszow - Pesticide Lab  | Rzeszow                  |        |                  |
| Poland      | InHort (Skierniewice) - Pesticide Lab   | Skierniewice             |        |                  |
| Poland      | IPP-NRI - Pesticide Lab (Sosnicowice)   | Sosnicowice              |        |                  |
| Poland      | WIW ZHW (Szczecin) - Pesticide Lab  | Szczecin                 |        |                  |
| Poland      | Wojewódzka Stacja Sanitarno-Epidemiologiczna w Warszawie  | Warszaw                  | NRL-CE |                  |
| Poland      | Zaklad Higieny Weterynaryjnej w Warszawie Pracownia Badan<br>Chemicznych                                    | Warszaw                  | THE CE |                  |
| Poland      | WIW ZHW (Wroclaw) - Pesticide Lab   | Wroclaw                  |        |                  |
| Portugal    | Laboratório Regional de Veterinária e Segurança Alimentar   | Funchal - Madeira Island |        |                  |
| Portugal    | INIAV Pesticide Lab (Vegetable & Animal Products)   | Vairão - Vila do Conde   | NRL-CE |                  |
| Romania     | LSVSA B-N   | Bistrita                 |        |                  |
| Romania     | IISPV   | Bucharest                | NRL-CF |                  |
| Romania     | National Phytosanitary Authority-Laboratory for Pesticide Residues Control in Plants and Vegetable Products | Bucharest                | NRL-CF |                  |
| Romania     | Sanitary Veterinary and Food Safety Directorate Bucharest   | Bucharest                |        |                  |
| Romania     | Sanitary Vet. and Food Safety Dir. Cluj Napoca  | Cluj Napoca              |        |                  |
| Romania     | Pesticides Residues Laboratory - Dolj   | Craiova                  |        |                  |
| Romania     | LRCRPPPV (Tirgu Mures) - Pesticide Lab  | Tirgu Mures              |        |                  |
| Slovakia    | Pesticide Lab of PHA SR - Bratislava  | Bratislava               |        |                  |
| Slovakia    | SLAPe OLC UKSUP   | Bratislava               |        |                  |
| Slovakia    | Veterinary and Food Institute in Bratislava   | Bratislava               | NRL-CF |                  |
| Slovenia    | Centralni laboratorij   | Ljubljana                | NRL-FE |                  |
| Slovenia    | Pesticide Lab - Maribor   | Maribor                  | NRL-CE |                  |
| Spain       | Laboratorio Agrario y Fitopatológico de Galicia   | Abegondo. A Coruña       |        |                  |
| Spain       | Instituto Tecnologico de Canarias, S. A. Laboratorio de Residuos. Departamento de Analisis Ambiental        | Agüimes, Gran Canaria    |        |                  |
| Spain       | Laboratorio Analítico Bioclínico - Spain, Almeria   | Almeria                  |        |                  |
| Spain       | Analytica Alimentaria GmbH - Almeria, Spain   | Almeria                  |        |                  |
| Spain       | Laboratorio Salud Pública de Badajoz  | Badajoz                  |        |                  |

| Country     | Institution   | City                            | NRL    | Results reported |
|-------------|---|---------------------------------|--------|------------------|
| Spain       | Laboratori Agència de Salut Pública de Barcelona      | Barcelona                       |        |                  |
| Spain       | Agrolab Ibérica SLU                                   | Burgos                          |        |                  |
| Spain       | LABORATORIO AGRARIO REGIONAL de CASTILLA y LEÓN       | Burgos                          |        |                  |
| Spain       | Labs & Technological Services AGQ - Burguillos        | Burguillos                      |        |                  |
| Spain       | Laboratori Agroalimentari - Generalitat Valenciana    | Burjassot, Valencia             |        |                  |
| Spain       | LAC - Generalitat de Catalunya                        | Cabrils                         |        |                  |
| Spain       | Laboratorio Agroalimentario de Extremadura            | Cáceres                         |        |                  |
| Spain       | LABORATORIO DE SALUD PUBLICA DE CUENCA-SPAIN          | Cuenca                          |        |                  |
| Spain       | Laboratorio Regional de la CCAA de La Rioja           | Logroño                         |        |                  |
| Spain       | EUROFINS ECOSUR, S.A.                                 | LORQUI - MURCIA                 |        |                  |
| Spain       | Laboratorio de Salud Pública de Galicia, Lugo         | Lugo                            |        |                  |
| Spain       | Laboratorio Arbitral Agroalimentario (MAPA, Spain)    | Madrid                          | NRL-CF |                  |
| Spain       | SALUD PUBLICA (LSP - MADRID SALUD)                    | Madrid                          |        |                  |
| Spain       | National Centre for Food (Majadahonda)                | Majadahonda                     | NRL-CF |                  |
| Spain       | LABORATORIO KUDAM, S.L.                               | Pilar de la Horadada (Alicante) |        |                  |
| Spain       | National Center for Technology and Food Safety (CNTA) | San Adrián (Navarra)            |        |                  |
| Spain       | LARAGA - Pesticide Lab (Toledo)                       | Toledo                          |        |                  |
| Spain       | Ainia (Valencia)                                      | Valencia                        |        |                  |
| Spain       | Nasertic - Spain, Villava                             | Villava                         |        |                  |
| Spain       | Laboratorio Agroambiental de Zaragoza                 | Zaragoza                        |        |                  |
| Sweden      | Eurofins Food and Feed Testing Sweden AB              | Lidköping                       |        |                  |
| Switzerland | Kantonales Laboratorium Zürich                        | Zürich                          |        |                  |

### Participating labs from EU candidate states and other non EU countries

| Country   | Institution  | City              | Report<br>data |
|-----------|--|-------------------|----------------|
| Australia | Symbio Laboratories - Australia, Eight Mile Plains   | Eight Mile Plains | No             |
| Brazil    | MAPA - Pesticide Lab - Brazil, Pedro Leopoldo  | Pedro Leopoldo    |                |
| Peru      | Bureau Veritas - Peru, Lima  | LIMA - CALLAO     |                |
| Serbia    | Field Test ltd.  | Belgrade          |                |
| Serbia    | Gradski zavod za javno zdravlje Beograd  | Belgrade          |                |
| Serbia    | Ministry of Agriculture, Forestry and Water Management - Directorate for national reference laboratories | Belgrade          |                |
| Thailand  | Central Laboratory - Pesticide Lab (Bangkok)   | Bangkok           | No             |
| UK        | AFBI   | Belfast           |                |
| UK        | Fera Science Limited   | York              |                |

### Appendix 2

Target Pesticide List

| Target Pesticide List  |              |
|--|--------------|
| Pesticides   | MRRL (mg/kg) |
| Compulsory Compounds (will be considered in Category A/B classification) |              |
| 2-phenylphenol   | 0.01         |
| Acephate   | 0.01         |
| Aclonifen  | 0.01         |
| Acetamiprid  | 0.01         |
| Acrinathrin  | 0.01         |
| Aldrin   | 0.005        |
| Ametoctradin   | 0.01         |
| Azinphos-methyl  | 0.005        |
| Azoxystrobin   | 0.01         |
| Bifenthrin   | 0.01         |
| Biphenyl   | 0.01         |
| Bitertanol   | 0.01         |
| Bixafen  | 0.01         |
| Boscalid   | 0.01         |
| Bromuconazole  | 0.01         |
| Buprofezin   | 0.01         |
| Cadusafos  | 0.005        |
| Carbaryl   | 0.005        |
| Carbendazim  | 0.01         |
| Carbofuran   | 0.005        |
| Carbofuran-3-hydroxy   | 0.005        |
| Carboxin   | 0.01         |
| Chlorantraniliprole  | 0.01         |
| Chlorfenapyr   | 0.01         |
| Chlorfenvinphos  | 0.01         |
| Chlorpropham   | 0.01         |
| Chlorpyrifos   | 0.005        |
| Chlorpyrifos-methyl  | 0.01         |
| Clothianidin   | 0.01         |
| Cyantraniliprole   | 0.01         |
| Cyazofamid   | 0.01         |
| Cyfluthrin   | 0.01         |
| Cymoxanil  | 0.01         |
| Cypermethrin   | 0.01         |
| Cyproconazole  | 0.01         |
| Cyprodinil   | 0.01         |
| Deltamethrin   | 0.01         |
| Demeton-S-methylsulfone  | 0.005        |
| Diazinon   | 0.005        |
| Dichlorvos   | 0.005        |
| Dieldrin   | 0.005        |
| Difenoconazole   | 0.01         |

| Pesticides              | MRRL (mg/kg) |
|-------------------------|--------------|
| Diflubenzuron           | 0.01         |
| Dimethoate              | 0.003        |
| Dimethomorph            | 0.01         |
| Diniconazole            | 0.01         |
| Endosulfan-alpha        | 0.01         |
| Endosulfan-beta         | 0.01         |
| Endosulfan-sulfate      | 0.01         |
| Epoxiconazole           | 0.01         |
| Ethion                  | 0.01         |
| Ethirimol               | 0.01         |
| Ethoprophos             | 0.005        |
| Etoxazole               | 0.01         |
| Famoxadone              | 0.01         |
| Fenbuconazole           | 0.005        |
| Fenhexamid              | 0.01         |
| Fenitrothion            | 0.01         |
| Fenpropathrin           | 0.01         |
| Fenpropidin             | 0.01         |
| Fenpropimorph           | 0.01         |
| Fenpyrazamine           | 0.01         |
| Fenpyroximate           | 0.01         |
| Fenthion                | 0.01         |
| Fenthion-oxon           | 0.01         |
| Fenthion-oxon-sulfone   | 0.01         |
| Fenthion-oxon-sulfoxide | 0.01         |
| Fenthion-sulfone        | 0.01         |
| Fenthion-sulfoxide      | 0.01         |
| Fenvalerate             | 0.01         |
| Fipronil                | 0.004        |
| Fipronil-sulfone        | 0.004        |
| Flonicamid              | 0.01         |
| Flubendiamide           | 0.01         |
| Fludioxonil             | 0.01         |
| Flufenoxuron            | 0.01         |
| Fluopicolide            | 0.01         |
| Fluopyram               | 0.01         |
| Fluquinconazole         | 0.01         |
| Flusilazole             | 0.01         |
| Flutolanil              | 0.01         |
| Flutriafol              | 0.01         |
| Fluxapyroxad            | 0.01         |
| Formetanate             | 0.01         |
| Hexaconazole            | 0.01         |
| Imazalil                | 0.005        |
| Imidacloprid            | 0.01         |

| Indoxacarb         0.0           Iprodione         0.0           Isocarbophos         0.0           Isoprothiolane         0.0           Isoproturon         0.0           Kresoxim-methyl         0.0           Lambda-cyhalothrin         0.0           Lindane         0.0           Linuron         0.0           Malaoxon         0.0           Malathion         0.0           Metaflumizone         0.0           Metaflumizone         0.0           Metanzifos         0.0           Methamidophos         0.0           Methamidophos         0.0           Methomyl         0.0           Metribuzin         0.0           Omethoate         0.0           Oxydemeton-methyl         0.0 | mg/kg)     |
|---|------------|
| Socarbophos   0.0   | )1         |
| Isoprothiolane         0.0           Isoproturon         0.0           Kresoxim-methyl         0.0           Lambda-cyhalothrin         0.0           Lindane         0.0           Linuron         0.0           Malaoxon         0.0           Malathion         0.0           Mandipropamid         0.0           Metaflumizone         0.0           Metaaxyl         0.0           Methoazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metribuzin         0.0           Omethoate         0.0  | 1          |
| Isoproturon         0.0           Kresoxim-methyl         0.0           Lambda-cyhalothrin         0.0           Lindane         0.0           Linuron         0.0           Malaoxon         0.0           Malathion         0.0           Mandipropamid         0.0           Metaflumizone         0.0           Metalaxyl         0.0           Metconazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.0   | 1          |
| Kresoxim-methyl         0.0           Lambda-cyhalothrin         0.0           Lindane         0.0           Linuron         0.0           Malaoxon         0.0           Malathion         0.0           Mandipropamid         0.0           Metaflumizone         0.0           Metalaxyl         0.0           Metconazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.0   | 1          |
| Lambda-cyhalothrin         0.0           Lindane         0.0           Linuron         0.0           Malaoxon         0.0           Malathion         0.0           Mandipropamid         0.0           Metaflumizone         0.0           Metalaxyl         0.0           Metconazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.00  | 1          |
| Lindane       0.0         Linuron       0.0         Malaoxon       0.0         Malathion       0.0         Mandipropamid       0.0         Metaflumizone       0.0         Metalaxyl       0.0         Metconazole       0.0         Methacrifos       0.0         Methamidophos       0.0         Methomyl       0.0         Metolachlor       0.0         Metrafenone       0.0         Metribuzin       0.0         Omethoate       0.00   | 1          |
| Linuron         0.0           Malaoxon         0.0           Malathion         0.0           Mandipropamid         0.0           Metaflumizone         0.0           Metalaxyl         0.0           Metconazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.00   | 1          |
| Malaoxon       0.0         Malathion       0.0         Mandipropamid       0.0         Metaflumizone       0.0         Metalaxyl       0.0         Metconazole       0.0         Methacrifos       0.0         Methamidophos       0.0         Methomyl       0.0         Metrolachlor       0.0         Metrafenone       0.0         Metribuzin       0.0         Omethoate       0.00  | 1          |
| Malathion       0.0         Mandipropamid       0.0         Metaflumizone       0.0         Metalaxyl       0.0         Metconazole       0.0         Methacrifos       0.0         Methamidophos       0.0         Methomyl       0.0         Metolachlor       0.0         Metrafenone       0.0         Metribuzin       0.0         Omethoate       0.00  | 1          |
| Mandipropamid       0.0         Metaflumizone       0.0         Metalaxyl       0.0         Metconazole       0.0         Methacrifos       0.0         Methamidophos       0.0         Methomyl       0.0         Metolachlor       0.0         Metrafenone       0.0         Metribuzin       0.0         Omethoate       0.00  | 1          |
| Metaflumizone         0.0           Metalaxyl         0.0           Metconazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.00  | 1          |
| Metalaxyl       0.0         Metconazole       0.0         Methacrifos       0.0         Methamidophos       0.0         Methomyl       0.0         Metolachlor       0.0         Metrafenone       0.0         Metribuzin       0.0         Omethoate       0.00  | 1          |
| Metconazole         0.0           Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.00  | )1         |
| Methacrifos         0.0           Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.00  | 1          |
| Methamidophos         0.0           Methomyl         0.0           Metolachlor         0.0           Metrafenone         0.0           Metribuzin         0.0           Omethoate         0.00  | 1          |
| Methomyl0.0Metolachlor0.0Metrafenone0.0Metribuzin0.0Omethoate0.00   | 1          |
| Metolachlor0.0Metrafenone0.0Metribuzin0.0Omethoate0.00  | )1         |
| Metrafenone0.0Metribuzin0.0Omethoate0.00  | 1          |
| Metribuzin 0.0 Omethoate 0.00   | )1         |
| Omethoate 0.00  | )1         |
|   | )1         |
| Oxydemeton-methyl 0.00  | 03         |
| 0.00  | <b>)</b> 5 |
| Paclobutrazol 0.0   | 1          |
| Parathion 0.0   | )1         |
| Penconazole 0.0   | 1          |
| Pencycuron 0.0  | )1         |
| Pendimethalin 0.0   | 1          |
| Permethrin 0.0  | )1         |
| Phosphamidon 0.0  | )1         |
| Pirimicarb 0.0  | 1          |
| Pirimicarb-desmethyl 0.0  | 1          |
| Pirimiphos-methyl 0.0   | )1         |
| Prochloraz 0.0  | )1         |
| Procymidone 0.0   | )1         |
| Profenofos 0.0  | 1          |
| Propamocarb 0.0   | )1         |
| Propiconazole 0.0   | 1          |
| Proquinazid 0.0   | )1         |
| Prosulfocarb 0.0  | 1          |
| Prothioconazole-desthio 0.0   | )1         |
| Prothiofos 0.0  | 1          |
| Pymetrozine 0.0   | )1         |
| Pyraclostrobin 0.0  |            |
| Pyridaben 0.0   | 1          |

| Pesticides  | MRRL (mg/kg) |
|---|--------------|
| Pyridalyl   | 0.01         |
| Pyrimethanil  | 0.01         |
| Pyriproxyfen  | 0.01         |
| Quinoxyfen  | 0.01         |
| Spinetoram  | 0.01         |
| Spirodiclofen   | 0.01         |
| Spiromesifen  | 0.01         |
| Spirotetramat   | 0.01         |
| Spirotetramat metabolite BYI08330-enol                                      | 0.01         |
| Spiroxamine   | 0.01         |
| Sulfoxaflor   | 0.01         |
| Tau-Fluvalinate   | 0.01         |
| Tebuconazole  | 0.01         |
| Tebufenozide  | 0.01         |
| Teflubenzuron   | 0.01         |
| Tefluthrin  | 0.01         |
| Terbuthylazine  | 0.01         |
| Tetraconazole   | 0.01         |
| Tetradifon  | 0.01         |
| Tetramethrin  | 0.01         |
| Thiabendazole   | 0.01         |
| Thiacloprid   | 0.01         |
| Thiamethoxam  | 0.01         |
| Thiodicarb  | 0.01         |
| Thiophanate-methyl  | 0.01         |
| Tolclofos-methyl  | 0.01         |
| Triadimefon   | 0.01         |
| Triadimenol   | 0.01         |
| Triflumizole  | 0.01         |
| Triflumizole metabolite (FM-6-1)  | 0.01         |
| Triazophos  | 0.005        |
| Tricyclazole  | 0.01         |
| Trifloxystrobin   | 0.01         |
| Trifluralin   | 0.01         |
| Triticonazole   | 0.01         |
| Vinclozolin   | 0.01         |
| Zoxamide  | 0.01         |
|   |              |
| Voluntary Compounds (will not be considered in Category A/B classification) |              |
| Benalaxyl (sum)   | 0.01         |
| Benzovindiflupyr  | 0.01         |
| Chlordane-cis   | 0.01         |
| Chlordane-oxy   | 0.01         |
| Chlordane-trans   | 0.01         |
| Chlorfluazuron  | 0.01         |
| Clomazone   | 0.01         |

| Pesticides             | MRRL (mg/kg) |
|------------------------|--------------|
| Cyflufenamid           | 0.01         |
| Cyhalofop-butyl        | 0.01         |
| DDD-pp                 | 0.01         |
| DDE-pp                 | 0.01         |
| DDT-op                 | 0.01         |
| DDT-pp                 | 0.01         |
| Dinotefuran            | 0.01         |
| Endrin                 | 0.01         |
| Endrin-ketone          | 0.01         |
| Fenobucarb             | 0.01         |
| Fenpicoxamid           | 0.01         |
| Florpyrauxyfen-benzyl  | 0.01         |
| Fluensulfone           | 0.01         |
| Flufenacet             | 0.01         |
| Flupyradifurone        | 0.01         |
| Flutianil              | 0.01         |
| HCH-alpha              | 0.01         |
| HCH-beta               | 0.01         |
| Heptachlor             | 0.01         |
| Heptachlorepoxid-cis   | 0.01         |
| Heptachlorepoxid-trans | 0.01         |
| Isofetamid             | 0.01         |
| Isopyrazam             | 0.01         |
| Isoxaflutole           | 0.01         |
| Mefentrifluconazole    | 0.01         |
| Metconazole            | 0.01         |
| Molinate               | 0.01         |
| Novaluron              | 0.01         |
| Oxadiargyl             | 0.01         |
| Oxathiapiprolin        | 0.01         |
| Oxyfluorfen            | 0.01         |
| Penflufen              | 0.01         |
| Pentachloro-aniline    | 0.01         |
| Penthiopyrad           | 0.01         |
| Picolinafen            | 0.01         |
| Propaquizafop          | 0.01         |
| Pyrethrins             | 0.01         |
| Pyridate               | 0.01         |
| Pyriofenone            | 0.01         |
| Quinalphos             | 0.01         |
| Quinoclamine           | 0.01         |
| Quintozene             | 0.01         |
| Rotenone               | 0.01         |
| Tolfenpyrad            | 0.01         |
| Tri-allate             | 0.01         |
| Tritosulfuron          | 0.01         |

Appendix 3 Homogeneity data

|            | Aclo      | nifen     | Azoxys    | strobin   | Benzovir  | ndiflupyr | Bose      | calid     |
|------------|-----------|-----------|-----------|-----------|-----------|-----------|-----------|-----------|
|            | mg        | /kg       | mg        | /kg       | mg        | /kg       | mg        | /kg       |
| Sample no. | Portion 1 | Portion 2 |
| 004        | 0.062     | 0.054     | 1.139     | 1.137     | 0.817     | 0.817     | 1.686     | 1.847     |
| 022        | 0.063     | 0.051     | 1.201     | 1.269     | 0.855     | 0.862     | 2.047     | 1.816     |
| 034        | 0.062     | 0.062     | 1.316     | 1.273     | 0.949     | 0.928     | 1.946     | 1.912     |
| 051        | 0.045     | 0.059     | 1.387     | 1.318     | 0.945     | 0.944     | 2.247     | 2.172     |
| 077        | 0.068     | 0.051     | 1.193     | 1.180     | 0.888     | 0.862     | 1.889     | 1.976     |
| 094        | 0.054     | 0.062     | 1.187     | 1.179     | 0.867     | 0.872     | 1.944     | 1.789     |
| 118        | 0.058     | 0.069     | 1.281     | 1.177     | 0.892     | 0.855     | 2.008     | 1.867     |
| 178        | 0.055     | 0.063     | 1.184     | 1.261     | 0.841     | 0.885     | 1.897     | 1.885     |
| 199        | 0.065     | 0.045     | 1.314     | 1.205     | 0.897     | 0.819     | 2.054     | 1.796     |
| 224        | 0.052     | 0.062     | 1.268     | 1.264     | 0.918     | 0.925     | 1.869     | 1.830     |
| 245        | 0.047     | 0.066     | 1.216     | 1.300     | 0.891     | 0.929     | 2.028     | 1.787     |

|            | Dieldrin  |           | Endosulfan-beta |           | Fenpropathrin |           | Fenpicoxamid |           |  |
|------------|-----------|-----------|-----------------|-----------|---------------|-----------|--------------|-----------|--|
|            | mg        | mg/kg     |                 | mg/kg     |               | mg/kg     |              | mg/kg     |  |
| Sample no. | Portion 1 | Portion 2 | Portion 1       | Portion 2 | Portion 1     | Portion 2 | Portion 1    | Portion 2 |  |
| 004        | 0.025     | 0.024     | 0.034           | 0.035     | 0.090         | 0.088     | 0.037        | 0.036     |  |
| 022        | 0.022     | 0.024     | 0.031           | 0.033     | 0.081         | 0.085     | 0.036        | 0.036     |  |
| 034        | 0.025     | 0.026     | 0.034           | 0.035     | 0.085         | 0.091     | 0.036        | 0.039     |  |
| 051        | 0.028     | 0.028     | 0.039           | 0.039     | 0.099         | 0.099     | 0.043        | 0.038     |  |
| 077        | 0.026     | 0.024     | 0.036           | 0.033     | 0.093         | 0.085     | 0.037        | 0.035     |  |
| 094        | 0.025     | 0.023     | 0.033           | 0.033     | 0.089         | 0.086     | 0.036        | 0.035     |  |
| 118        | 0.027     | 0.027     | 0.038           | 0.036     | 0.102         | 0.099     | 0.040        | 0.038     |  |
| 178        | 0.022     | 0.025     | 0.031           | 0.032     | 0.077         | 0.086     | 0.031        | 0.033     |  |
| 199        | 0.024     | 0.022     | 0.032           | 0.030     | 0.084         | 0.078     | 0.039        | 0.033     |  |
| 224        | 0.025     | 0.024     | 0.034           | 0.033     | 0.087         | 0.085     | 0.034        | 0.036     |  |
| 245        | 0.027     | 0.026     | 0.036           | 0.035     | 0.092         | 0.092     | 0.037        | 0.035     |  |

|            | Fluopyram |           | Fluxapyroxad |           | Lamba-cyhalothrin |           | Lindane   |           |
|------------|-----------|-----------|--------------|-----------|-------------------|-----------|-----------|-----------|
|            | mg        | /kg       | mg/kg        |           | mg/kg             |           | mg/kg     |           |
| Sample no. | Portion 1 | Portion 2 | Portion 1    | Portion 2 | Portion 1         | Portion 2 | Portion 1 | Portion 2 |
| 004        | 1.345     | 1.436     | 0.661        | 0.716     | 0.068             | 0.069     | 0.038     | 0.037     |
| 022        | 1.469     | 1.379     | 0.786        | 0.687     | 0.071             | 0.075     | 0.034     | 0.035     |
| 034        | 1.448     | 1.466     | 0.751        | 0.726     | 0.075             | 0.075     | 0.037     | 0.037     |
| 051        | 1.540     | 1.513     | 0.834        | 0.811     | 0.086             | 0.089     | 0.041     | 0.038     |
| 077        | 1.371     | 1.413     | 0.717        | 0.739     | 0.076             | 0.072     | 0.038     | 0.035     |
| 094        | 1.398     | 1.408     | 0.739        | 0.681     | 0.073             | 0.068     | 0.038     | 0.036     |
| 118        | 1.438     | 1.328     | 0.756        | 0.703     | 0.073             | 0.075     | 0.042     | 0.040     |
| 178        | 1.429     | 1.321     | 0.722        | 0.700     | 0.070             | 0.077     | 0.030     | 0.038     |
| 199        | 1.486     | 1.375     | 0.771        | 0.696     | 0.070             | 0.065     | 0.036     | 0.033     |
| 224        | 1.348     | 1.381     | 0.716        | 0.703     | 0.074             | 0.076     | 0.036     | 0.038     |
| 245        | 1.445     | 1.350     | 0.772        | 0.679     | 0.079             | 0.078     | 0.037     | 0.038     |

|            | Mefentrifluconazole Metco |           | nazole    | Proquinazid |           | Prothioconazole-desthio |           |           |
|------------|---------------------------|-----------|-----------|-------------|-----------|-------------------------|-----------|-----------|
|            | mg                        | /kg       | mg        | /kg         | mg/kg     |                         | mg/kg     |           |
| Sample no. | Portion 1                 | Portion 2 | Portion 1 | Portion 2   | Portion 1 | Portion 2               | Portion 1 | Portion 2 |
| 004        | 0.705                     | 0.734     | 0.038     | 0.037       | 0.193     | 0.190                   | 1.101     | 1.220     |
| 022        | 0.734                     | 0.775     | 0.031     | 0.035       | 0.204     | 0.206                   | 1.252     | 1.176     |
| 034        | 0.710                     | 0.827     | 0.029     | 0.037       | 0.216     | 0.213                   | 1.253     | 1.284     |
| 051        | 0.871                     | 0.791     | 0.042     | 0.036       | 0.225     | 0.231                   | 1.345     | 1.339     |
| 077        | 0.741                     | 0.773     | 0.036     | 0.034       | 0.206     | 0.204                   | 1.195     | 1.310     |
| 094        | 0.747                     | 0.756     | 0.035     | 0.035       | 0.206     | 0.206                   | 1.219     | 1.228     |
| 118        | 0.714                     | 0.716     | 0.039     | 0.039       | 0.203     | 0.205                   | 1.231     | 1.140     |
| 178        | 0.734                     | 0.727     | 0.030     | 0.031       | 0.208     | 0.201                   | 1.229     | 1.138     |
| 199        | 0.717                     | 0.689     | 0.032     | 0.030       | 0.219     | 0.219                   | 1.253     | 1.117     |
| 224        | 0.759                     | 0.733     | 0.032     | 0.031       | 0.222     | 0.226                   | 1.186     | 1.218     |
| 245        | 0.783                     | 0.758     | 0.038     | 0.035       | 0.220     |                         | 1.274     | 1.137     |

|            | Pyraclo   | strobin   | Tau-Flu   | valinate  |
|------------|-----------|-----------|-----------|-----------|
|            | mg        | /kg       | mg        | /kg       |
| Sample no. | Portion 1 | Portion 2 | Portion 1 | Portion 2 |
| 004        | 4.476     | 4.568     | 0.299     | 0.298     |
| 022        | 4.657     | 4.733     | 0.320     | 0.325     |
| 034        | 4.754     | 5.103     | 0.342     | 0.316     |
| 051        | 5.369     | 4.930     | 0.392     | 0.382     |
| 077        | 4.603     | 4.668     | 0.339     | 0.327     |
| 094        | 4.549     | 4.796     | 0.315     | 0.316     |
| 118        | 4.726     | 4.382     | 0.328     | 0.317     |
| 178        | 4.554     | 4.405     | 0.318     | 0.344     |
| 199        | 4.801     | 4.415     | 0.303     | 0.306     |
| 224        | 4.556     | 4.473     | 0.360     | 0.356     |
| 245        | 4.730     | 4.633     | 0.314     | 0.285     |


#### Appendix 4 Stability figures


The stability test was performed according to ISO 13528 Annex B [5]. Two different storage temperatures were used; room temperature and -18 °C.

The dates of testing were as follows:

Day 1: 28 March 2022 Day 2: 13 April 2022 Day 3: 2 May 2022


All pesticides passed the test at -18 °C see **1.6 Stability test**. At room temperature indoxacarb and pirimicarb did not pass the test when stored for 11 weeks.

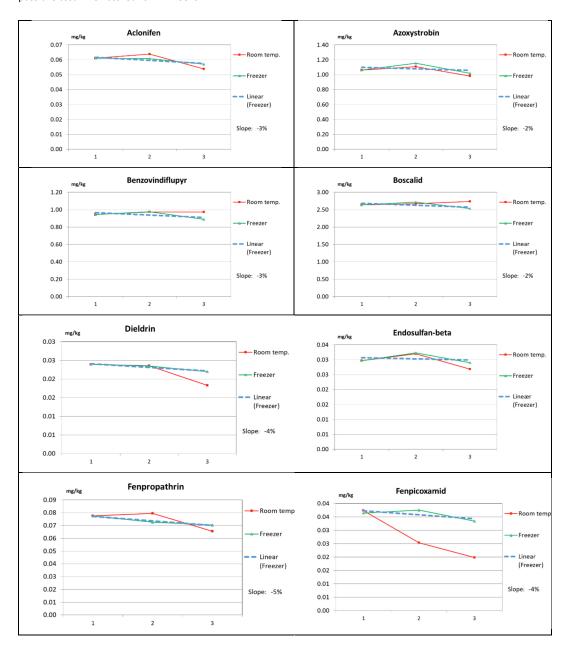

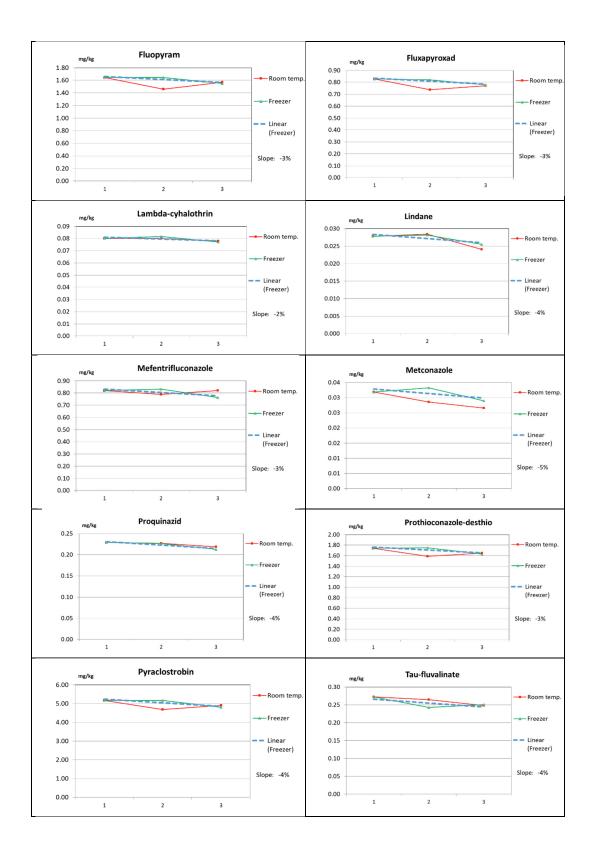


Appendix 5 Graphical presentation of z-scores




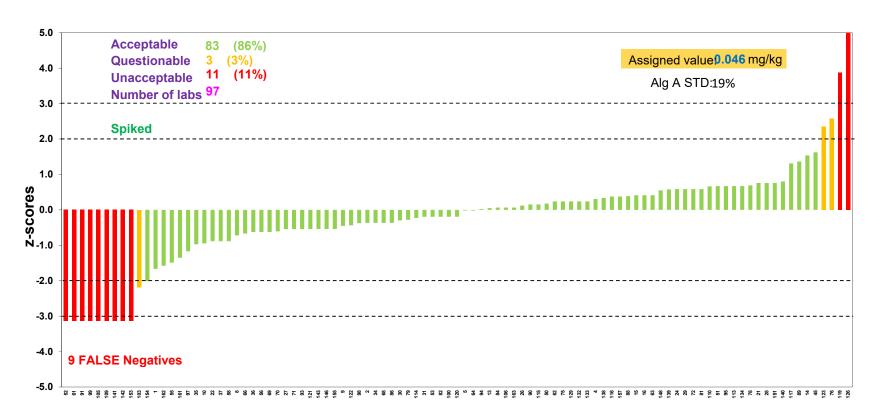


# Azoxystrobin

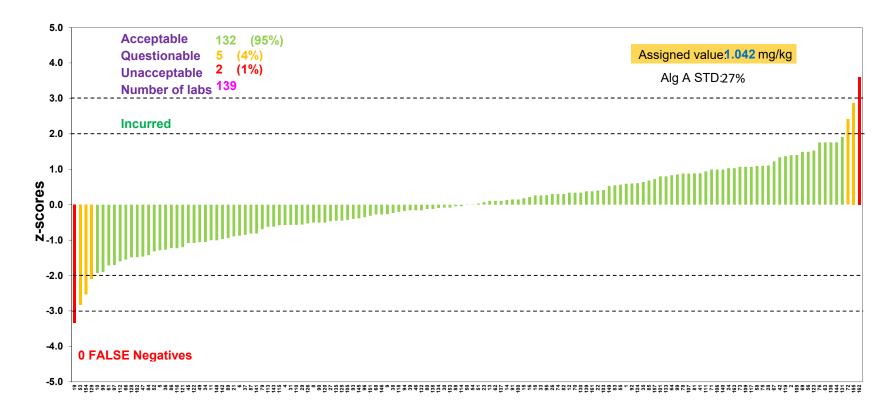

# Benzov in diflupyr

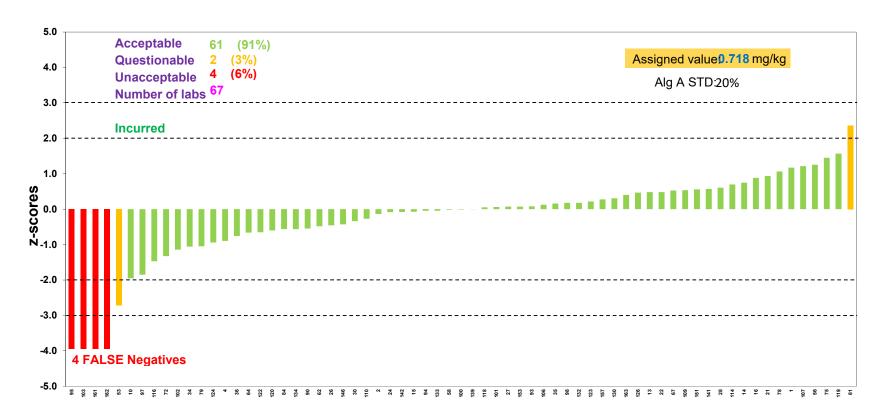

## **Boscalid**

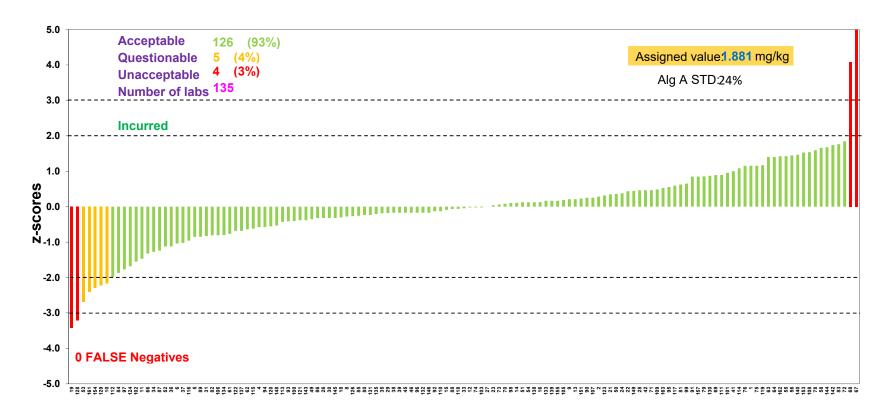

### Dieldrin




### **Endosulfan-beta**






# **Fenpicoxamid**

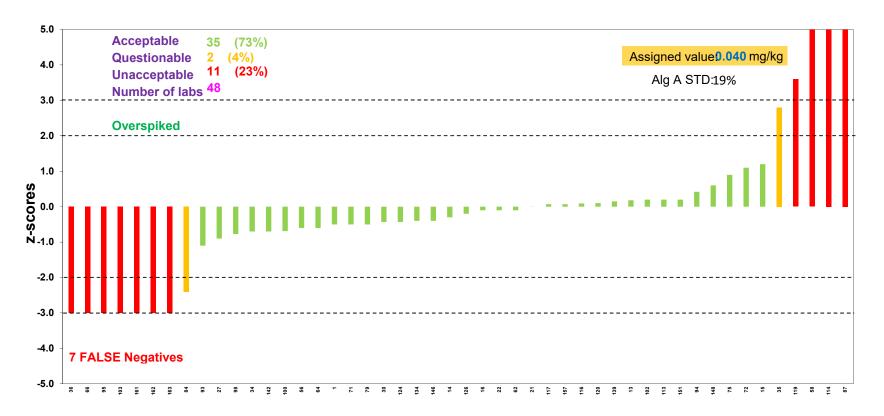

# **Fenpropathrin**



# Fluopyram

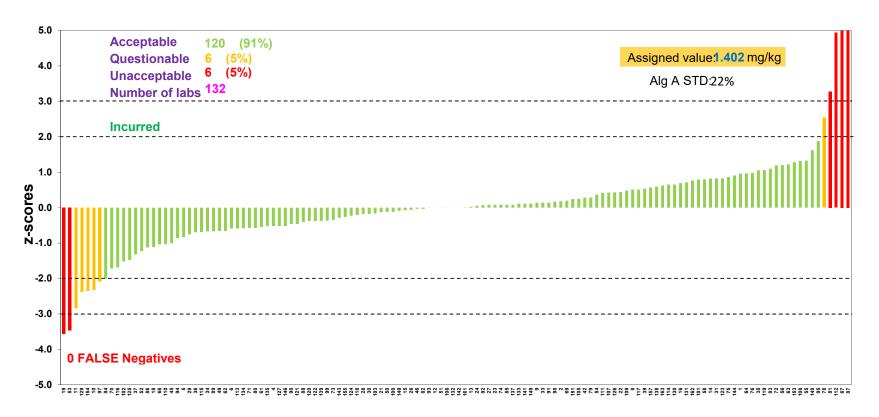

# Fluxapyroxad




# Lambda-cyhalothrin



### Lindane



## Mefentrifluconazole

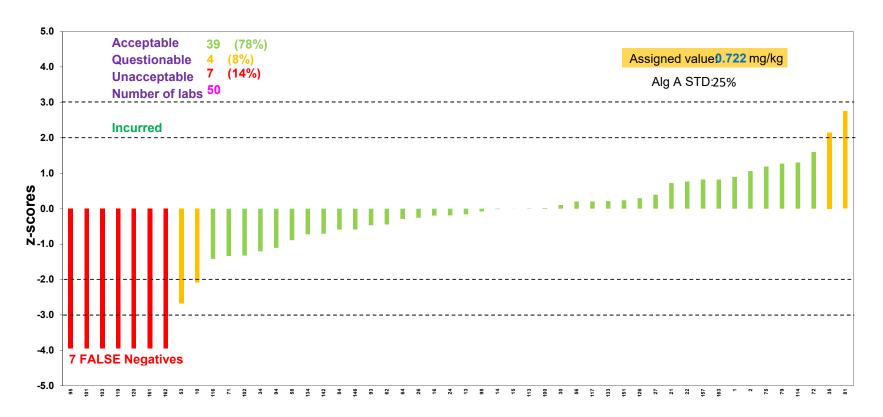

### Metconazole



# **Proquinazid**



### Prothioconazole-desthio EU and EFTA Laboratories




# **Pyraclostrobin**

#### **EU and EFTA Laboratories**



# Tau-Fluvalinate

#### **EU and EFTA Laboratories**

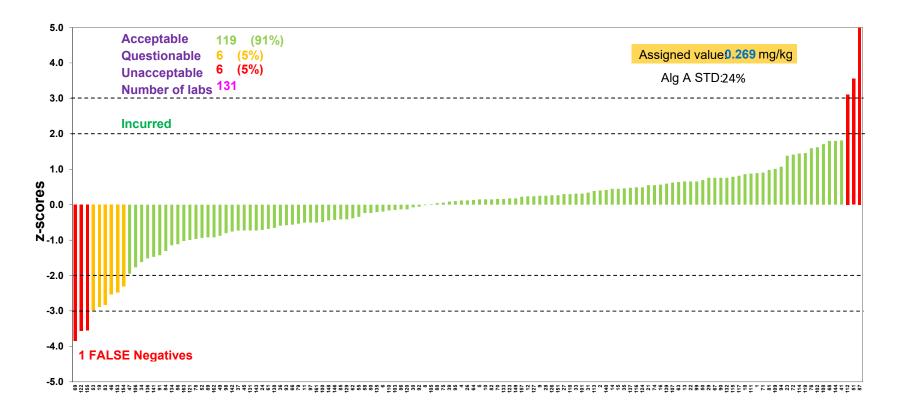

# Teflubenzuron

#### **EU and EFTA Laboratories**





### **GENERAL PROTOCOL**

# for EU Proficiency Tests on Pesticide Residues in Food and Feed

#### Introduction

This protocol contains general procedures valid for all European Union Proficiency Tests (EUPTs) organised on behalf of the European Commission, DG-SANTE<sup>1</sup> by the four European Union Reference Laboratories (EURLs) responsible for pesticide residues in food and feed. These EUPTs are directed at laboratories belonging to the Network<sup>2</sup> of National Reference Laboratories (NRLs) and Official Laboratories (OfLs) of the EU Member States. OfLs from EFTA countries and EU-Candidate countries are also welcome to participate in the EUPTs. OfLs from Third countries may be permitted to participate on a case-by-case basis.

The following four EURLs for pesticide residues were appointed by DG-SANTE based on regulation 882/2004/EC that was repealed by regulation 625/2017/EC<sup>3</sup>:

- EURL for Fruits and Vegetables (EURL-FV),
- EURL for Cereals and Feedingstuffs (EURL-CF),
- EURL for Food of Animal Origin and Commodities with High Fat Content (EURL-AO) and
- EURL for pesticides requiring Single Residue Methods (EURL-SRM).

The aim of these EUPTs is to obtain information regarding the quality, accuracy and comparability of pesticide residue data in food and feed reported to the European Union within the framework of the national control programmes and the EU multiannual co-ordinated control programme<sup>4</sup>. Participating laboratories will be provided with an assessment of their analytical performance that

<sup>&</sup>lt;sup>1</sup> DG-SANTE = European Commission, Health and Food Safety Directorate-General

<sup>&</sup>lt;sup>2</sup> For more information about the EURL/NRL/OfL-Network please refer to the EURL-Web-portal under: "<a href="http://www.eurl-pesticides.eu">http://www.eurl-pesticides.eu</a>"

<sup>&</sup>lt;sup>3</sup> Regulation (EU) 2017/625 of the European Parliament and of the Council on official controls and other official activities performed to ensure the application of food and feed law, rules on animal health and welfare, plant health and plant protection products.. Published at OJ of the EU L95 of 07.04.2017

<sup>&</sup>lt;sup>4</sup> European Commission Proficiency Tests for Pesticide Residues in Fruits and Vegetables, Trends in Analytical Chemistry, 2010, 29 (1), 70 – 83.



they can use to demonstrate their analytical performance and compare themselves with other participating laboratories.

#### **EUPT-Organisers and Scientific Committee**

EUPTs are organised by individual EURLs, or by more than one EURL, in collaboration.

An **Organising Team** (in the following named Organisers) is appointed by the EURL(s) in charge. This team is responsible for all administrative and technical matters concerning the organisation of the PT, e.g. the PT-announcement, the production of the PT-material (Test Item), the undertaking of homogeneity and stability tests, the packing and shipment of the PT-materials, the handling and evaluation of the results and method information submitted by the participants, the drafting of the preliminary and final reports as well as generation and distribution of EUPT-participation certificates.

To complement the internal expertise of the EURLs, a group of external consultants forming the **EUPT-Scientific Committee** (EUPT-SC)<sup>5</sup> has been established and approved by DG-SANTE. The EUPT-SC consists of expert scientists with many years of experience in PTs and/or pesticide residue analysis. The actual composition of the EUPT-SC and the affiliation of each of its members is shown on the EURL-Website. The members of the EUPT-SC are also listed in the Specific Protocol and the Final Report of each EUPT.

The EUPT-SC is made up of the following two subgroups:

- a) An independent Quality Control Group (EUPT-QCG) and
- b) An Advisory Group (EUPT-AG).

The EUPT-SC's role is to help the Organisers make decisions regarding the EUPT design: the selection of the commodity, the selection of pesticides to be included in the Target Pesticide List (see below), the establishment of the Minimum Required Reporting Levels (MRRLs), the statistical treatment and evaluation of the participants' results (in anonymous form), and the drafting and updating of documents, such as the General and Specific PT Protocols and the Final EUPT-Reports.


<sup>&</sup>lt;sup>5</sup> Link to the List of current members of the EUPT Scientific Committee: http://www.eurl-pesticides.eu/library/docs/allcrl/EUPT-SC.pdf



The EUPT-QCG has the additional function of supervising the quality of EUPTs and of assisting the EURLs in confidential aspects such as the choice of the pesticides to be present in the Test Item and the approximate concentrations at which they should be present.

The EUPT-SC typically meets once a year, after the EUPTs of all four pesticide EURLs have been conducted, to discuss the evaluation of the EUPT-results and to assist the EURLs in their decision making. Upcoming EUPTs are also planned during these meetings.

The EUPT-Organising Team and the EUPT-SC together form the **EUPT-Panel**.

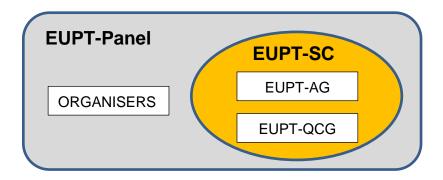

The decisions of the EUPT-Panel will be documented.

This present EUPT General Protocol was jointly drafted by the EUPT-SC and the EURLs.

#### **EUPT Participants**

Within the European Union all NRLs operating in the same area as the organising EURL, as well as all OfLs whose scope overlaps with that of the EUPT, are legally obliged to participate in EUPTs. The legal obligation of NRLs and OfLs to participate in EUPTs arises from:

- Art 38 (b) of Reg. 625/2017/EC and Art. 28 of Reg. 396/2005/EC<sup>6</sup> (for all OfLs analysing for pesticide residues within the framework of official controls<sup>7</sup> of food or feed)
- Art. 101 (1)(a) of Reg. 625/2017/EC (for all NRLs)

<sup>6</sup> Regulation (EC) No 396/2005, published at OJ of the EU L70 of 16.03.2005, as last amended by Regulation 839/2008 published at OJ of the EU L234 of 30.08.2008.

81

<sup>&</sup>lt;sup>7</sup> Official controls in the sense of Reg. 625/2017/EC. This includes labs involved in controls within the framework of national and/or EU-controlled programmes as well as labs involved in import controls according to Regulation 669/2009/EC.



The four EURLs will annually issue and distribute, via the EURL-website, a joint list of all OfLs that must participate in each of the EUPTs to be conducted within a given year. The list of obliged labs will be updated every year to take account of any changes in the lab profiles. Interim updates will be issued to eliminate any possible errors.

NRLs are responsible for checking whether all relevant OfLs within their network are included in the list of obligated laboratories with their actual commodity-scopes and contact information.

OfLs are furthermore urged to keep their own profiles within the EURL-DataPool up-to-date, especially their commodity and pesticide scopes and their contact information.

Labs that are obliged to participate in a given EUPT, and that are not able to participate, must provide the reasons for their non-participation This also applies to any participating laboratories that fail to report results.

OfLs not paying the EUPT sample delivery fee will be initially warned that their participation in subsequent EUPTs could be denied. In case of a repetitive non-payment, the EUPT organisers will inform the corresponding NRL to take action.

#### **Confidentiality and Communication**

The proprietor of all EUPT data is DG-SANTE and as such has access to all information.

For each EUPT, the laboratories are given a unique code (lab code), initially only known to themselves and the Organisers. In the final EUPT-Report, the names of participating laboratories will not be linked to their laboratory codes. It should be noted, however, that the Organisers, at the request by DG-SANTE, may present the EUPT-results on a country-by-country basis. It may therefore be possible that a link between codes and laboratories could be made, especially for those countries where only one laboratory has participated. Furthermore, the EURLs reserve the right to share EUPT results and codes amongst themselves: for example, for the purpose of evaluating overall lab or country performance as requested by DG-SANTE.

As laid down in Regulation 625/2017/EC, NRLs are responsible for evaluating and improving their own OfL-Network. On request from the NRLs, the EURLs will provide them with the PT-codes of the participating OfLs belonging to their OfL-Network. This will allow NRLs to follow the participation and performance of the laboratories within their network.



Communication between participating laboratories during the test, on matters concerning a PT exercise, is not permitted from the start of the PT exercise until the distribution of the preliminary report.

For each EUPT the organising EURL prepares a specific EUPT-Website where all PT-relevant documents in their latest version are linked. In case of important modifications on any of these documents, the participating laboratories will be informed via e-mail. In any case, as soon as the PT-period starts the participants are encouraged to visit the particular EUPT-Website, to make sure that they are using the latest versions of all PT-relevant documents.

The official language used in all EUPTs is English.

#### **Announcement / Invitation Letter**

At least 3 months before the distribution of the Test Item the EURLs will publish an Announcement/Invitation letter on the EURL-web-portal and distribute it via e-mail to the NRL/OfL mailing list available to the EURLs. This letter will inform about the commodity to be used as Test Item, as well as links to the tentative EUPT-Target Pesticide List and the tentative EUPT-Calendar.

#### **Target Pesticide List**

This list contains all analytes (pesticides and metabolites) to be sought for, along with the Minimum Required Reporting Levels (MRRLs) valid for the specific EUPT. The MRRLs are typically based upon the lowest MRLs found either in Regulation 396/2005/EC or Commission Directive 2006/125/EC (Baby Food Directive).

Labs must express their results as stated in the Target Pesticides List.

#### Specific Protocol

For each EUPT the organizing EURL will publish a Specific Protocol at least 2 weeks before the Test Item is distributed to the participating laboratories. The Specific Protocol will contain all the information previously included in the Invitation Letter but in its final version, information on payment and delivery, instructions on how to handle the Test Item upon receipt and on how to submit results, as well as any other relevant information.



#### Homogeneity of the Test Item

The Test Item will be tested for homogeneity typically before distribution to participants. The homogeneity tests usually involve the analysis of two replicate analytical portions, taken from at least ten randomly chosen units of treated Test Item. Both, sample preparation and measurements should be conducted in random order.

The homogeneity test data are statistically evaluated according to ISO 13528, Annex B or to the International Harmonized Protocols jointly published by ISO, AOAC and IUPAC. The results of all homogeneity tests are presented to the EUPT-SC. In special cases, where the above homogeneity test criteria are not met, the EUPT-Panel, considering all relevant aspects (e.g. the homogeneity results of other pesticides spiked at the same time, the overall distribution of the participants' results (CV\*), the analytical difficulties faced during the test, knowledge of the analytical behaviour of the pesticide question), may decide to overrule the test. The reasons of this overruling have to be transparently explained in the Final EUPT-Report. For certain analytes with comparable properties, an equivalent distribution within the sample can be expected if they were spiked/used at simultaneously. The homogeneity test, of one or more of these analytes, may thus be skipped or simplified. If, however, the distribution of participants' results for an analyte that was not or not fully tested for homogeneity, is found to be atypically broad, compared to the tested analytes, the EUPT-SC may decide that a homogeneity test should be performed a posteriori by the EURL.

#### Stability of the analytes contained in the Test Item

The Test Items will also be tested for stability - according to ISO 13528, Annex B. The time delay between the first and the last stability test must exceed the period of the EUPT-exercise. Typically the first analysis is carried out shortly before the shipment of the Test Items and the last one shortly after the deadline for submission of results. To better recognise trends and gain additional certainty one or more additional tests may be conducted by the Organisers. At least 6 sub-samples (analytical portions) should be analysed on each test day (e.g. 2 analytical portions withdrawn from three randomly chosen containers OR 6 portions withdrawn from a single container). In principle all pesticides contained in the Test Item should be checked for stability. However, in individual cases, where sufficient knowledge exists that the stability of a certain analyte is very unlikely to be significantly affected during storage (e.g. based on experience from past stability tests or knowledge of its physicochemical properties), the Organisers, after consultation with the EUPT-QCG, may decide to omit a specific stability test. The EUPT-Panel will finally decide whether





analytes for which the stability test was not undertaken will be included in the Final EUPT-Report, considering all relevant aspects such as the distribution of the participant's results (CV\*).

A pesticide is considered to be adequately stable if  $|y_i-y| \le 0.3 \times \sigma_{pt}$ , with  $y_i$  being the mean value of the results of the last phase of the stability test, y being the mean value of the results of the first phase of the stability test and  $\sigma_{pt}$  being the standard deviation used for proficiency assessment (typically 25 % of the assigned value).

The results of all stability tests are presented to the EUPT-SC. In special cases where the above stability test criteria are not met, the EUPT-SC considering all relevant aspects (e.g. the past experience with the stability of the compound, the overall distribution the participants' results, the measurement variability, analytical difficulties faced during the test and knowledge about the analytical behaviour of the pesticide question) may decide to overrule the test. The reasons of this overruling will be transparently explained in the Final EUPT-Report.

The Organisers may also decide to conduct additional stability tests at different storage conditions than those recommended to the participants e.g. at ambient temperature.

**Stability during shipment:** Considering knowledge about the expected susceptibility of pesticides in the Test Item to possible losses, the Organisers will choose the shipment conditions to be such that pesticide losses are minimised (e.g. shipment of frozen samples, addition of dry ice). As shipment time can differ between labs/countries it is recommended that the Organisers keep track of the shipment duration and then decide whether it is reasonable to conduct additional stability tests at conditions simulating shipment. Should critical losses be detected for certain pesticides, the EUPT-SC will be informed (or the EUPT-QCG before or during the test). Case-by-case decisions may be taken by the EUPT-Panel considering all relevant aspects including the duration and conditions of the shipment to the laboratory as well as the feedback by the laboratory.

#### Methodologies to be used by the participants

Participating laboratories are instructed to use the analytical procedure(s) that they would routinely employ in official control activities (monitoring etc.). Where an analytical method has not yet been established routinely this should be stated.



#### General procedures for reporting results

Participating laboratories are responsible for reporting their own <u>quantitative results</u> to the Organiser within the stipulated deadline. Any pesticide that was targeted by a participating laboratory should be reported as "analysed". Each laboratory will be able to report only <u>one</u> result for each analyte detected in the Test Item. The concentrations of the pesticides detected should be expressed in 'mg/kg' unless indicated otherwise in the specific protocol. Laboratories should not report results below their reporting limits.

#### Correction of results for recovery

Correction of results for recovery is recommended if the average recovery rate significantly deviates from 100 % (typically if outside the 80–120% range). Approaches for recovery correction explicitly stated in the DG-SANTE document are

- a) the use of recovery correction factors,
- b) the use of stable isotope labelled analogues of the target analytes as Internal Standards (ILISs),
- c) the 'procedural calibration' approach as well as
- d) the approach of 'standard addition' with additions of analyte(s) being made to analytical portions.

Results may be corrected for recovery only in cases where this correction is applied in routine practice (including cases of MRL-violations). Laboratories are required to report whether their results were adjusted for recovery and, if a recovery factor was used, the recovery rate (in percentage) must also be reported. If one or more of the approaches b), c) and d) were employed, in which correction for recovery is inherent to the procedures, the apparent recovery figures obtained during validation experiments are not mandatory, and the approached followed are to be reported in the appropriate fields within the data submission tool.

#### **Methodology information**

All laboratories are requested to provide information on the analytical method(s) they have used. A compilation of the methodology information submitted by all participants is presented in an Annex of the Final EUPT-Report or in a separate report. Where necessary the methods are evaluated and discussed, especially in those cases where the result distribution is not unimodal or very broad (e.g.  $CV^* > 35$  %). If no sufficient information on the methodology used is provided, the Organisers reserve the right not to accept the analytical results reported by the participants concerned or even refuse participation in the following PT.



#### **Results evaluation**

The procedures used for the treatment and assessment of results are described below.

#### False Positive results

These are results of pesticides from the Target Pesticides List, that are reported, at or above, their respective MRRL although they were: (i) not detected by the Organiser, even after repeated analyses, and/or (ii) not detected by the overwhelming majority (e.g. > 95 %) of the participating laboratories that had targeted the specific pesticides. In certain instances, case-by-case decisions by the EUPT-Panel may be necessary.

Any results reported lower than the MRRL will not be considered as false positives, even though these results should not have been reported.

#### False Negative results

These are results for pesticides reported by the laboratories as 'analysed' but without reporting numerical values although they were: a) used by the Organiser to treat the Test Item and b) detected by the Organiser as well as the majority of the participants that had targeted these specific pesticides at or above the respective MRRLs. Results reported as '< RL' (RL= Reporting Limit of the laboratory) will be considered as not detected and will be judged as false negatives. In certain instances, case-by-case decisions by the EUPT-Panel may be necessary.

In cases of the assigned value being less than a factor of 3 times the MRRL, false negatives will typically not be assigned. The EUPT-Panel may decide to take case-by-case decisions in this respect after considering all relevant factors such as the result distribution and the reporting limits of the affected labs.

#### Estimation of the assigned value (x<sub>pt</sub>)

In order to minimise the influence of out-lying results on the statistical evaluation, the assigned value  $x_{pt}$  (= consensus concentration) will typically be estimated using the robust estimate of the



participant's mean ( $x^*$ ) as described in ISO 13528:2015<sup>8</sup>, taking into account the results reported by EU and EFTA countries laboratories only. In special justifiable cases, the EUPT-Panel may decide to eliminate certain results traceably associated with gross errors (see "Omission or Exclusion of results" below) or to use only the results of a subgroup consisting of laboratories that have repeatedly demonstrated good performance for the specific or similar compounds in the past.

#### Omission or Exclusion of results

Before estimating the assigned value, results associated with obvious mistakes have to be examined to decide whether they should be removed from the population. Such gross errors may include incorrect recording (e.g. due to transcription errors by the participant, decimal point faults or transposed digits, incorrect unit), calculation errors (e.g. missing factors), analysis of a wrong sample/extract (e.g. a spiked blank), use of wrong concentrations of standard solutions, incorrect data processing (e.g. integration of wrong peak), inappropriate storage or transport conditions (in case of susceptible compounds), and the use of inappropriate analytical steps or procedures that demonstrably lead to significantly biased results (e.g. employing inappropriate internal standards or analytical steps or conditions leading to considerable losses, due to degradations, adsorptions, incomplete extractions, partitioning etc.). Where the Organisers (e.g. after the publication of the preliminary report) receive information of such gross errors, having a significant impact on a generated result, the affected results will be examined on a case-by-case basis to decide whether, or not, they should be excluded from the population used for robust statistics. Results may also be omitted e.g. if an inappropriate method has been used even if they are not outliers. All decisions to omit/exclude results will be discussed with the EUPT-SC and the reasoning for the omission of each result clearly stated in the Final EUPT-Report. However, z scores will be calculated for all results irrespective of the fact that they were omitted from the calculation of the assigned value.

Omitted results might be interesting as they might give indications about possible source(s) of errors. The Organisers will thus ask the relevant lab(s) to provide feedback on possible sources of errors (see also "follow-up activities").

Results reported by laboratories from non EU member states are typically excluded from the population that is used to derive the assigned value (see also "Estimation of the assigned value").

\_

<sup>&</sup>lt;sup>8</sup> DIN ISO 13528:2015, Statistical methods for use in proficiency testing by interlaboratory comparisons, International Organization for Standardization. Therein a specific robust method for determination of the consensus mean and standard deviation without the need for removal of deviating results is described (Algorithm A in Annex C).



#### - Uncertainty of the assigned value

The uncertainty of the assigned values  $u(x_{pt})$  is calculated according to ISO 13528:2015 as:

$$u\left(x_{pt}\right) = 1,25 \times \frac{s^*}{\sqrt{p}}$$

where  $s^*$  is the robust standard deviation and p is the number of results.

In certain cases, and considering all relevant factors (e.g. the result distribution, multimodality, the number of submitted results, information regarding analyte homogeneity/stability, information regarding the use of methodologies that might produce a bias that were used by the participants), the EUPT-Panel may consider the assigned value of a specific analyte to be too uncertain and decide that the results should not be evaluated, or only evaluated for informative purposes. The provisions of ISO 13528:2015 concerning the uncertainty of the assigned value will be taken into account.

#### - Standard deviation of the assigned value (target standard deviation)

The target standard deviation of the assigned value (FFP- $\sigma_{pt}$ ) will be calculated using a Fit-For-Purpose approach with a fixed Relative Standard Deviation (FFP-RSD).

Based on experience from previous EUPTs9, a percentage FFP-RSD of 25 % is currently used for all analyte-matrix combination, with the target standard deviation being calculated as follows:

FFP-
$$\sigma_{pt} = 0.25 \times x_{pt}$$

The EUPT-Panel reserves the right to also employ other FFP-RSDs or other approaches for setting the assigned value on a case-by-case basis, considering analytical difficulties and experience gained from previous proficiency tests.

For informative purposes the robust relative standard deviation (CV\*) of the participants results is calculated according to ISO 13528:2015; Chapter 7.7 following Algorithm A in Annex C (so called "consensus approach").

\_

Omparative Study of the Main Top-down Approaches for the Estimation of Measurement Uncertainty in Multiresidue Analysis of Pesticides in Fruits and Vegetables. J. Agric. Food Chem., 2011, 59(14), 7609-7619.



#### z scores

This parameter is calculated using the following formula:

$$z_i = \frac{\left(x_i - x_{pt}\right)}{FFP - \sigma_{nt}}$$

where  $x_i$  is the value reported by the laboratory,  $x_{pt}$  is the assigned value, and FFP- $\sigma_{pt}$  is the standard deviation using the FFP approach. Z scores will be rounded to one decimal place. For the calculation of combined z scores (see below) the original z scores will be used and the combined z-scores will be rounded to one decimal place after calculation.

Any z scores > 5 will be typically reported as '> 5' and a value of '5' will be used to calculate combined z scores (see below).

Z scores will be interpreted in the following way, as is set in the ISO 17043:2010<sup>10</sup>:

$$|z| \le 2.0$$
 Acceptable  $2.0 < |z| < 3.0$  Questionable  $|z| \ge 3.0$  Unacceptable

For results considered as false negatives, z scores will be calculated using the MRRL or RL (the laboratory's Reporting Limit) if RL < MRRL. Where, using this approach, the calculated z scores for false negatives are > -3 (still questionable), they will be fixed at -3.5 to underline that these are unacceptable results. These z-scores will typically appear in the z-score histograms and used in the calculation of combined z-scores.

#### Collection of measurement uncertainty (MU) figures

The participating labs will be asked to report the MU figure they would routinely report with each EUPT result. The EUPT-Panel will decide whether and how to evaluate these figures and whether indications will be made to the laboratories in this respect.

\_

<sup>&</sup>lt;sup>10</sup> ISO/IEC 17043:2010. Conformity assessment – General requirements for proficiency testing



#### Category classification

The EUPT-Panel will decide if and how to classify the laboratories into categories based on their scope and/or performance. Currently a scope-based classification into Category A and Category B is employed. Laboratories that a) are able to analyse at least 90% of the compulsory pesticides in the target pesticides list, b) have correctly detected and quantified a sufficiently high percentage of the pesticides present in the Test Item (at least 90 %) and c) reported no false positives, will have demonstrated 'sufficient scope' and will be therefore classified into Category A. For the 90% criterion the number of pesticides needed to be correctly analysed to have sufficient scope will be calculated by multiplying the number of compulsory pesticides from the Target Pesticides List by 0.9 and rounding to the nearest full number with 0.5 decimals being rounded downwards (see some examples in Table 1).

Table 1. No. of pesticides from the Target Pesticides List needed to be targeted or pesticides present in the Test Item that need to be correctly detected and quantified to have sufficient scope.

| No. of compulsory pesticides present in the Test Item / Target Pesticides List (N) | 90 % | No. of pesticides needed to be correctly detected and quantified / targeted to have sufficient scope (n) | n     |  |
|--|------|--|-------|--|
| 3  | 2.7  | 3  | N     |  |
| 4  | 3.6  | 4  | IN    |  |
| 5  | 4.5  | 4  |       |  |
| 6  | 5.4  | 5  |       |  |
| 7  | 6.3  | 6  |       |  |
| 8  | 7.2  | 7  |       |  |
| 9  | 8.1  | 8  | NI 1  |  |
| 10   | 9.0  | 9  | N - 1 |  |
| 11   | 9.9  | 10   |       |  |
| 12   | 10.8 | 11   |       |  |
| 13   | 11.7 | 12   |       |  |
| 14   | 12.6 | 13   |       |  |
| 15   | 13.5 | 13   |       |  |
| 16   | 14.4 | 14   |       |  |
| 17   | 15.3 | 15   |       |  |
| 18   | 16.2 | 16   |       |  |
| 19   | 17.1 | 17   | N O   |  |
| 20   | 18   | 18   | N - 2 |  |
| 21   | 18.9 | 19   |       |  |
| 22   | 19.8 | 20   |       |  |
| 23   | 20.7 | 21   |       |  |
| 24   | 21.6 | 22   |       |  |
| 25   | 22.5 | 22   | NI O  |  |
| 26   | 23.4 | 23   | N - 3 |  |

The EUPT-Panel reserves the right to develop and apply alternative classification rules.



#### Overall performance of laboratories - combined z scores

For evaluation of the overall performance of laboratories within Category A, the Average of the Squared z score  $(AZ^2)^{11,12}$  (see below) will be used. The  $AZ^2$  is calculated as follows:

$$AZ^2 = \frac{\sum_{i=1}^n z_i^2}{n}$$

Where n is the number of z scores to be considered in the calculation. In the calculation of the  $AZ^2$ , z scores higher than 5 will be set as 5. Based on the  $AZ^2$  achieved, the laboratories are classified as follows:

$$AZ^2 \le 2.0$$
 Good  
 $2.0 < AZ^2 < 3.0$  Satisfactory  
 $AZ^2 \ge 3.0$  Unsatisfactory

Combined z scores are considered to be of lesser importance than individual z scores. The EUPT-Panel retains the right not to calculate  $AZ^2$  if it is considered as not being useful or if the number of results reported by any participant is considered to be too low.

In the case of EUPT-SRMs, where only a few results per lab may be available, the Average of the Absolute z scores (AAZ) may be calculated for informative purposes, but only for labs that have reported enough results to obtain 5 or more z scores. For the calculation of the AAZ, z scores higher than 5 will also be set as 5. The z-scores appointed to false negatives will be also included in the calculation of the combined z-scores.

Laboratories within Category B will be typically ranked according to the total number of pesticides they correctly reported to be present in the Test Item. The number of acceptable z scores achieved will be presented, too. The EURL-Panel retains the right to calculate combined z scores (see above) also for labs within Category B, e.g. for informative purposes, provided that a minimum number of results (z scores) have been reported.

<sup>&</sup>lt;sup>11</sup> Formerly named "Sum of squared z scores (SZ<sup>2</sup>)"

<sup>&</sup>lt;sup>12</sup> Laboratory assessment by combined z score values in proficiency tests: experience gained through the EUPT for pesticide residues in fruits and vegetables. Anal. Bioanal. Chem., 2010, 397, 3061–3070.



#### **Publication of results**

The EURLs will publish a preliminary report, containing tentative assigned values and z score values for all pesticides present in the Test Item, within 2 months of the deadline for result submission.

The Final EUPT-Report will be published after the EUPT-Panel has discussed the results. Taking into account that the EUPT-Panel meets normally only once a year (typically in late summer or autumn) to discuss the results of all EUPTs organised by the EURLs earlier in the year, the Final EUPT-Report may be published up to 10 months after the deadline for results submission. Results submitted by non-EU/EFTA laboratories might not always be used in the tables or figures in the Final EUPT-Report.

#### **Certificates of participation**

Together with the Final EUPT-Report, the EURL Organiser will deliver a Certificate of Participation to each participating laboratory showing the z scores achieved for each individual pesticide, the combined z scores calculated (if any), and the classification into Categories.

#### **Feedback**

At any time before, during or after the PT participants have the possibility to contact the Organisers and make suggestions or indicate errors. After the distribution of the Final EUPT-Report, participating laboratories will be given the opportunity to give their feedback to the Organisers and make suggestions for future improvements.

#### Correction of errors

Should errors be discovered in any of the documents issued prior to the EUPT (Calendar, Target Pesticides List, Specific Protocol, General Protocol) the corrected documents will be uploaded onto the website and in the case of substantial errors the participants will be informed. **Before starting** the exercise, participants should make sure to download the latest version of these documents.

If substantial errors are discovered in the Preliminary EUPT-Report the Organisers will distribute a new corrected version, where it will be stated that the previous version is no longer valid.



Where substantial errors are discovered in the Final EUPT-Report the EUPT-Panel will decide whether a corrigendum will be issued and how this should look like. The online version of the Final EUPT report will be replaced by the new one and all affected labs will be contacted.

Where errors are discovered in EUPT-Certificates the relevant laboratories will be sent new corrected ones. Where necessary the laboratories will be asked to return the old ones.

#### Follow-up activities

Laboratories are expected to undertake follow-up activities to trace back the sources of erroneous or strongly deviating results (typically those with |z| > 2.0) - including all false positives. In exceptional cases, follow-up activities may even be indicated for results within  $|z| \le 2.0$  (e.g. where two errors with opposed tendency cancel each other leading to acceptable results).

Upon request, the laboratory's corresponding NRL and EURL are to be informed of the outcome of any investigative activities for false positives, false negatives and for results with  $|z| \ge 3.0$ . Concerning z scores between 2.0 and 3.0 the communication of the outcome of follow-up activities is optional but highly encouraged where the source of deviation could be identified and could be of interest to other labs.

According to instructions from DG-SANTE, the "Protocol for management of underperformance in comparative testing and/or lack of collaboration of National Reference Laboratories (NRLs) with EU Reference Laboratories (EURLs) activities" is to be followed.

NRLs will be considered as **underperforming in relation to scope** if in at least two of the last four EUPTs falling within their responsibility area they: a) haven't participated, or b) targeted less than 90% of the compulsory pesticides in the target lists (80% for SRM-compounds), or c) detected less than 90% of the compulsory compounds present in the test items (80% for SRM-compounds). Additionally, NRLs that obtained AZ<sup>2</sup> higher than 3 (AAZ higher than 1.3 for SRM-compounds) in two consecutive EUPTs of the last four EUPTs, will be considered as **underperforming in accuracy**. A two-step protocol established by DG-SANTE will be applied as soon as underperformance of an NRL is detected<sup>13</sup>:

#### Phase 1:

• Identifying the origin of the bad results (failure in EUPTs).

<sup>&</sup>lt;sup>13</sup> Article 101 of Regulation (EC) 625/2017





 Actions: On the spot visits and training if necessary and repetition of the comparative test if feasible and close the assessment of results by the EURL.

#### Phase 2:

- If the results still reveal underperformance the Commission shall be informed officially by the EURL including a report of the main findings and corrective actions.
- The Commission shall inform the Competent Authority and require that appropriate actions are taken.

Underperformance rules for the OfLs will be established at a later stage.

#### **Disclaimer**

The EUPT-Panel retains the right to change any parts of this EUPT – General Protocol based on new scientific or technical information. Any changes will be communicated in due course.





# SPECIFIC PROTOCOL

# for the EU Proficiency Test for Pesticide Residues in Cereals/Feeding stuff using Multi Residue Methods, EUPT-CF16 (2022)

(last updated: 10 March 2022)

#### Introduction

This protocol is complementary to the <u>General Protocol for EU Proficiency Tests for Pesticide Residues in Food and Feed</u> (9th Edition). The current proficiency test covers pesticides that are determined by Multi Residue Methods. This EUPT is to be performed by all National Reference Laboratories for Cereals and/or Feeding stuffs (NRL-CFs) as well as by all official EU laboratories (OfLs) responsible for official pesticide residue controls on feeding stuff, as far as their scope overlaps with that of the EUPT-CF15.

#### **Test Item (Test Material)**

This proficiency test concerns the analysis of pesticide residues in barley kernels. The barley was grown in Denmark and pesticides were applied in the field.

The Organiser, will check the Test Items for sufficient homogeneity and for stability at conditions reproducing sample shipment and storage during the duration of the test, according to ISO 13528, Annex B. All these tests will be conducted by the organiser, the EURL-CF which is (ISO 17025 accredited).

#### **Analytical Parameters**

The Test Item contains several pesticides from the **Target Pesticides List**.

Laboratories must report their results as stated in the Target Pesticides List.

#### **Amount of Test Item**

The participants will receive:

approximately 100 g of barley kernels Test Item with incurred and spiked pesticides

Blank material will not be distributed to the participants.

#### **Shipment of Test Items**

The Test Items are planned to be shipped on 28 March 2022.

Test Items will be shipped frozen and packed in thermo-boxes together with a freezer block. The organiser will aim to ensure that all participating laboratories will receive their shipments on the same day. Prior to shipment a reminder will be sent to the participating laboratories by e-mail.

Laboratories must make their own arrangements for the receipt of the package. They should inform the Organiser of any public holidays in their country/city during the week of the shipment, and must make the necessary arrangements to receive the shipment, even if the laboratory is closed.

#### **Instructions on Test Item Handling**

Once received, the Test Items should be stored deep-frozen (at -18°C or below) before analysis to avoid any possible deterioration/spoilage and to minimize pesticide losses. The test Item should be milled before analysis. After milling, mix the flour thoroughly, before taking the analytical portion(s).

All participants should use their own routine standard operating procedures for milling, extraction, cleanup and analytical measurement and their own reference standards for identification and quantification purposes.

The homogeneity test is conducted using 5 g of milled Test Item in all cases. As sub-sampling variability increases with decreasing analytical portion size, sufficient homogeneity can only be guaranteed where participants employ sample portions that are equal to or larger than the ones stated above.

#### **EUPT Webtool and Deadlines**

To select pesticide scope and report results and method information, the participants should log in to the **EUPT Webtool** using the username and password send by email. Please, save the credentials, as it will be valid for the EUPTs next year.

<u>Selection/deselection of scope</u>: The analytical scope must be selected prior to the shipment of the samples. This is done via the **EUPT Webtool**. The scope selection subpage will be open from 14 March to 28 March 2022. As default all mandatory pesticides are preselected.

<u>Results and method submission</u>: The **EUPT Webtool** will be accessible from 29 March 2022 for sample receipt acknowledgement and submission results and method information.

The deadline for submission is 2 May 2022 at 23.00 CET.

**IMPORTANT:** After the final submission it will NOT be possible to edit the results. Participants will receive an email confirming the submission of their results. Attached to the email will be an excel file with their submitted data.

<u>Test Item Receipt and Acceptance:</u> Once the laboratory has received the Test Items it must report to the organiser, via the **EUPT Webtool**, the date of receipt, and its acceptance. If the laboratory does not respond by 1 April 2022 at 12.00 CET, the Organiser will assume that the Test Items have been received and accepted.

If participants have not received the Test Items by **the 1 April 2022 at noon**, they must inform the Organiser immediately by e-mail to <a href="mailto:eurl-cf@food.dtu.dk">eurl-cf@food.dtu.dk</a>.

#### **Reporting Quantitative Results:**

Results should not be reported where a pesticide

- a) was not detected,
- b) was detected below the RL (Reporting Limit) of the laboratory, or

#### Significant Figures:

Residue levels < 0.010 mg/kg;

to be expressed by two significant figures (e.g. 0.0058 mg/kg).

Residue levels ≥ 0.010 mg/kg;

- to be expressed by three significant figures, e.g. 0.156, 1.64, 10.3 mg/kg.

<u>Reporting Analytical method</u>: The laboratory must to report details of the analytical methods they used. If not it will not be possible to submit results.

#### Reporting of supplementary information in case of false negative results

In case of false negative results, the affected laboratories will be asked to provide details on the methodology used after the deadline for result submission. This has also to be done by accessing EUPT Webtool. Deadline for this is 11 May 2022.

#### Follow-up actions

In accordance with Art. 32 1b of Regulation (EC) No 2017/625, underperformance of any NRL-CF in comparative testing will be followed by EURL-CF.

#### **Documents**

All documents related to EUPT-CF16 can be found on EUPT-CF15 website.

#### https://www.eurl-

<u>pesticides.eu/docs/public/tmplt\_article.asp?LabID=400&CntID=1196&Theme\_ID=1&Pdf=False&Lang=EN</u>

#### Calendar

| Activity   | Dates                   |
|--|-------------------------|
| Announcement Calendar Target Pesticide List                  | December 2021           |
| EUPT-Registration Website open                               | December 2021           |
| Deadline for registration                                    | 28 February 2022        |
| Specific Protocol published                                  | 14 March 2022           |
| Website for selecting pesticide scope open                   | 14 March 2022           |
| Website for selecting pesticide scope closed                 | 28 March 2022           |
| Distribution of Test items                                   | 28 March 2022           |
| Deadline for receipt and acceptance of Test Materials        | within 24 hr on receipt |
| Deadline for Result Submission                               | 2 May 2022              |
| Deadine for Result Submission                                | at 23.00 CET            |
| Deadline for submission of additional method information for | 11 May 2022             |
| false negative results                                       | at 24.00 CET            |
| Preliminary Report (only compilation of results) published   | 11 July 2022            |
| Final Report published                                       | December 2022           |

#### **Participation Fees**

For participating laboratories from the EU, EU-candidate states and EFTA states the participation fee will be:

• 200€

The participation fees for laboratories from third countries will be:

• 350€

For further information, visit <u>www.eurl-pesticides.eu.</u>

#### **Delays in Payment**

The participants will receive an invoice from DTU. The terms of payment are 30 days net. After this deadline reminders will be sent. From the second reminder onwards an administration fee of DKK 100.00 excluding VAT (ca. 13 €) will be charged per reminder.

If the participant ask DTU to issue a new invoice because additional/new information are needed on the invoice, or just want a copy of the original invoice, that may add additional cost due to the administrative workload.

Any questions concerning invoices must be directed to Anne Linda Bisp, <a href="mailto:abi@adm.dtu.dk">abi@adm.dtu.dk</a> at the financial department of DTU.

#### Contact information:

DTU Food

National Food Institute





#### **Mette Erecius Poulsen**

Head of EURL Cereals and Feeding stuff

National Food Institute Technical University of Denmark Kemitorvet, Building 202 DK-2800 Lyngby

Phone: +45-3588-7463 E-Mail: <u>eurl-cf@food.dtu.dk</u> http://www.eurl-pesticides.eu

#### **Organising Team:**

| Elena Hakme, Chemist      | EURL for Cereals and Feeding stuff |
|---------------------------|------------------------------------|
| Ederina Ninga, Chemist    | EURL for Cereals and Feeding stuff |
| Ban M. Kadhum, Technician | EURL for Cereals and Feeding stuff |
| Susanne Pless, Technician | EURL for Cereals and Feeding stuff |

#### **Quality Control Group:**

| Dr. Antonio Valverde | University of Almería, Spain       |
|----------------------|------------------------------------|
| Dr. Paula Medina     | European Food Safety Agency, Italy |

#### **Advisory Group**

| Advisory Group                 |  |
|--------------------------------|--|
| Prof. Amadeo R. Fernández-Alba | University of Almeria, Spain   |
| Dr. Carmen Ferrer Amate        | University of Almeria, Spain   |
| Dr. André de Kok               | Former Wargeningen Food Safety Reseach, Wargeningen, The Netherlands                 |
| Mr. Ralf Lippold               | Chemisches und Veterinäruntersuchungsamt (CVUA) Freiburg, Germany                    |
| Dr. Michelangelo Anastassiades | Chemisches und Veterinäruntersuchungsamt (CVUA) Stuttgart, Germany                   |
| Dr. Hermann Unterluggauer      | AGES Competence Center for Residues of Plant Protection Products, Innsbruck, Austria |
| Dr. Tuija Pihlström            | National Food Administration, Uppsala, Sweden  |
| Dr. Magnus Jezussek            | Bavarian Authority of Health an Food Safety, Erlangen, Germany                       |
| Mr. Finbarr o'Regan            | Pesticide Control Laboratory, Celbridge, Ireland                                     |
| Dr. Patrizia Pelosi            | Istituto Superiore di Sanità, Roma, Italy  |
| Dr. Hans Mol                   | Wargeningen Food Safety Reseach, Wargeningen, The Netherlands                        |
| Mr. Bjoern Hardenbusch         | Chemisches und Veterinäruntersuchungsamt (CVUA) Freiburg, Germany                    |
| Mr. Radim Stepan               | Czech Agriculture and Food Inspection Authority                                      |

National Food Institute
Technical University of Denmark
Building 202
Kemitorvet
2800 Lyngby
Denmark

Tel: +45 35 88 70 00

ISBN: 978-87-7586-005-0

www.food.dtu.dk